


Unison UK-Phenyl Unison UK-C1
Unison UK-C8 Cadenza CX-C18
Unison UK-C18 Cadenza CD-C18

100 x 4.6 mm


HPLC COLUMN APPLICATION

2,7-Bis [N- (1-naphthyl) anilino] -9,9'-spirobi [9H-fluorene] (Organic light-emitting diode material)

2,7-ビス [N -(1-ナフチル) アニリノ] -9,9'-スピロビ [9H-フルオレン]
(有機EL材料)

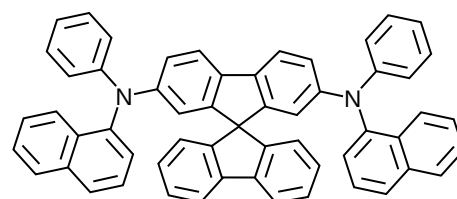

2,7-bis [N- (1-naphthyl) anilino]
-9,9'-spirobi [9H-fluorene]

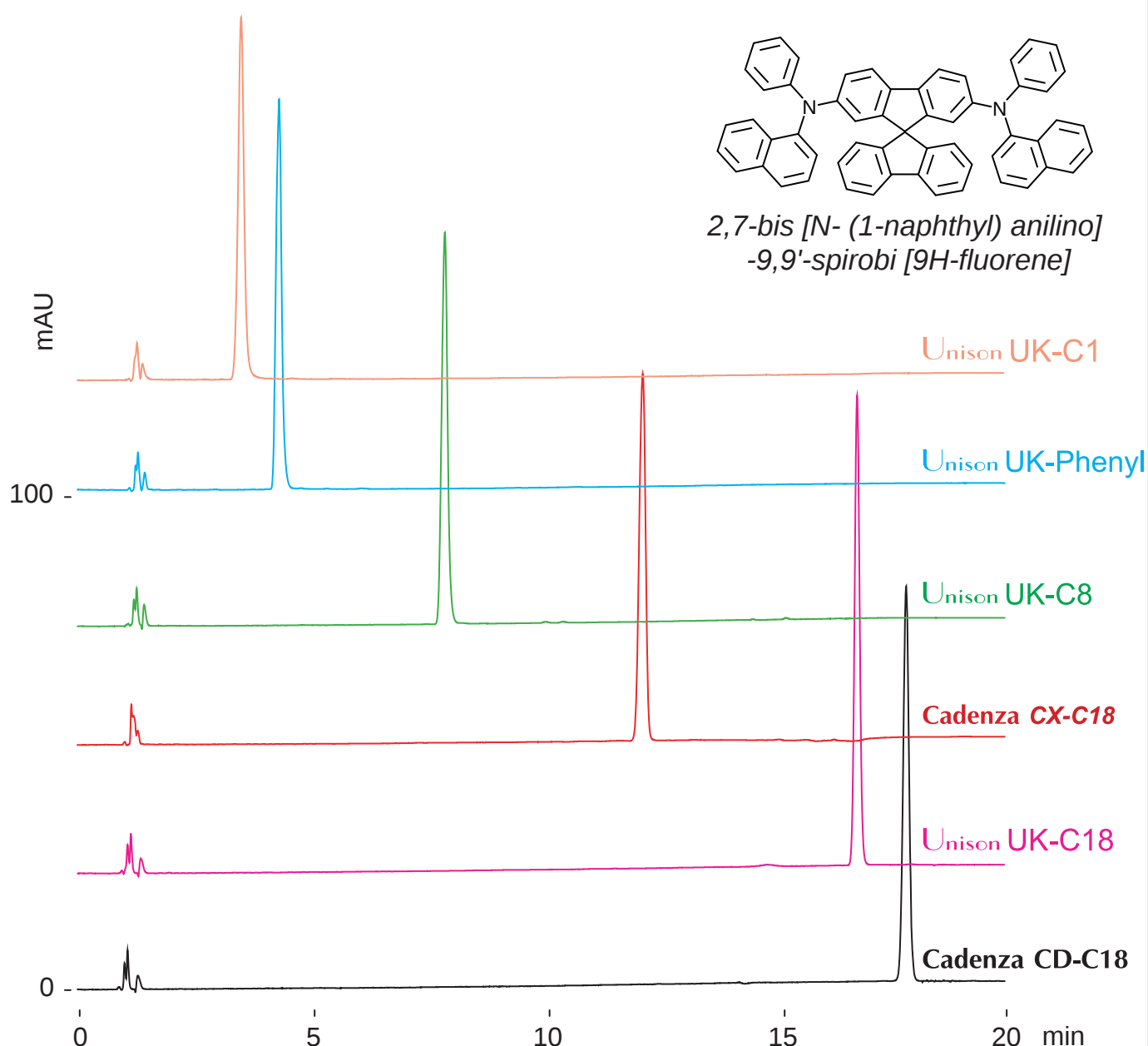

100 x 4.6 mm

A: water / acetonitrile / formic acid = 50 / 50 / 0.1

B: acetonitrile / formic acid = 100 / 0.1

60-100%B (0-10min), 100%B (10-20min)

1.0mL/min (6MPa), 37°C, 260nm

10μL (0.05mg/mL, THF / ACN / HCOOH = 50 / 50 / 1)