

Nardis ND-RX


Nardis ND-NX

150 x 3 mm


HPLC COLUMN APPLICATION

LC-MS: trimethylamine, trimethylamine N-oxide
(Choline metabolites)



トリメチルアミン, トリメチルアミン N-オキシド (コリン代謝物) の LC - MS 分析

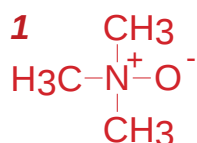

trimethylamine N-oxide
pKa 4.7

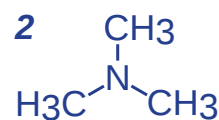

trimethylamine
pKa 9.8

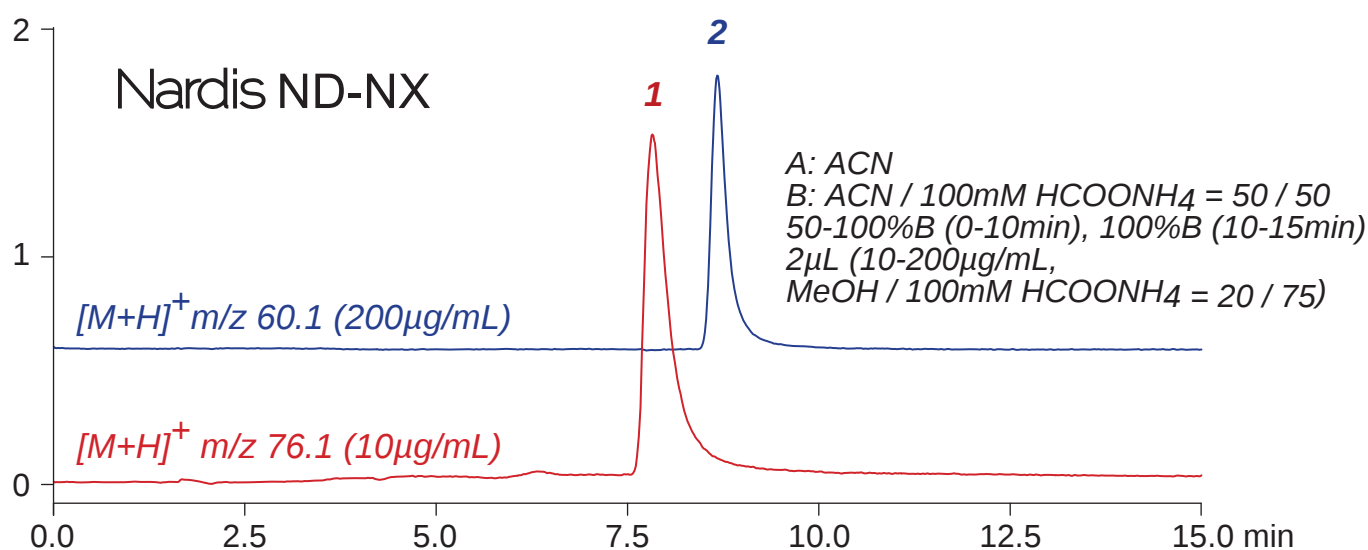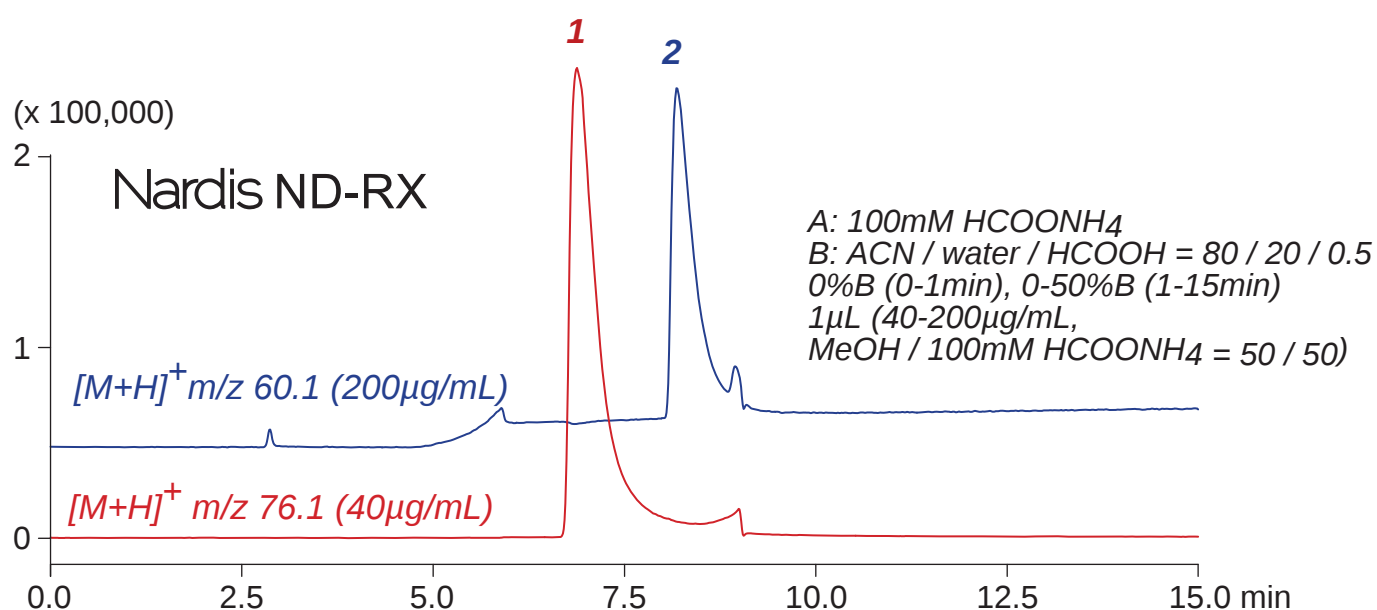

150 x 3 mm

0.4mL/min (6-9MPa), 37°C, Single Quad. MS (ESI)