

Nardis ND-NX




Unison UK-Amino

150 x 3 mm

HPLC COLUMN APPLICATION

LC-MS: trimethylamine, trimethylamine N-oxide
(Choline metabolites)

トリメチルアミン, トリメチルアミン N-オキシド(コリン代謝物)の LC-MS 分析

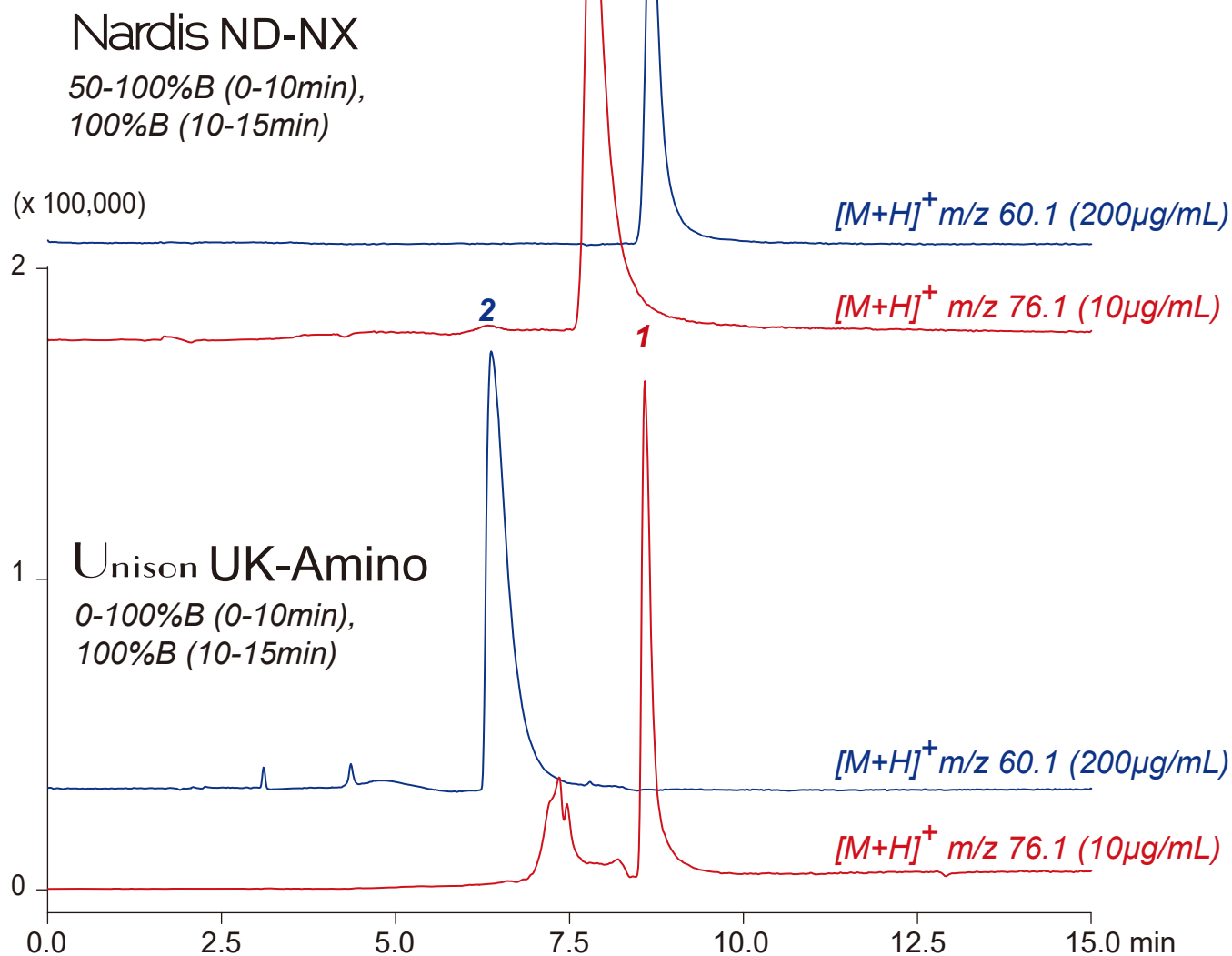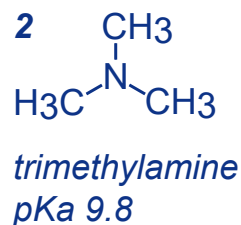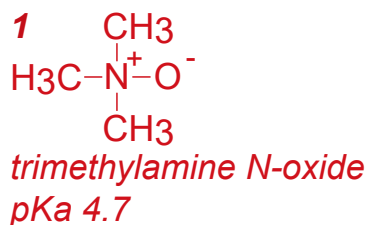

150 x 3 mm

A: acetonitrile

B: acetonitrile / 100mM ammonium formate = 50 / 50

0.4mL/min (5-6MPa), 37°C, Single Quad. MS (ESI)

2 μ L (10-200 μ g/mL, MeOH / 100mM HCOONH₄ = 25 / 75)