


Unison UK-Phenyl Unison UK-C1
Unison UK-C8 Cadenza CX-C18
Unison UK-C18 Cadenza CD-C18

100 x 4.6 mm


HPLC COLUMN APPLICATION

2,4,6-Triphenyl-1,3,5-triazine (Organic light-emitting diode material)

2,4,6-トリフェニル-1,3,5-トリアジン (有機EL材料)

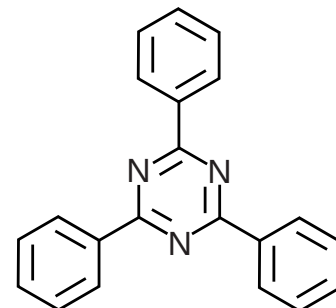

2,4,6-triphenyl-1,3,5-triazine

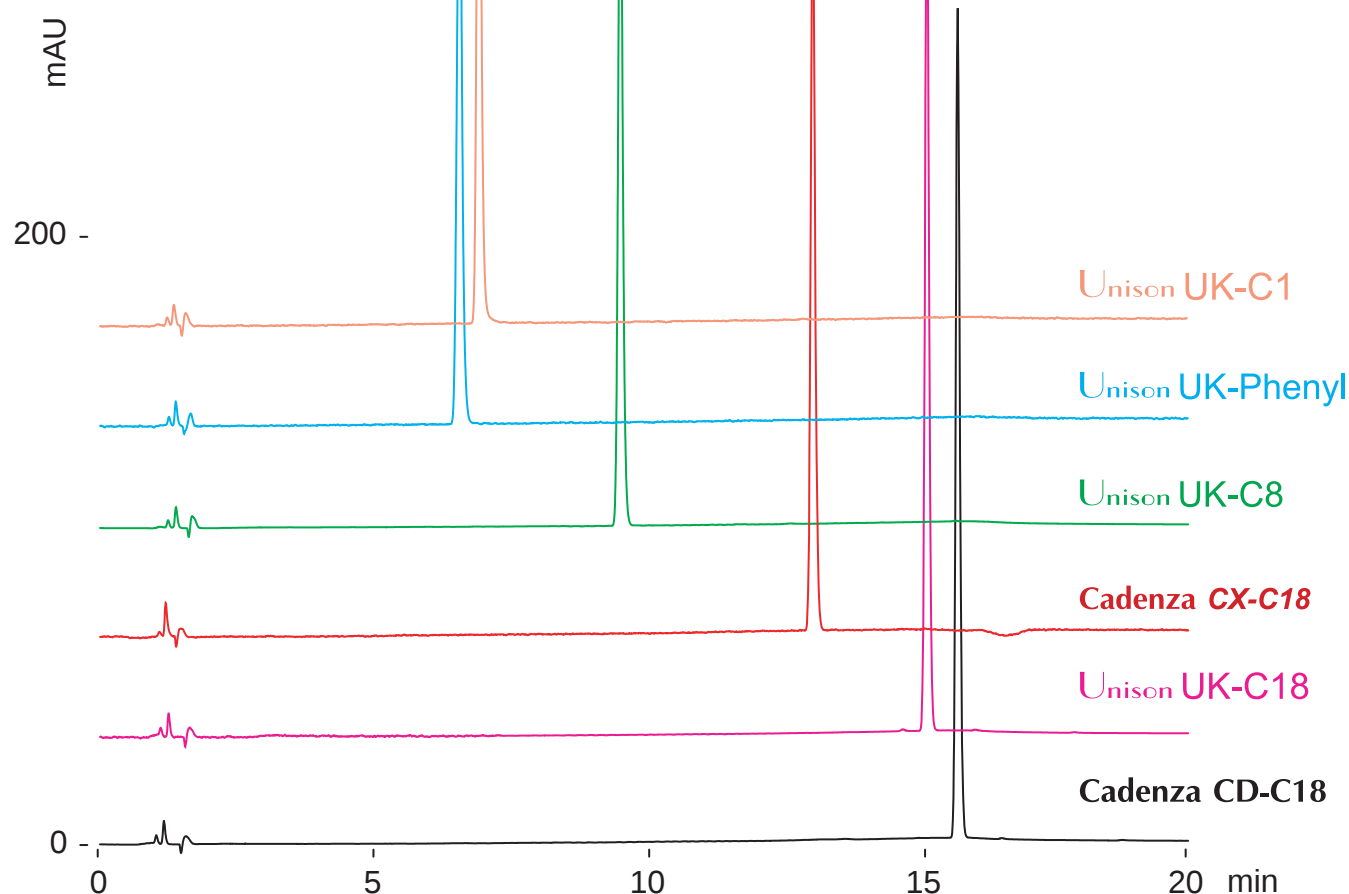

100 x 4.6 mm

A: water / acetonitrile / formic acid = 50 / 50 / 0.1

B: acetonitrile / formic acid = 100 / 0.1

0-100%B (0-15min), 100%B (15-20min)

1.0mL/min (9MPa), 37°C, 260nm

10μL (0.025mg/mL, THF / ACN / HCOOH = 50 / 50 / 1)