


Unison UK-Phenyl Unison UK-C1
Unison UK-C8 Cadenza CX-C18
Unison UK-C18 Cadenza CD-C18

100 x 4.6 mm


HPLC COLUMN APPLICATION

1,3-Bis [5- (4-tert-butylphenyl) -2- [1,3,4] oxadiazolyl] benzene (Organic light-emitting diode material)

1,3ビス [5-(4-tert-ブチルフェニル)2-[1,3,4]オキサジアゾリル]ベンゼン
(有機EL材料)

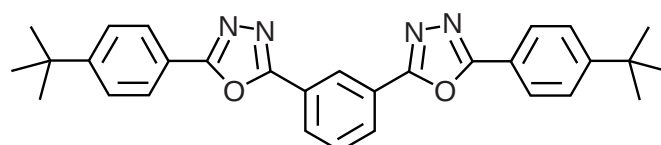

1,3-bis [5- (4-tert-butylphenyl) -2-
[1,3,4] oxadiazolyl] benzene

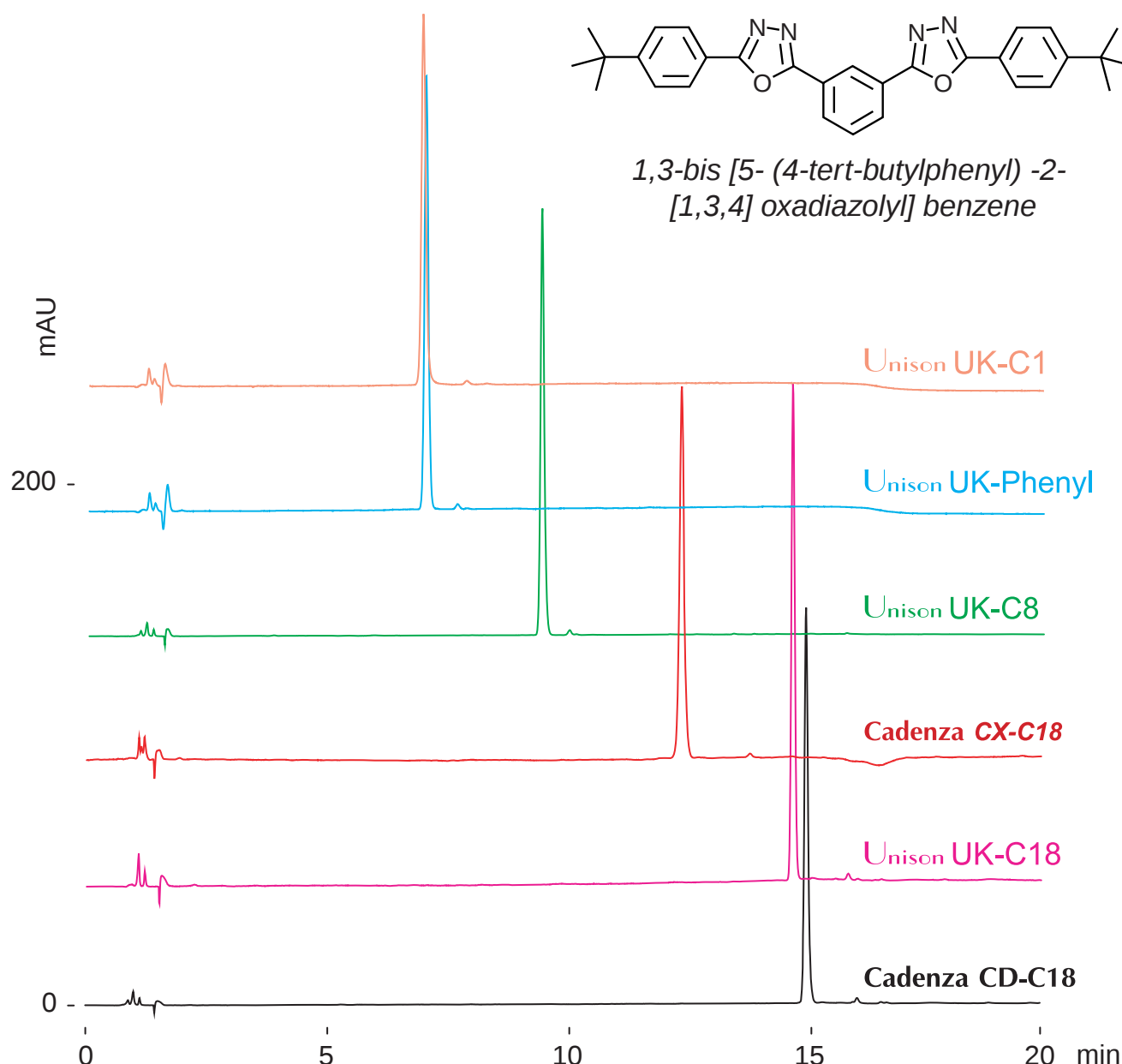

100 x 4.6 mm

A: water / acetonitrile / formic acid = 50 / 50 / 0.1

B: acetonitrile / formic acid = 100 / 0.1

0-100%B (0-15min), 100%B (15-20min)

1.0mL/min (9MPa), 37°C, 290nm

10μL (0.025mg/mL, THF / ACN / HCOOH = 50 / 50 / 1)