
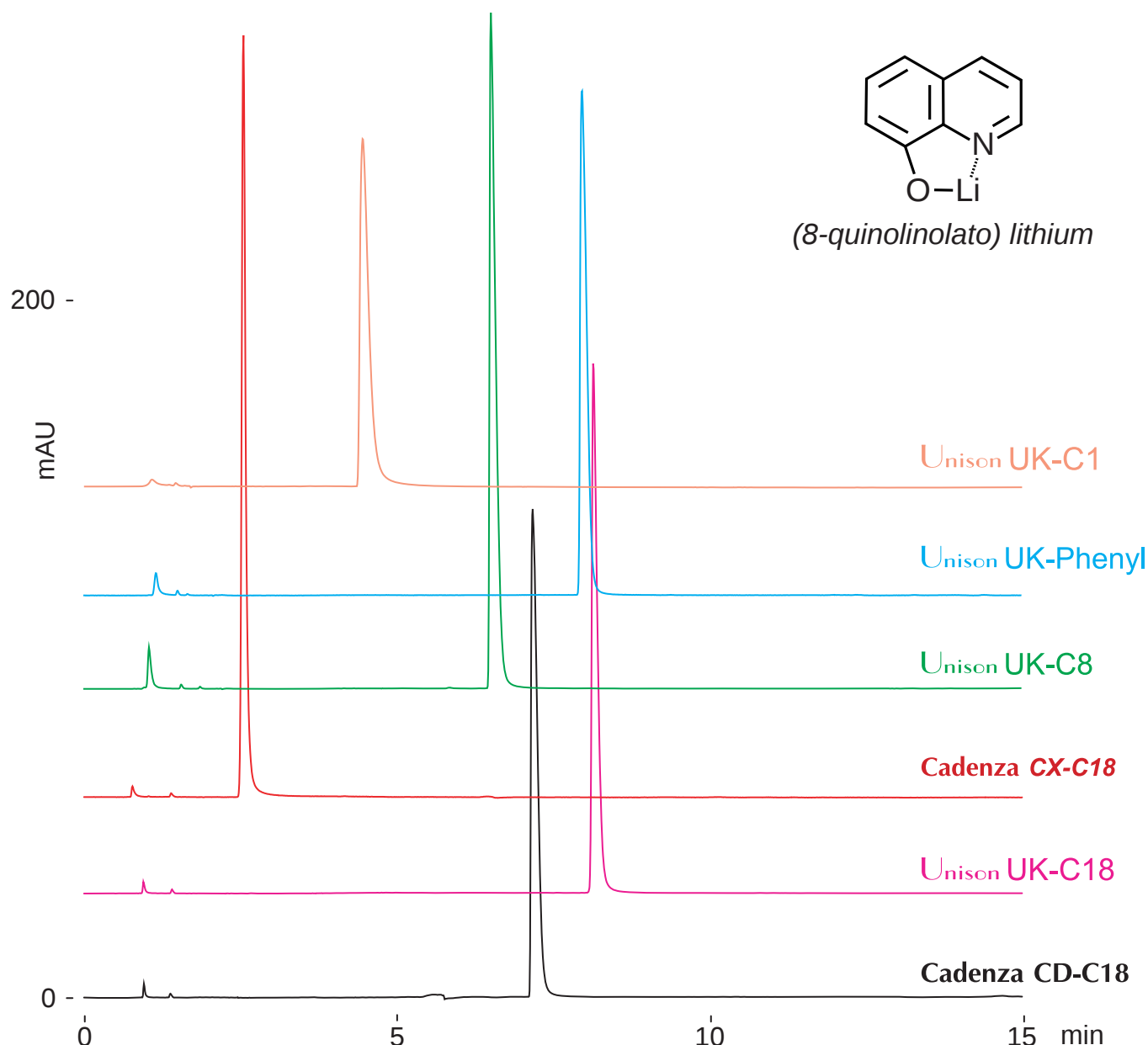


Unison UK-Phenyl Unison UK-C1
Unison UK-C8 Cadenza CX-C18
Unison UK-C18 Cadenza CD-C18

100 x 4.6 mm

HPLC COLUMN APPLICATION

(8-Quinolinolato) lithium
(Organic light-emitting diode material)
(8キノリノラト) リチウム
(有機 E L 材料)

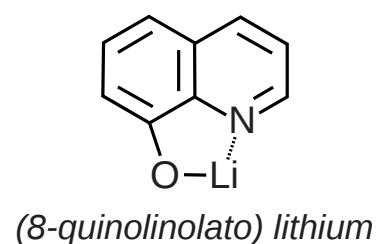

100 x 4.6 mm

A: water / formic acid = 100 / 0.5

B: water / acetonitrile / formic acid = 50 / 50 / 0.5

0%B (0-1min), 0-25%B (1-15min)

1.0mL/min (9MPa), 37°C, 260nm

10μL (0.05mg/mL, water / formic acid = 100 / 1)