


Unison UK-C8

100 x 4.6 mm

HPLC COLUMN APPLICATION

Organic light-emitting electron transport materials

有機エレクトロルミネッセンス電子輸送層材料

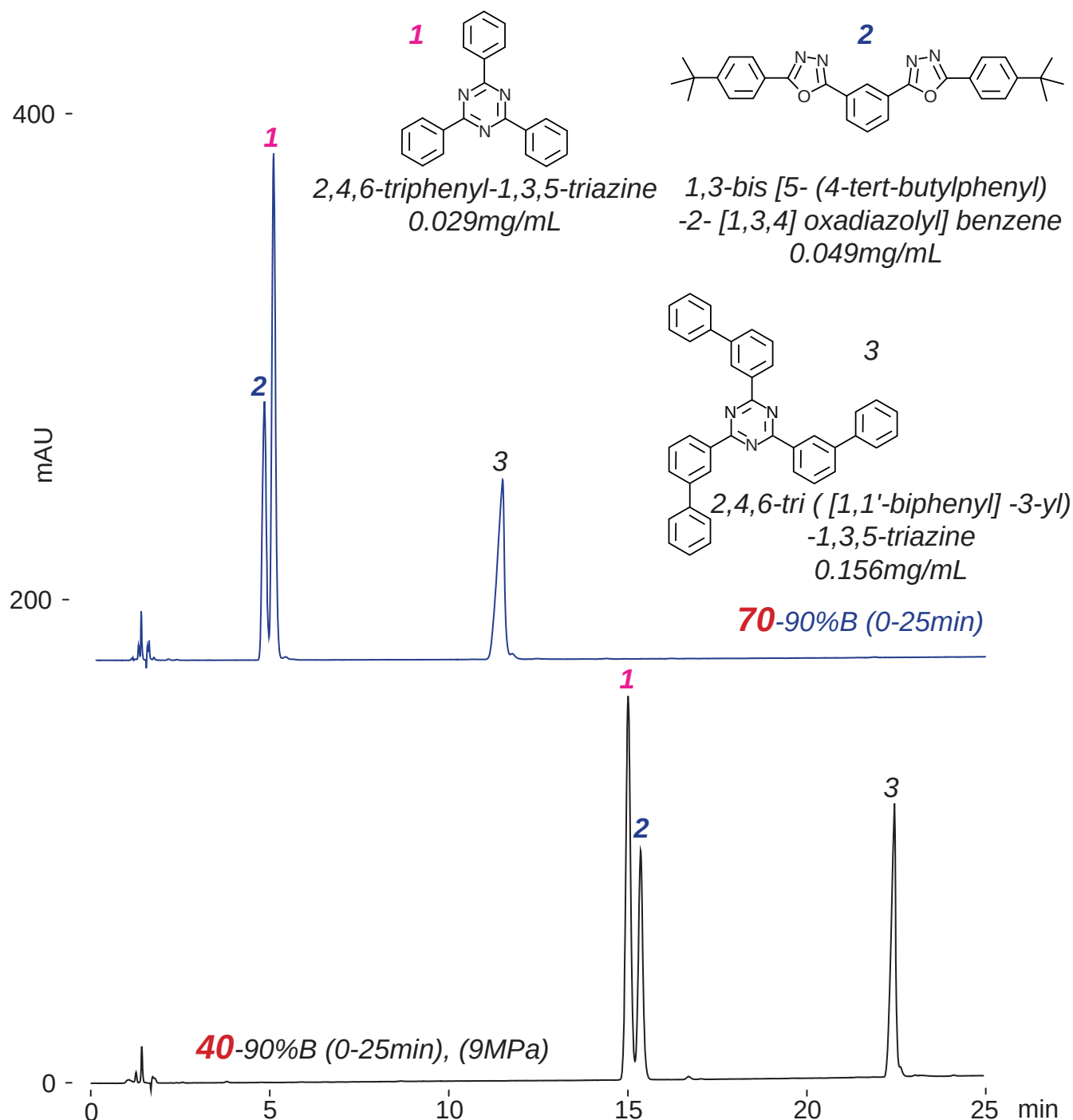

Unison UK-C8, 100 x 4.6 mm

A: water / acetonitrile / formic acid = 50 / 50 / 0.1

B: acetonitrile / formic acid = 100 / 0.1

1.0mL/min (6-9MPa), 37°C, 260nm

10μL (0.029-0.156mg/mL, THF / ACN / HCOOH = 50 / 50 / 1)