

Unison UK-C8

Cadenza CD-C18


Unison UK-C18

150 x 4.6 mm


HPLC COLUMN APPLICATION

PSI-6206 (HCV(hepatitis C virus) protease inhibitor)

PSI-6206 (HCV (C型肝炎ウイルス) プロテアーゼ阻害剤)

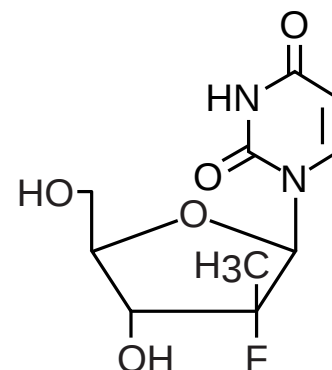

PSI-6206
((2'R)-2'-deoxy-2'-fluoro-2'-methyluridine)

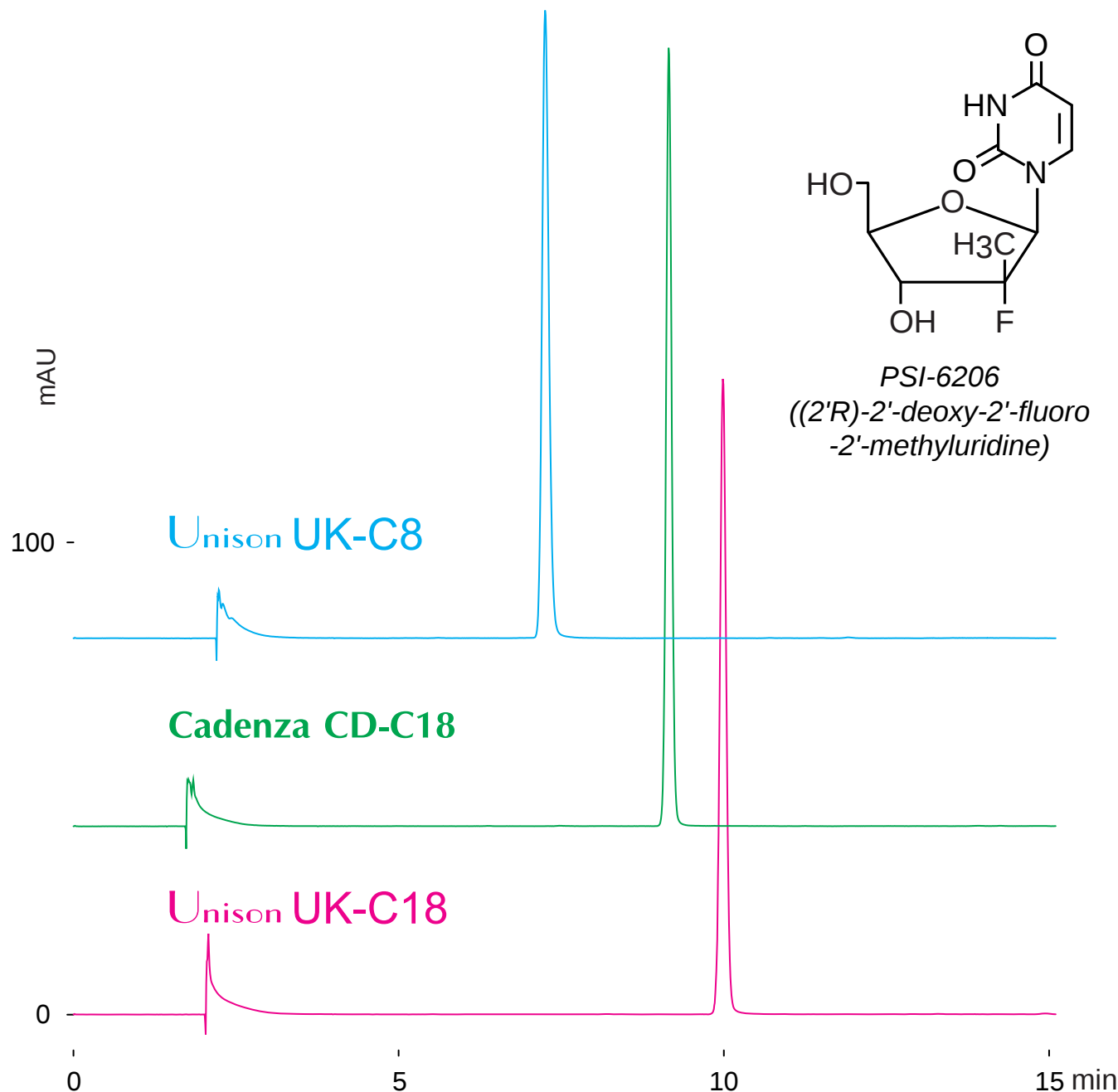

150 x 4.6 mm

A: 100mM ammonium formate

B: 100mM ammonium formate / acetonitrile = 50 / 50

0-15%B (0-15min)

1.0mL/min (9-10MPa), 37°C, 260nm

10μL (0.1mg/mL, dimethyl sulfoxide)