



Unison UK-C8
Unison UK-Phenyl
Unison UK-C18
Cadenza CD-C18

150 x 4.6 mm

HPLC COLUMN APPLICATION

Entecavir (Antiviral)

エンテカビル (抗ウイルス薬)

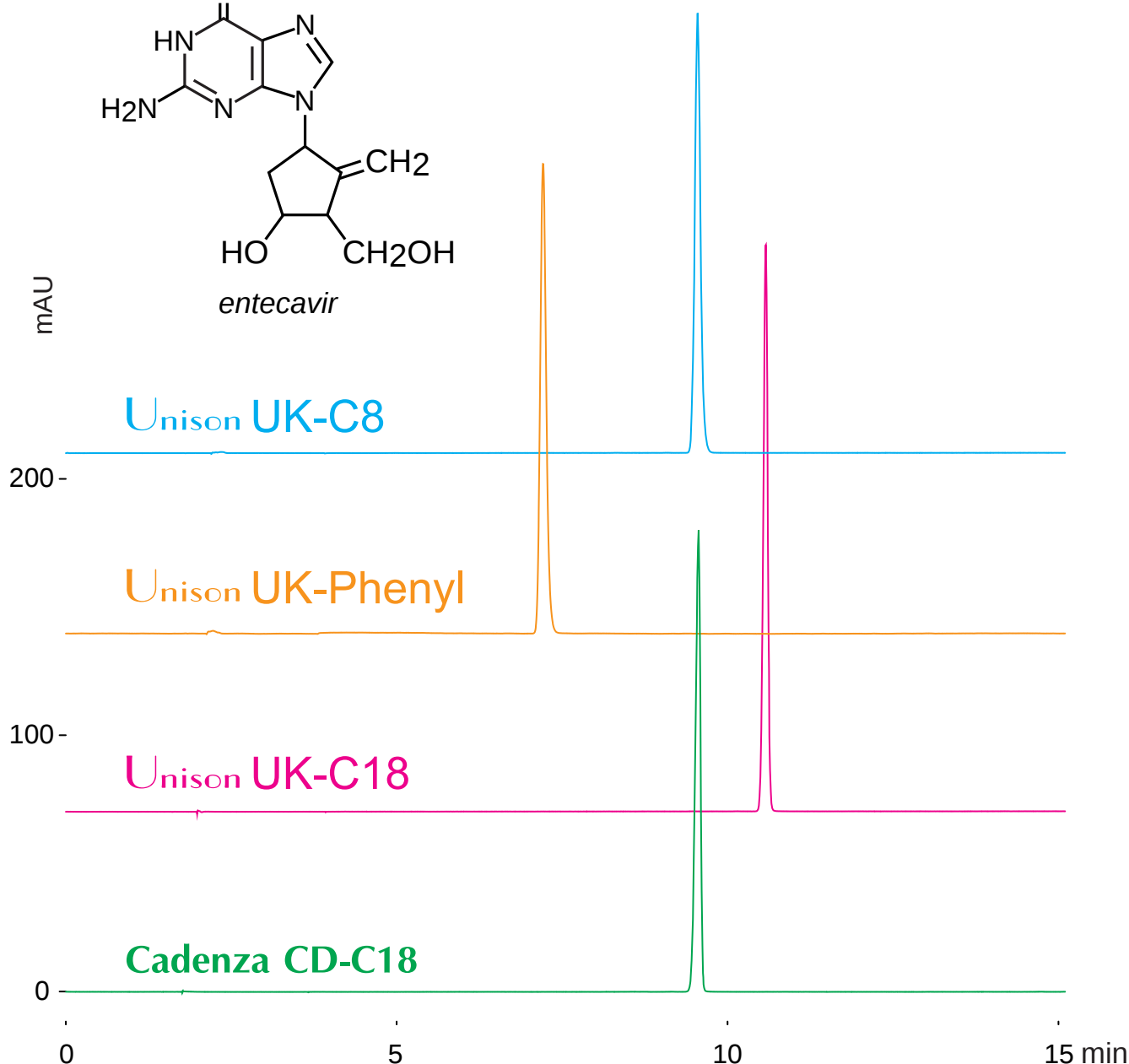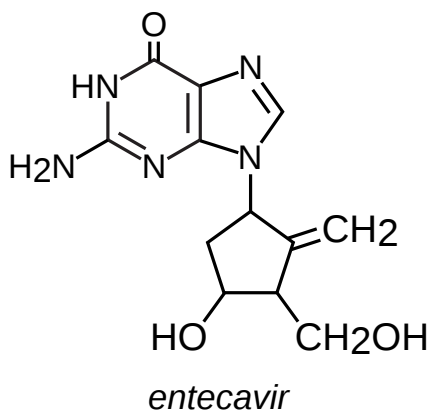

150 x 4.6 mm

A: 100mM ammonium formate

B: 100mM ammonium formate / acetonitrile = 50 / 50

0-25%B (0-15min)

1.0mL/min (10-11MPa), 37°C, 260nm

10μL (0.1mg/mL, methanol/100mM ammonium formate=50/50)