



Cadenza CL-C18  
Cadenza CD-C18  
Unison UK-C18

250 x 4.6 mm


HPLC COLUMN APPLICATION

n-octyl-beta-D-thioglucopyranoside,  
n-octyl beta-D-glucopyranoside(Membrane protein solubilizer)

n-オクチル -Dチオグルコピラノシド、n-オクチル -Dグルコピラノシド  
(膜タンパク質可溶化剤)

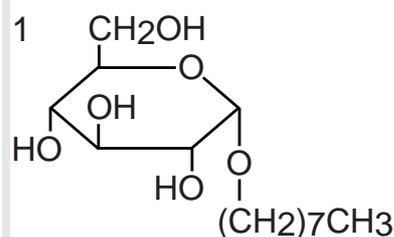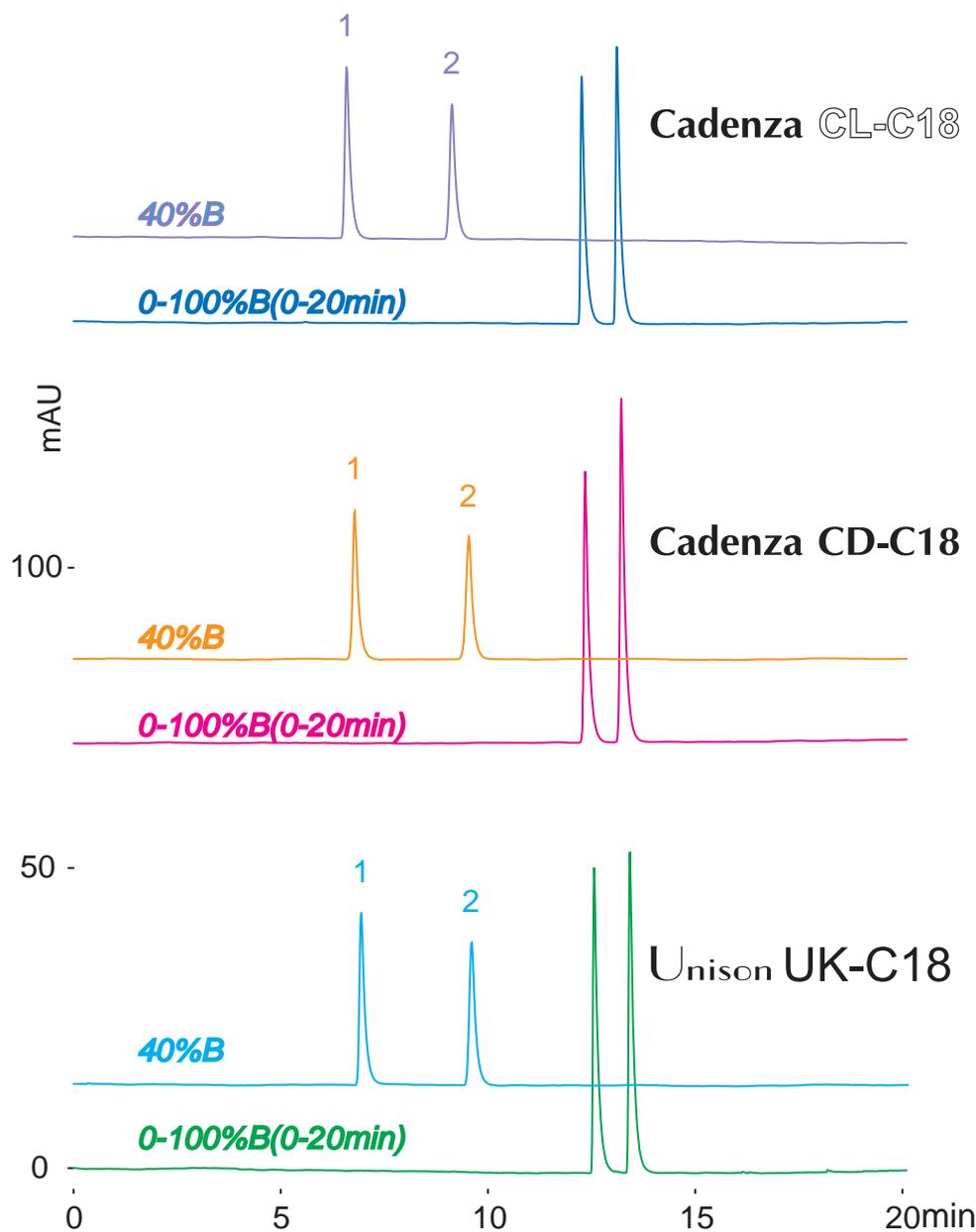

n-octyl beta-D-glucopyranoside

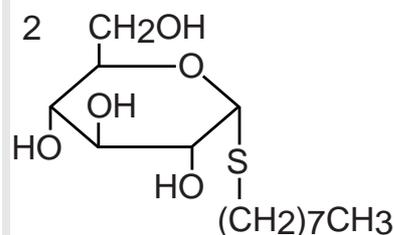

n-octyl beta-D-thioglucopyranoside

250 x 4.6 mm

A: water / acetic acid = 100 / 0.1

B: acetonitrile / acetic acid = 100 / 0.1

1.0mL/min (16-18MPa), 37°C

ELSD (spray chamber 70°C, drift chamber 100°C)

10µL (1mg/mL, water)