



Unison UK-C1
 Unison UK-Phenyl Cadenza CX-C18
 Unison UK-C8 Cadenza CW-C18
 Unison UK-C18 Cadenza CD-C18

100 x 4.6 mm

HPLC COLUMN APPLICATION

Ceftazidime (Cephem antibiotics)

セフトジジム (セフェム系抗生物質)

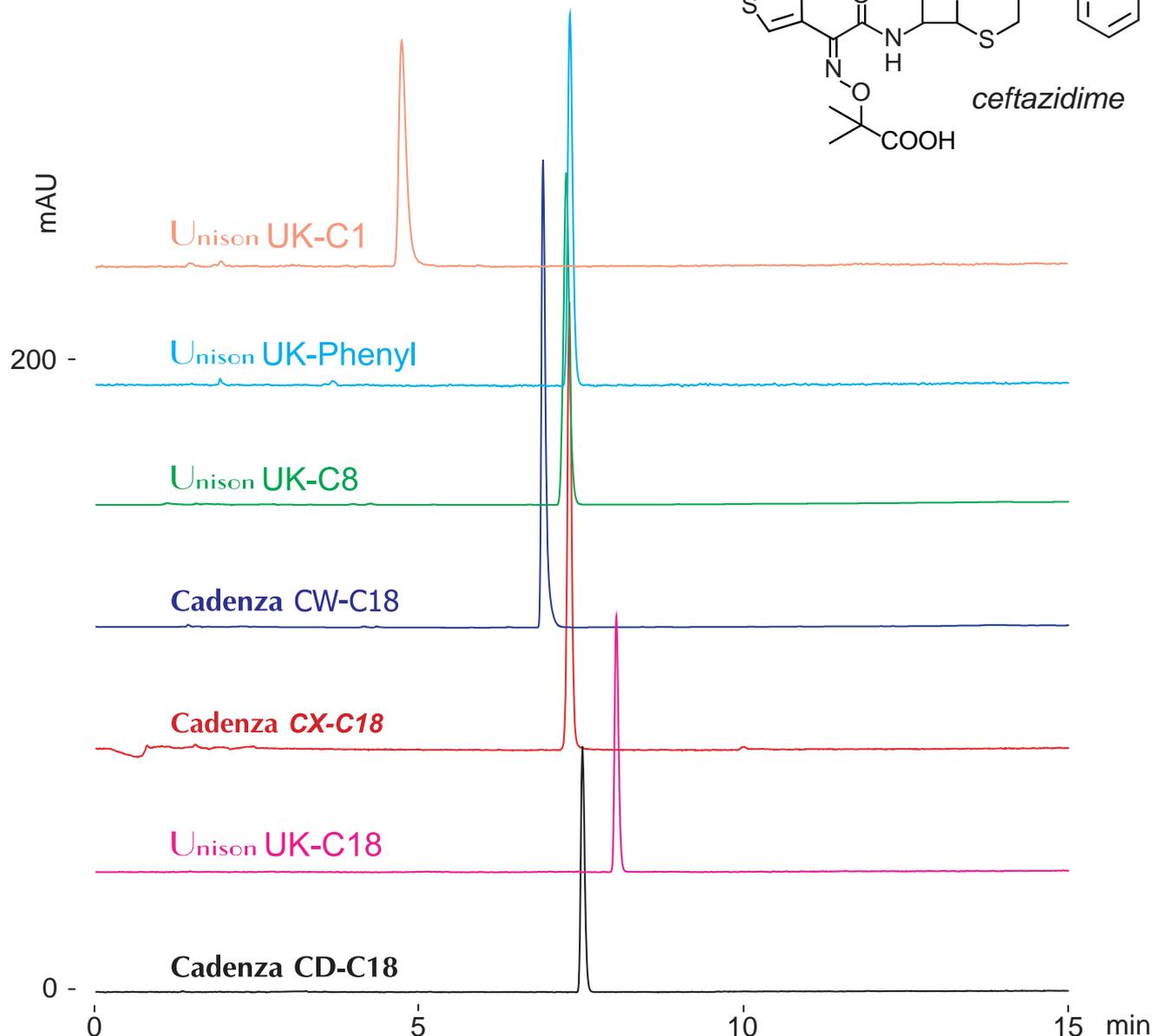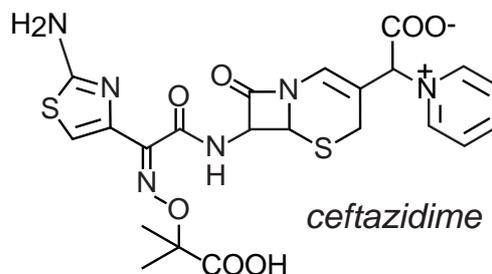

100 x 4.6 mm

A: water / formic acid = 100 / 0.1

B: water / acetonitrile / formic acid = 50 / 50 / 0.1

0-50%B (0-15min)

1.0mL/min (10MPa), 37°C, 260nm

10µL (0.05mg/mL, water / formic acid = 100 / 1)