




Scherzo SS-C18
 Scherzo SM-C18
 Scherzo SW-C18

100 x 4.6 mm

HPLC COLUMN APPLICATION

Dimesitylammonium pentafluorobenzenesulfonate
 (Condensation organocatalysts)

ジメシチルアンモニウム ペンタフルオロベンゼンスルホナート
 (縮合有機触媒)

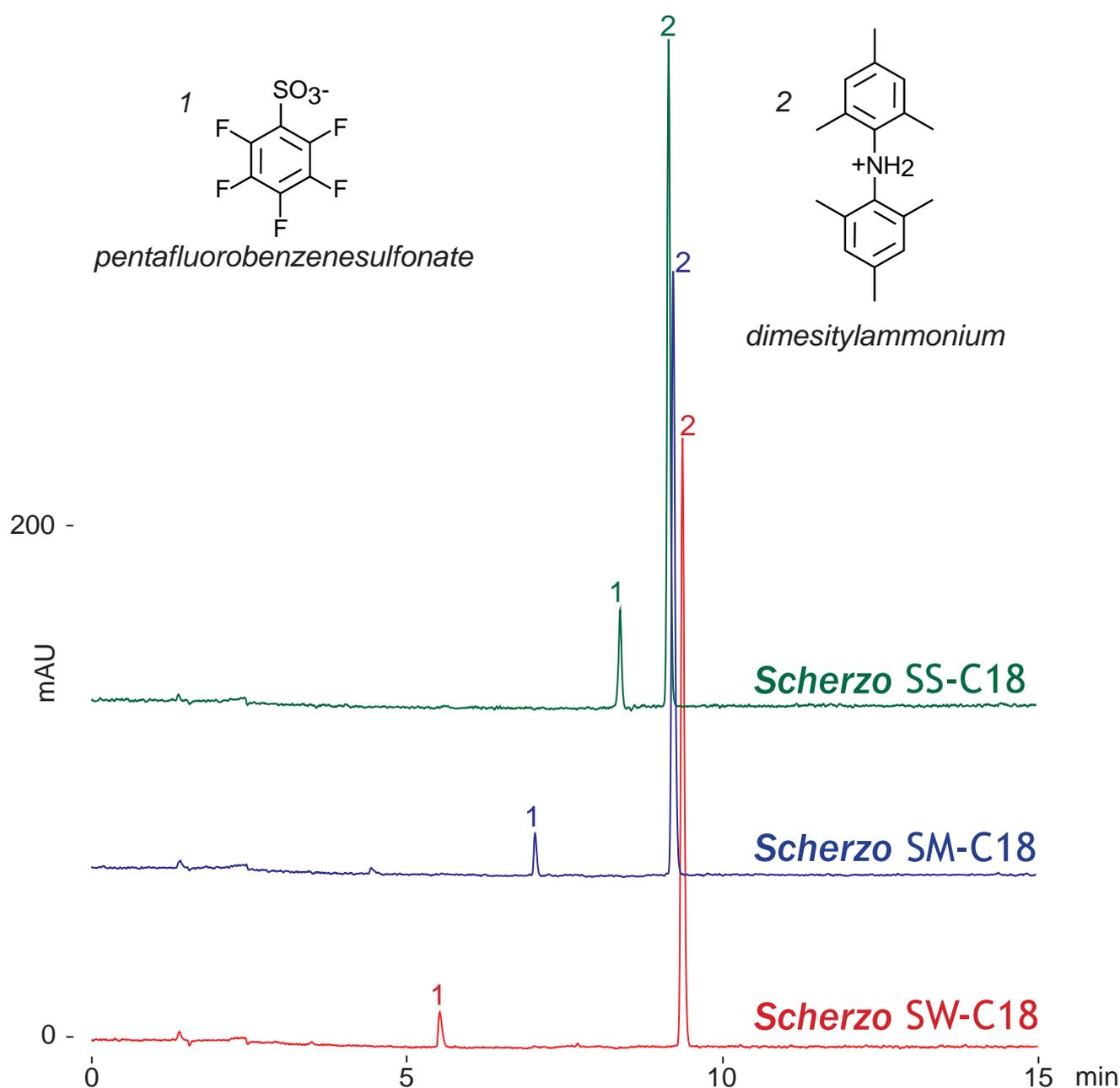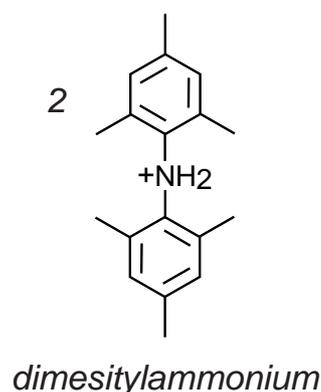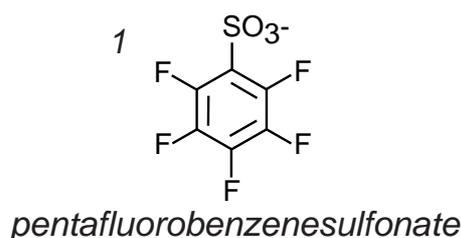

100 x 4.6 mm

A: water / formic acid = 100 / 0.1

B: 100mM ammonium formate / acetonitrile = 10 / 90

0-100%B (0-7min), 100%B (7-15min)

1.0mL/min (11.2MPa), 37°C, 260nm

10µL (0.1mg/mL, water / acetonitrile / formic acid= 25 / 75 / 1)