


Nardis ND-RX

100 x 4.6 mm

HPLC COLUMN APPLICATION

Dimesitylammonium pentafluorobenzenesulfonate (Condensation organocatalysts)

ジメシチルアンモニウム ペンタフルオロベンゼンスルホナート
 (縮合有機触媒)

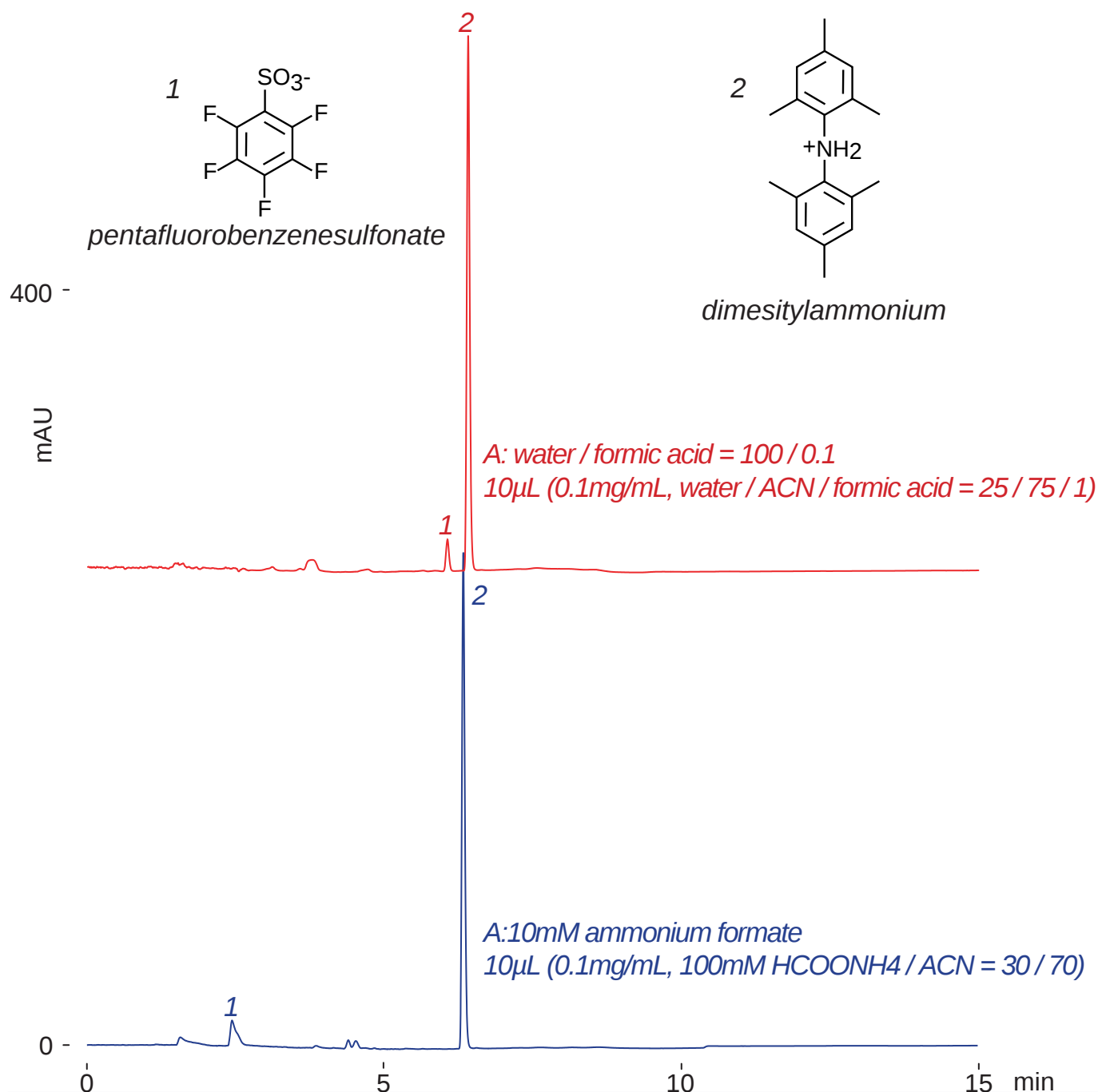

Nardis ND-RX, 100 x 4.6 mm

B: 100mM ammonium formate / acetonitrile = 10 / 90

0-100%B (0-7min), 100%B (7-15min)

1.0mL/min (8.7MPa), 37°C, 260nm