

Cadenza CW-C18


Presto FF-C18

250 x 4.6 mm

HPLC COLUMN APPLICATION

Sorbitan monopalmitate(Span 40)

ソルビタンモノパルミタート(スパン40)

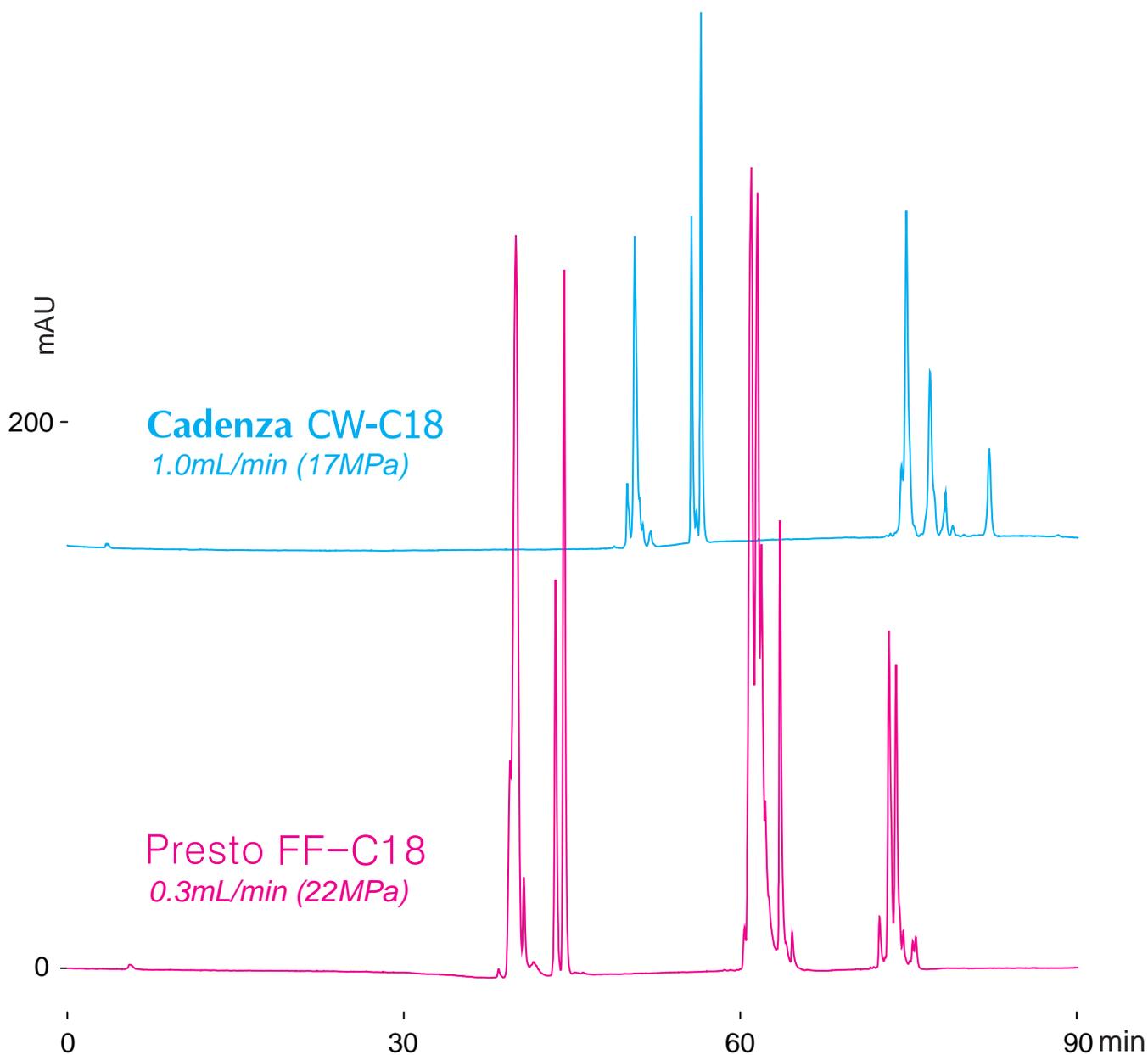

250 x 4.6 mm

A: water / acetic acid = 100 / 0.1

B: acetonitrile / acetic acid = 100 / 0.1

0-100%B (0-60min), 100%B (60-90min)

37°C

ELSD (spray chamber 70°C, drift chamber 100°C)

10 μ L (20mg/mL, ethanol)