


Unison UK-C18

150 x 2 mm

HPLC COLUMN APPLICATION

LC-MS: Veterinary drugs (Hormones)

動物用医薬品（ホルモン剤）のLC-MS分析

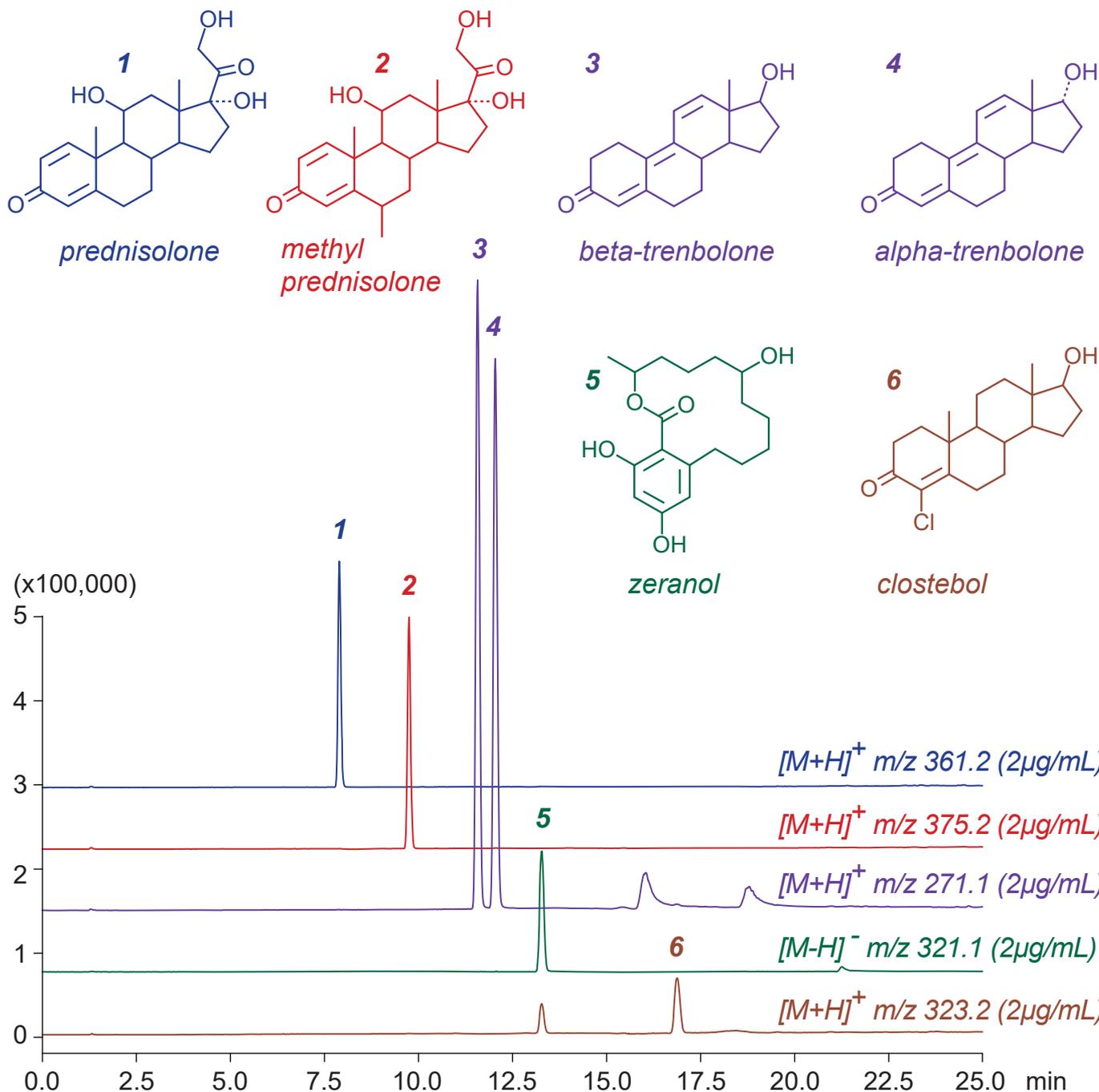

Unison UK-C18, 150 x 2 mm

A: acetonitrile / water / formic acid = 10 / 90 / 0.1

B: acetonitrile / water / formic acid = 90 / 10 / 0.1

15-100%B (0-25min)

0.2mL/min (10MPa), 37°C, Single Quad. MS (ESI)

2 μL (2 $\mu\text{g/mL}$, ACN / water / HCOOH = 50 / 50 / 1)