

Unison UK-Phenyl

Unison UK-C8

Unison UK-C18


Cadenza CD-C18

100 x 4.6 mm


HPLC COLUMN APPLICATION

Andrographolide(diterpene lactone)

アンドログラホリド (ジテルペンラクトン)

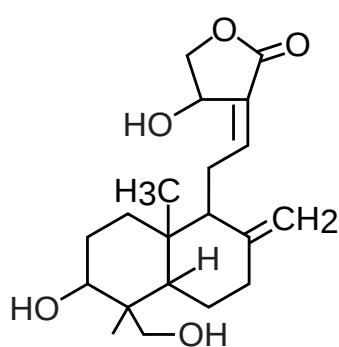

andrographolide

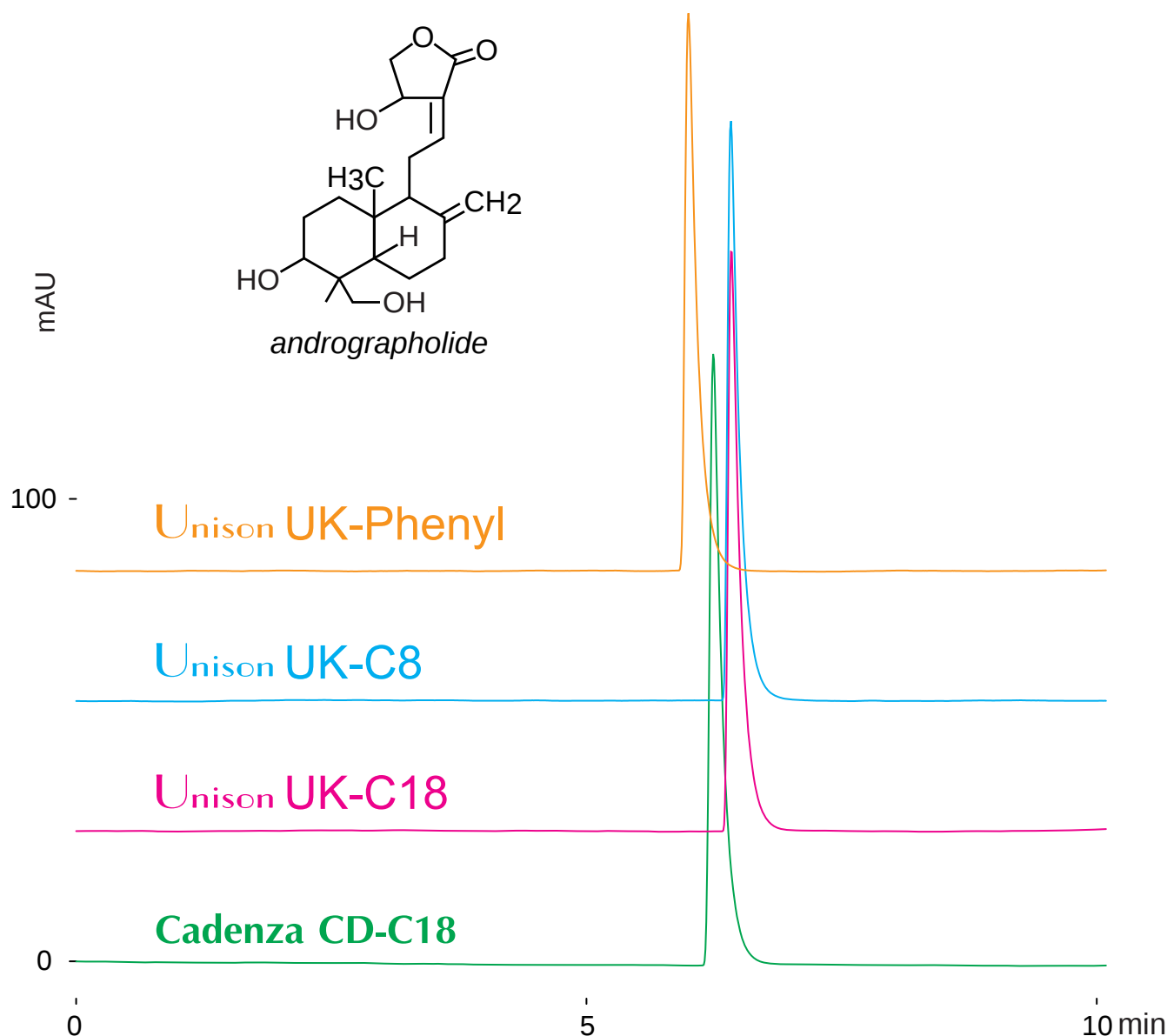

100 x 4.6 mm

A: water / formic acid = 100 / 0.1

B: acetonitrile / formic acid = 100 / 0.1

0-80%B (0-10min)

1.0mL/min (9-13MPa), 37°C

ELSD (spray chamber 70°C, drift chamber 100°C)

10μL (2mg/mL, dimethyl sulfoxide)