

Unison UK-C1

Unison UK-Phenyl

Unison UK-C8

Unison UK-C18

Cadenza CX-C18

Cadenza CW-C18


Cadenza CD-C18

100 x 4.6 mm


HPLC COLUMN APPLICATION

2- [(Diphenylphosphoryl) methyl] phenol (Condensation organocatalysts)

2-[(ジフェニルホスホリル) メチル] フェノール
(縮合有機触媒)

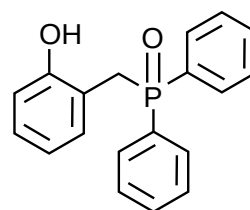

2- [(diphenylphosphoryl) methyl] phenol

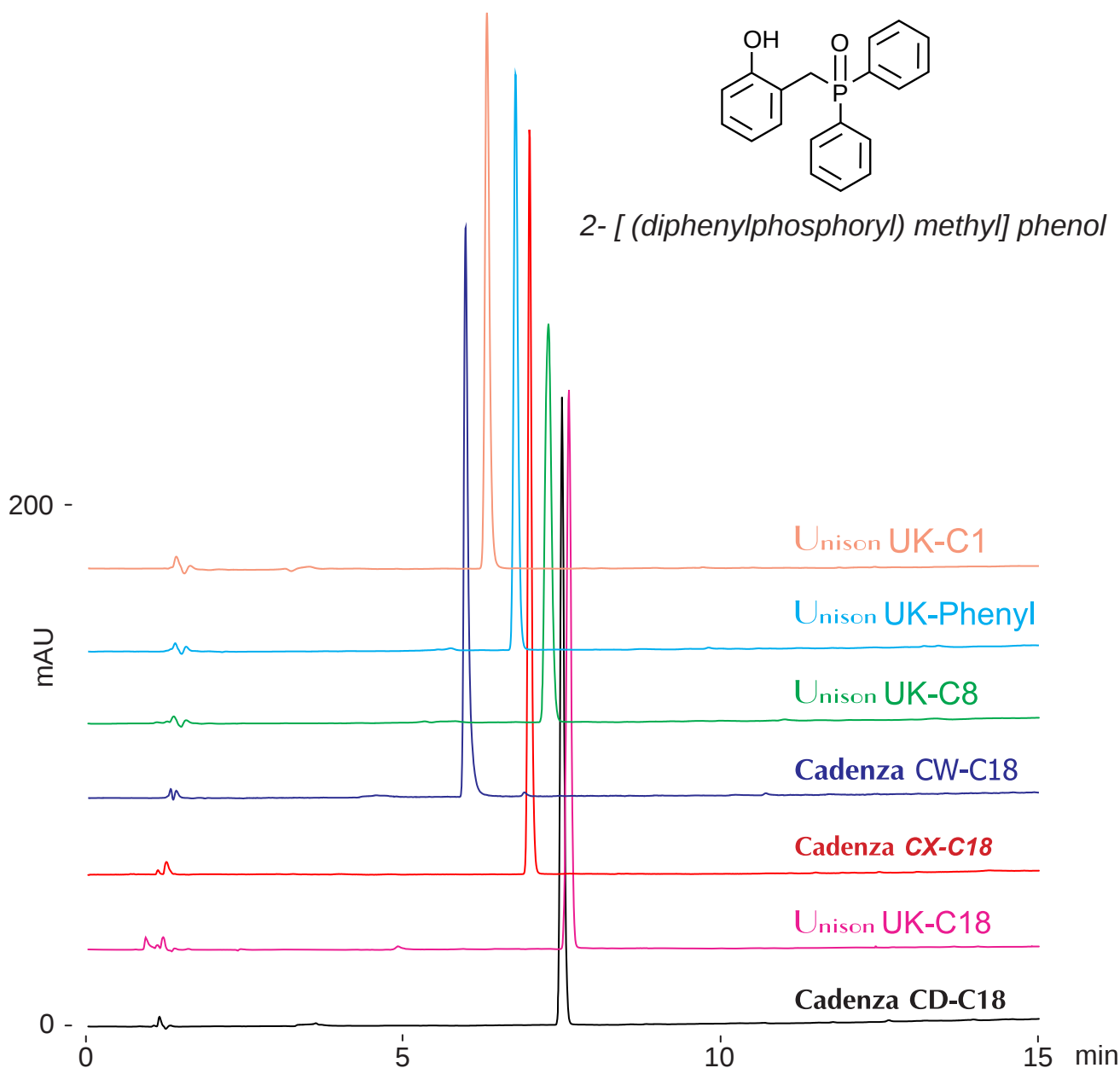

100 x 4.6 mm

A: water / acetonitrile / formic acid = 90 / 10 / 0.1

B: water / acetonitrile / formic acid = 10 / 90 / 0.1

15-100%B (0-15min)

1.0mL/min (10MPaq), 37°C, 260nm

10μL (0.5mg/mL, water / acetonitrile / formic acid = 50 / 50 / 1)