

Unison UK-C1
 Unison UK-Phenyl Cadenza CX-C18
 Unison UK-C8 Cadenza CW-C18
 Unison UK-C18 Cadenza CD-C18


100 x 4.6 mm

HPLC COLUMN APPLICATION

Pyridinium 3-nitrobenzenesulfonate

(Condensation organocatalysts)

ピリジニウム 3-ニトロベンゼンスルホナート
 (縮合有機触媒)

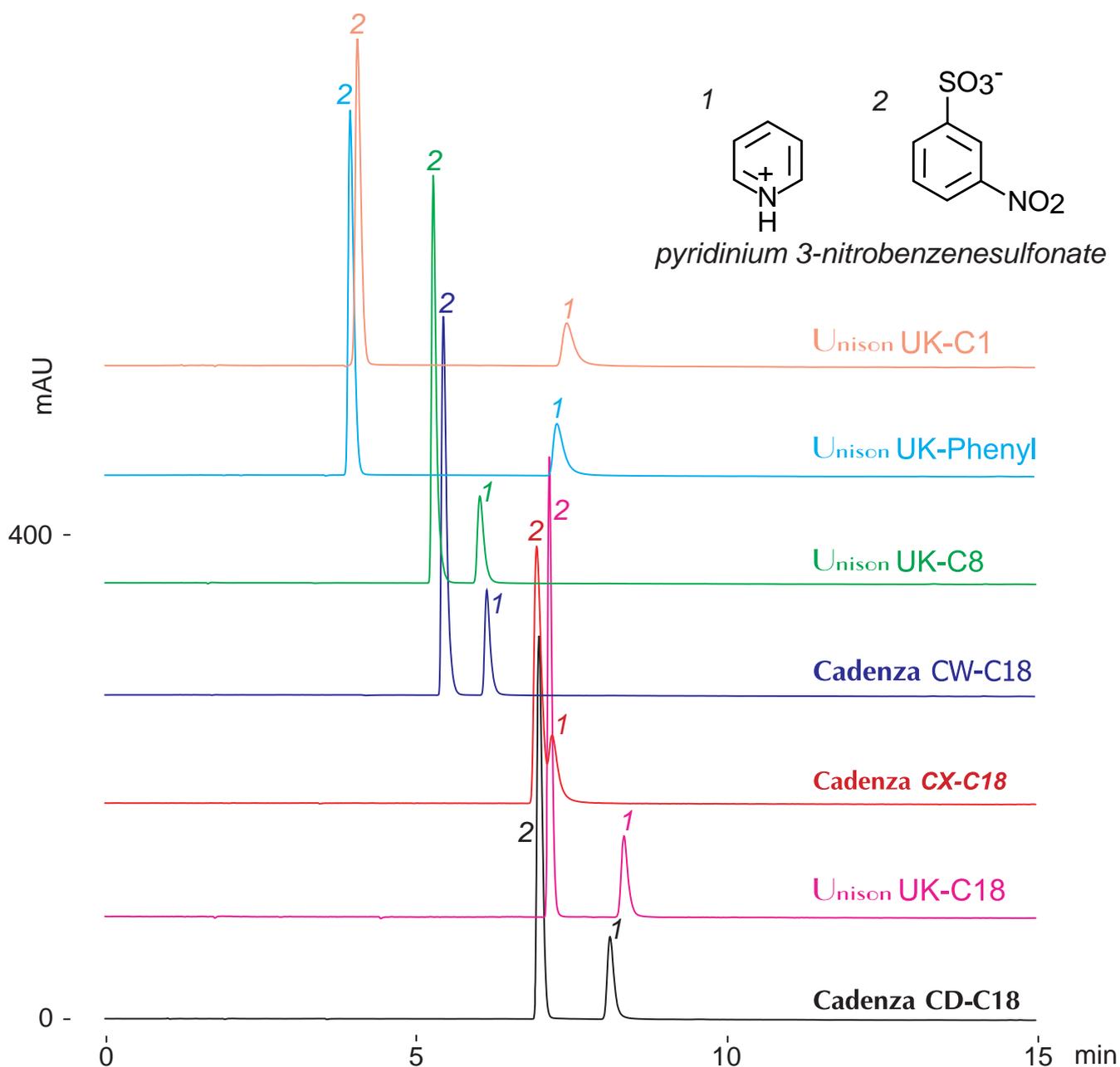

100 x 4.6 mm

A: 100mM ammonium formate

B: 100mM ammonium formate / acetonitrile = 50 / 50

0%B (0-1min), 0-50%B (1-15min)

1.0mL/min (12MPa), 37°C, 260nm

10µL (0.25mg/mL, 100mM ammonium formate)