

Scherzo SM-C18 Nardis ND-RX

Scherzo SW-C18 Scherzo SS-C18 100 x 4.6 mm


HPLC COLUMN APPLICATION

Pyridinium 3-nitrobenzenesulfonate

(Condensation organocatalysts)

ピリジニウム 3-ニトロベンゼンスルホナート

(縮合有機触媒)

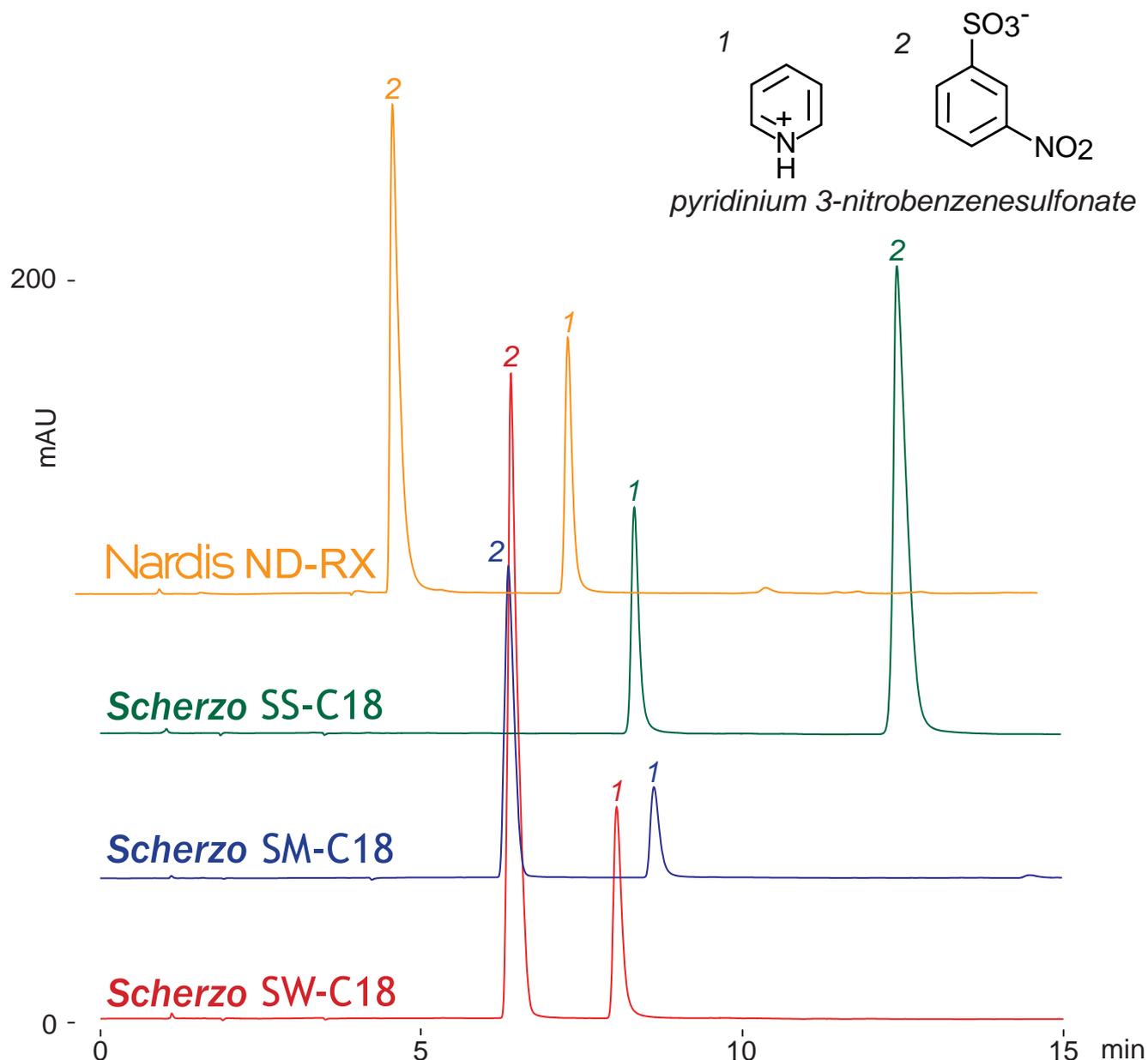

100 x 4.6 mm

A: 10mM ammonium formate

B: 100mM ammonium formate / acetonitrile = 50 / 50

0%B (0-1min), 0-50%B (1-15min)

1.0mL/min (12MPa), 37°C, 260nm

10µL (0.25mg/mL, 100mM ammonium formate)