

Scherzo SS-C18
 Scherzo SM-C18
 Scherzo SW-C18


100 x 4.6 mm

HPLC COLUMN APPLICATION

Pyridinium 3-nitrobenzenesulfonate

(Condensation organocatalysts)

ピリジニウム 3-ニトロベンゼンスルホナート
 (縮合有機触媒)

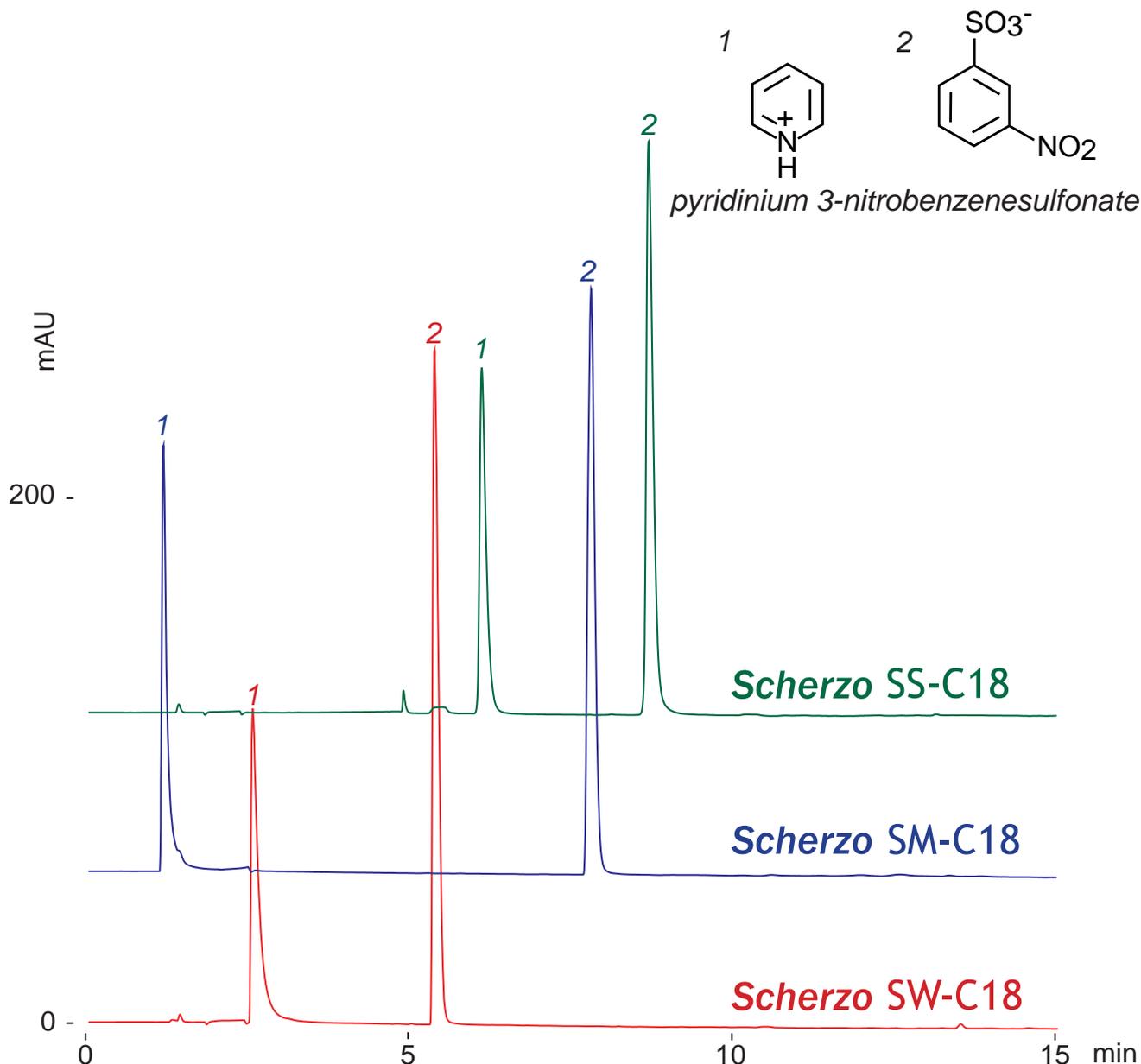

100 x 4.6 mm

A: water / formic acid = 100 / 0.1

B: 100mM ammonium formate / acetonitrile = 50 / 50

0-100%B (0-12min), 100%B(12-15min)

1.0mL/min (12MPa), 37°C, 260nm

10µL (0.25mg/mL, water / formic acid = 100 / 1)