


Scherzo SS-C18  
 Scherzo SM-C18  
 Scherzo SW-C18

100 x 4.6 mm


HPLC COLUMN APPLICATION

3,4,5-Trifluorophenylboronic acid  
 (Condensation organocatalysts)

3,4,5-トリフルオロフェニルボロン酸  
 (縮合有機触媒)

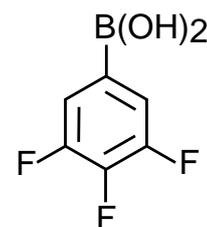

3,4,5-Trifluorophenyl  
boronic acid

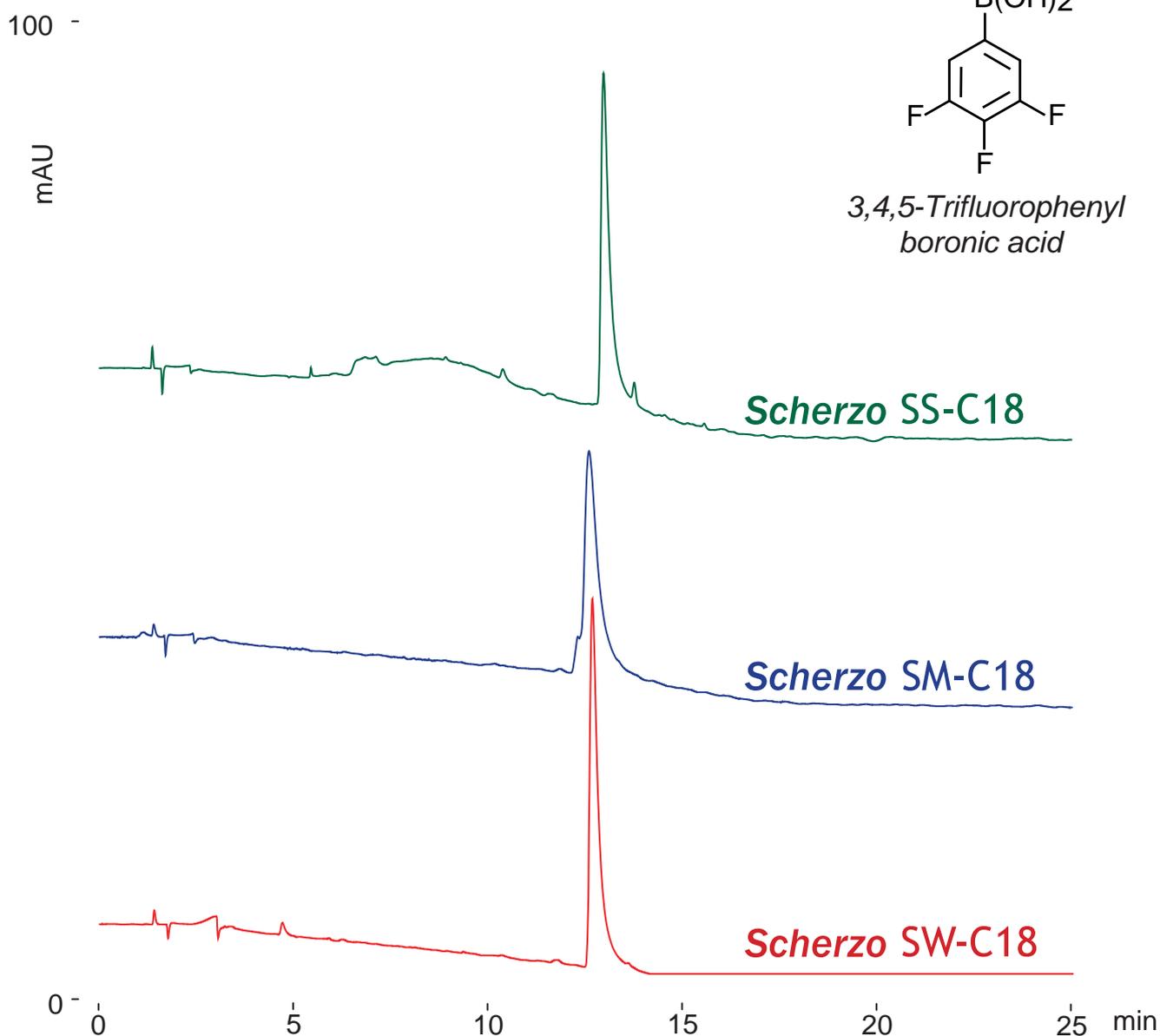

100 x 4.6 mm  
 A: water / formic acid = 100 / 0.5  
 B: 100mM ammonium formate / acetonitrile = 50 / 50  
 0-100%B (0-15min), 100%B (15-25min)  
 1.0mL/min (12MPa), 37°C, 260nm  
 10µL (1mg/mL, water / acetonitrile / formic acid = 80 / 20 / 1)