


Unison UK-C1
 Unison UK-Phenyl Cadenza CX-C18
 Unison UK-C8 Cadenza CW-C18
 Unison UK-C18 Cadenza CD-C18

100 x 4.6 mm


HPLC COLUMN APPLICATION

Crotamiton (Antipruritic)

クロタミトン (鎮痒薬)

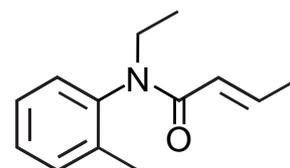

crotamiton

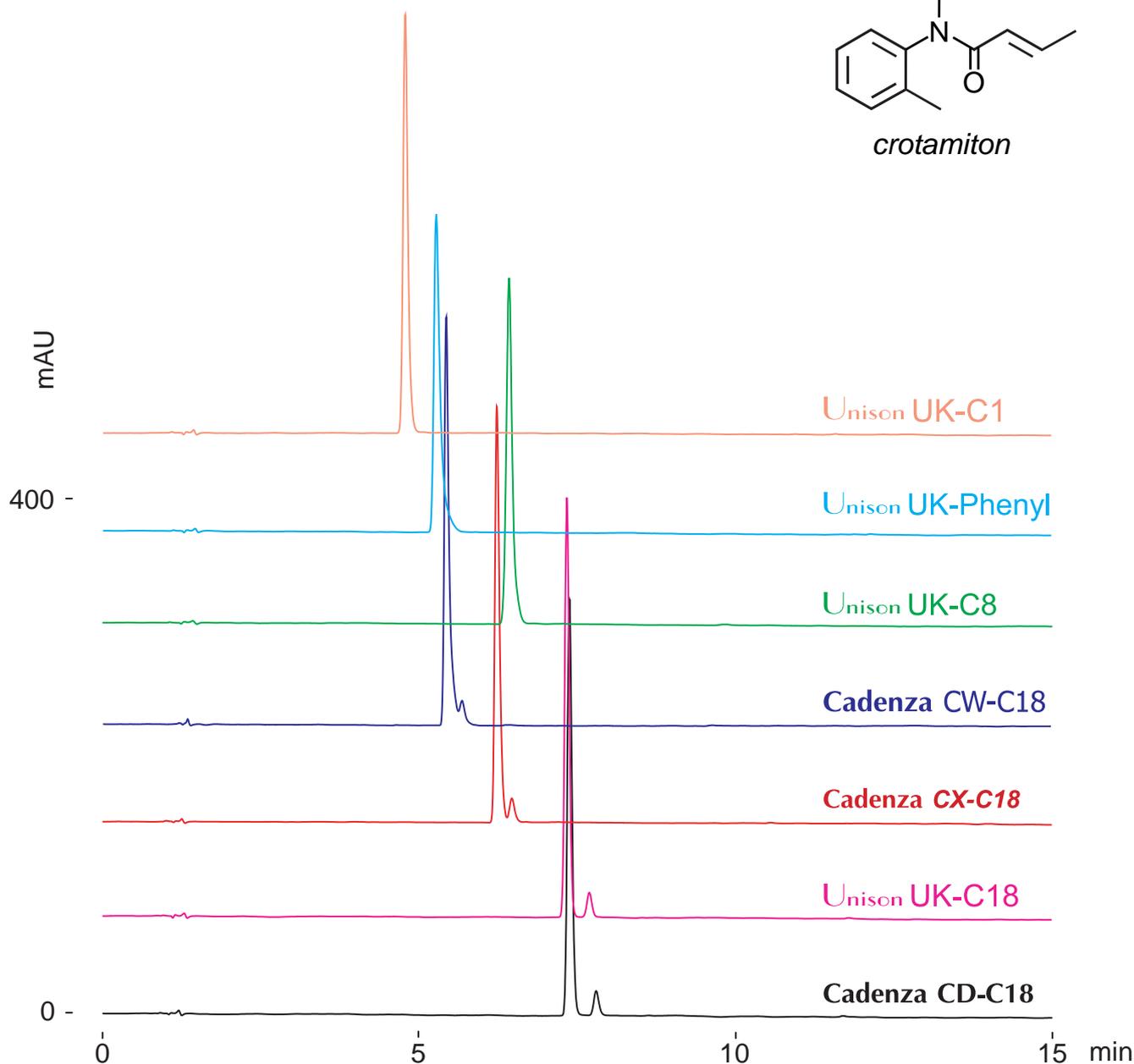

100 x 4.6 mm

A: 100mM ammonium formate / acetonitrile = 90 / 10

B: 100mM ammonium formate / acetonitrile = 10 / 90

25-100%B (0-15min)

1.0mL/min(13MPa), 37°C, 260nm

10µL (0.15mg/mL, 100mM ammonium formate / acetonitrile = 80 / 20)