



Unison UK-C1
 Unison UK-Phenyl Cadenza CX-C18
 Unison UK-C8 Cadenza CW-C18
 Unison UK-C18 Cadenza CD-C18

100 x 4.6 mm

HPLC COLUMN APPLICATION

Cromolyn (Antiallergic agent)

クロモグリク酸 (抗アレルギー薬)

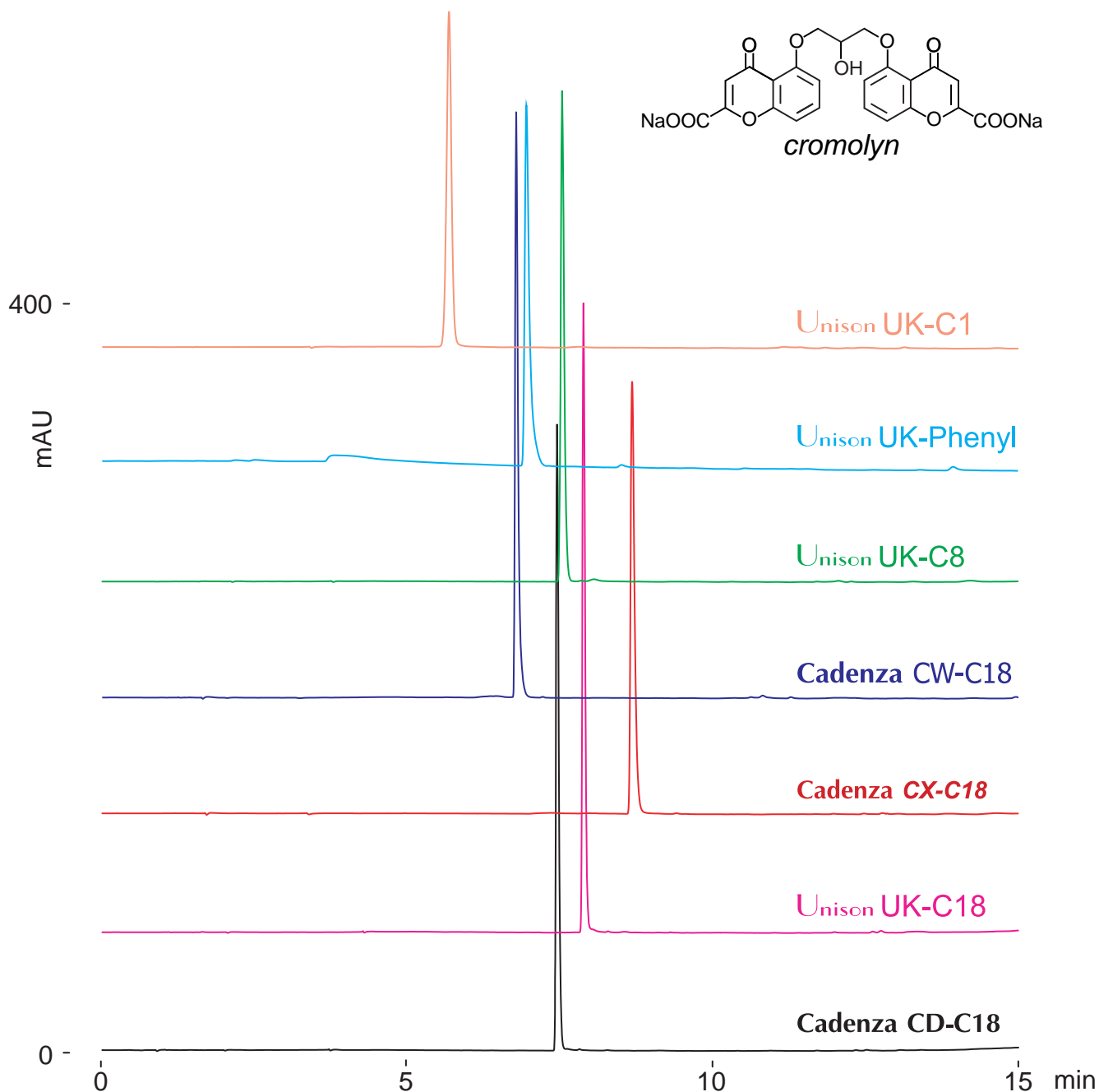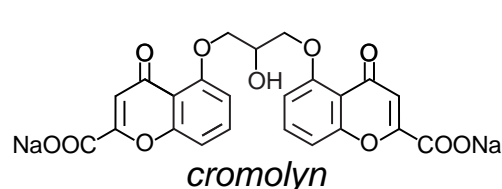

100 x 4.6 mm
 A: 100mM ammonium formate
 B: 100mM ammonium formate / acetonitrile = 50 / 50
 0%B (0-1min), 0-70%B (1-15min)
 1.0mL/min (12MPa), 37°C, 260nm
 10µL (0.1mg/mL, 100mM ammonium formate)