


Unison UK-C1
 Unison UK-Phenyl
 Unison UK-C8
 Unison UK-C18

Cadenza CX-C18
 Cadenza CW-C18
 Cadenza CD-C18


100 x 4.6 mm

HPLC COLUMN APPLICATION

Bupropion (Antidepressant) ブロピオン (抗うつ薬)

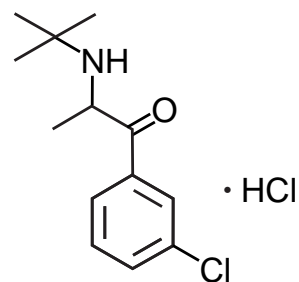

bupropion hydrochloride

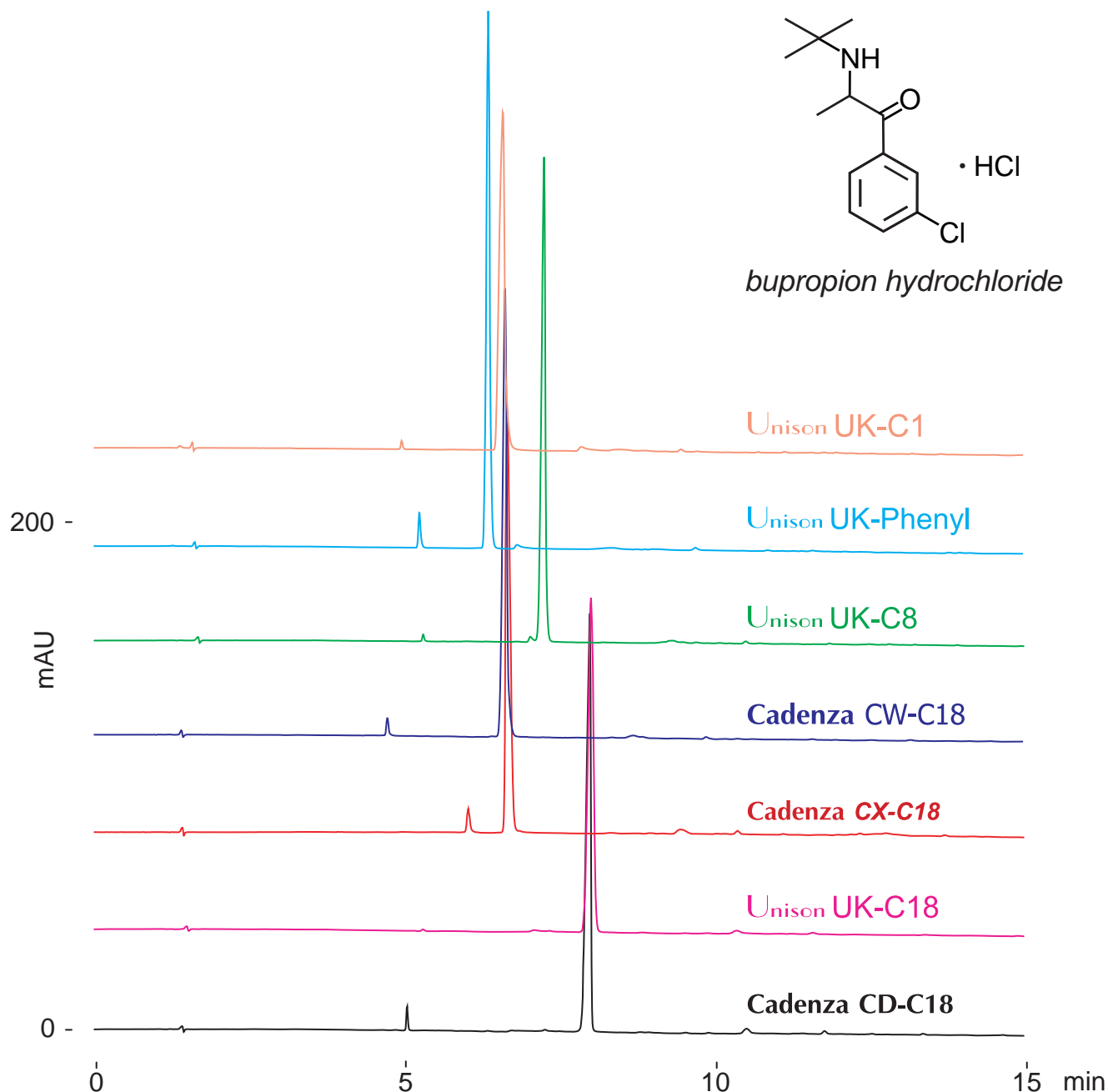

100 x 4.6 mm

A: 100mM ammonium formate / acetonitrile = 90 / 10

B: 100mM ammonium formate / acetonitrile = 10 / 90

0-100%B (0-15min)

1.0mL/min (10MPa), 37°C, 260nm

10µL (0.15mg/mL, 100mM ammonium formate)