



Unison UK-C1
 Unison UK-Phenyl Cadenza CX-C18
 Unison UK-C8 Cadenza CW-C18
 Unison UK-C18 Cadenza CD-C18

100 x 4.6 mm

HPLC COLUMN APPLICATION

Dobutamine (Cardiotonic)

ドブタミン (強心薬)

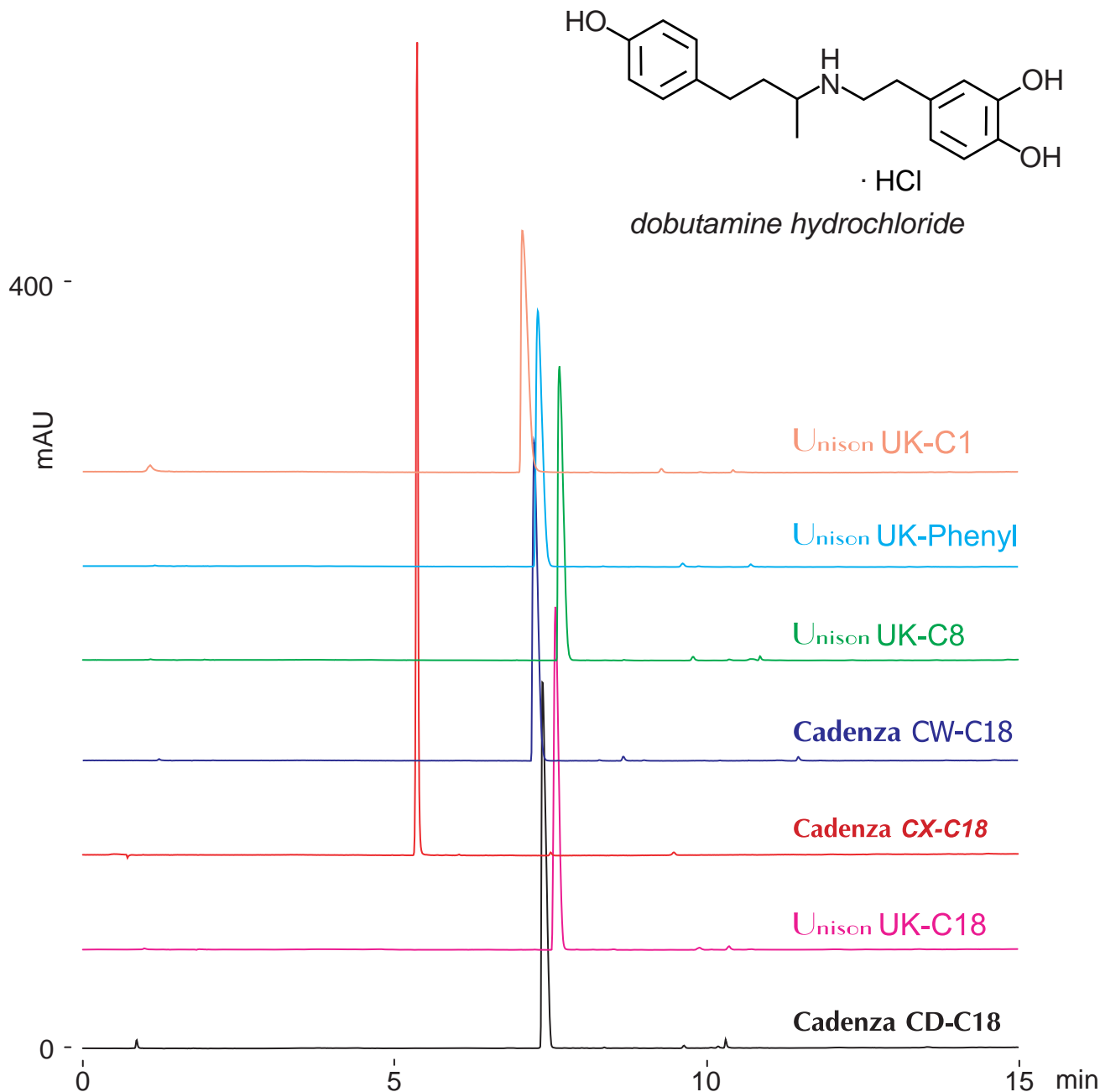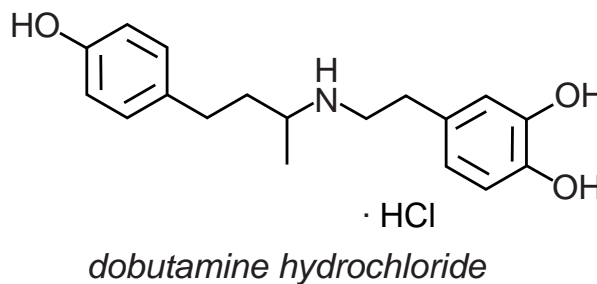

100 x 4.6 mm

A: water / formic acid = 100 / 0.5

B: water / acetonitrile / formic acid = 10 / 90 / 0.5

0%B (0-1min), 0-60%B (1-15min)

1.0mL/min (10MPa), 37°C, 280nm

10µL (0.33mg/mL, water / formic acid = 100 / 1)