


Unison UK-C1  
 Unison UK-Phenyl Cadenza CX-C18  
 Unison UK-C8 Cadenza CW-C18  
 Unison UK-C18 Cadenza CD-C18

100 x 4.6 mm

HPLC COLUMN APPLICATION

## Dobutamine (Cardiotonic)

ドブタミン (強心薬)

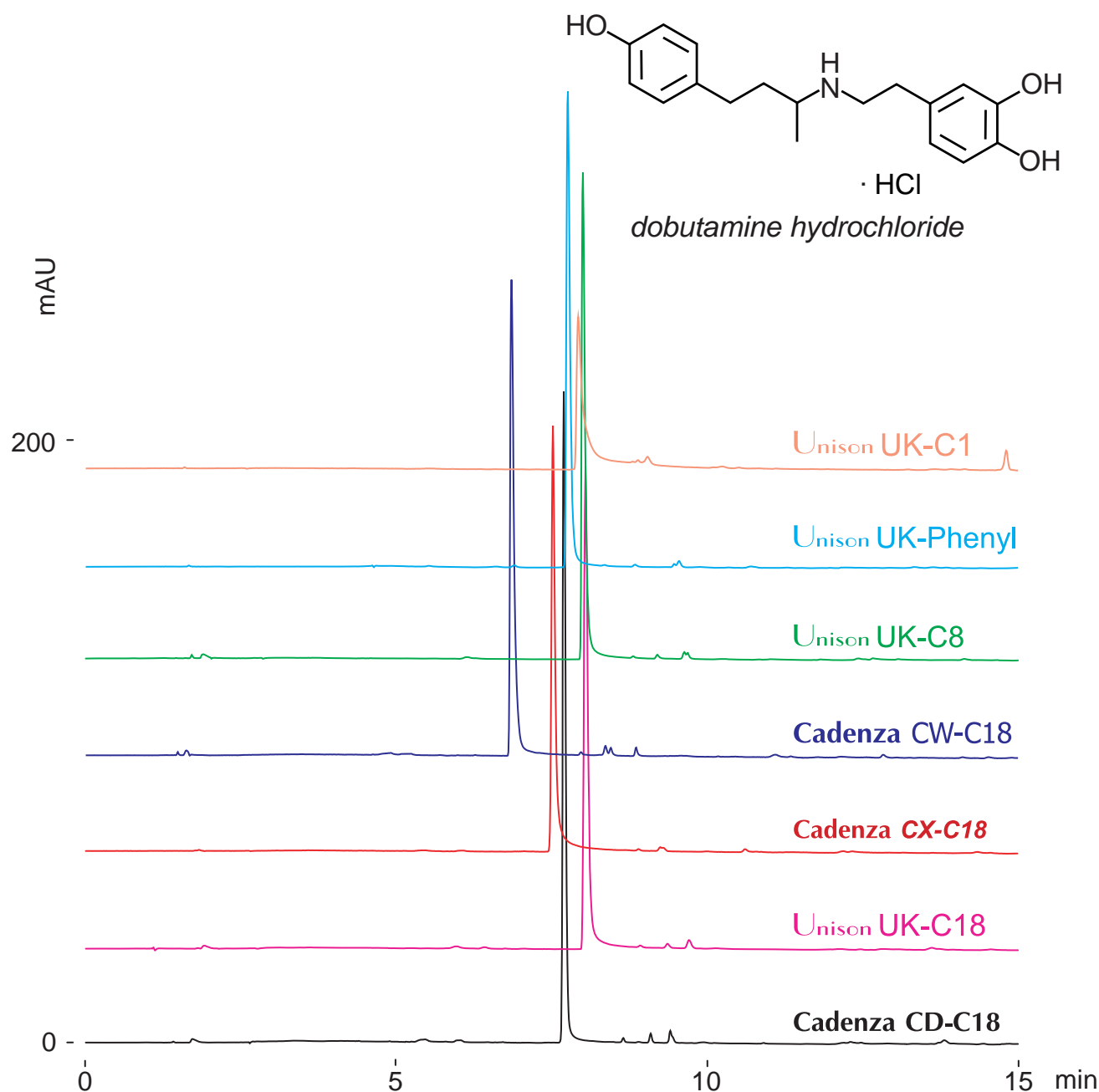

100 x 4.6 mm  
 A: 100mM ammonium formate  
 B: 100mM ammonium formate / acetonitrile = 50 / 50  
 0-100%B (0-15min)  
 1.0mL/min (12MPa), 37°C, 280nm  
 10µL (0.33mg/mL, 100mM ammonium formate)