

Unison UK-Amino  
Nardis ND-NX


150 x 4.6 mm

HPLC COLUMN APPLICATION

Emtricitabine

(Nucleoside Analogue Reverse Transcriptase Inhibitor:NRTI)

エムトリシタビン (核酸系逆転写酵素阻害剤)

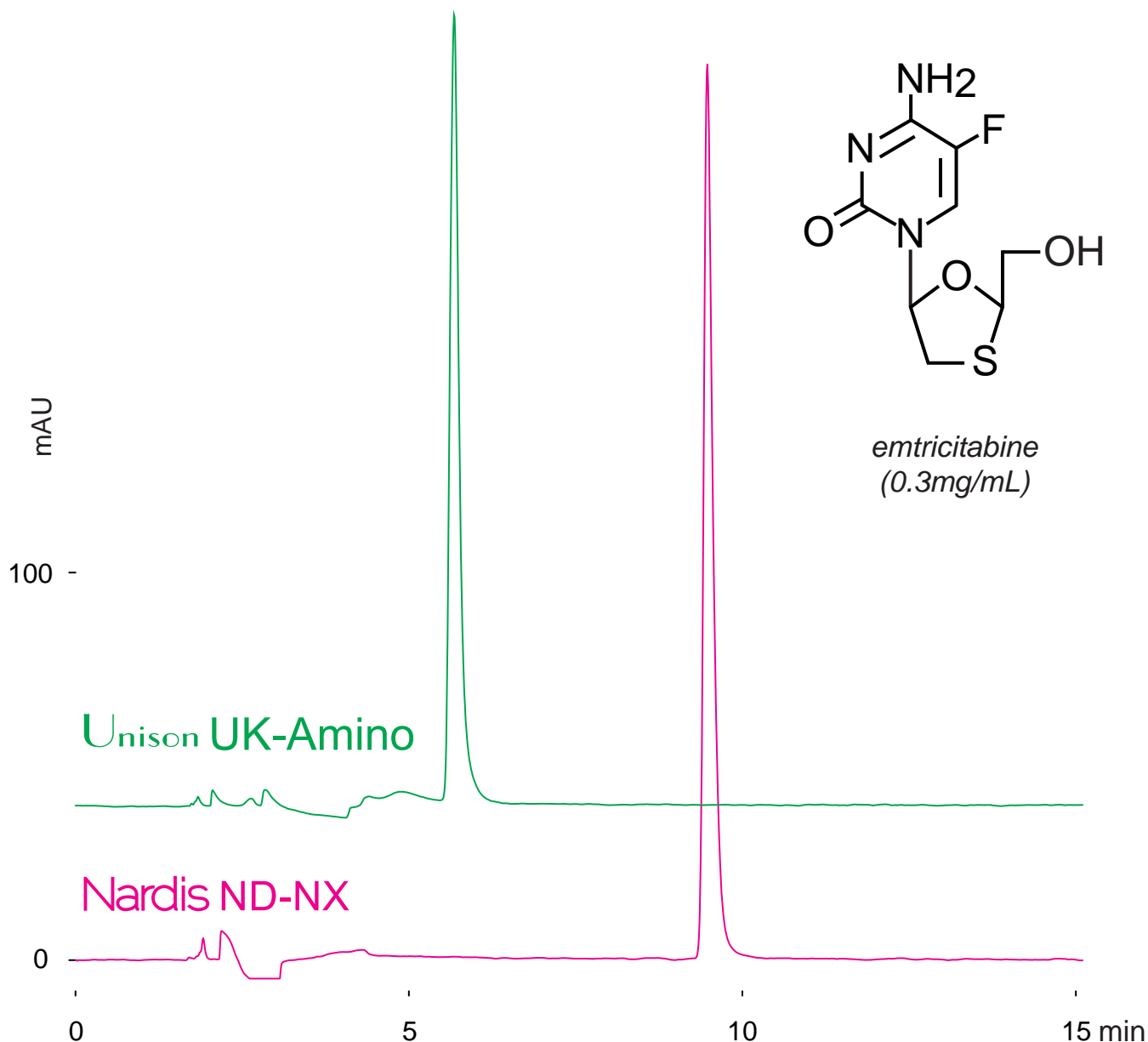

150 x 4.6 mm

A: acetonitrile / acetic acid = 100 / 0.5

B: water / acetic acid = 100 / 0.5

0-20%B (0-15min)

1.0mL/min (5-11MPa), 37°C, 260nm

10µL (0.3mg/mL, methanol)