


Unison UK-C1  
 Unison UK-Phenyl Cadenza CX-C18  
 Unison UK-C8 Cadenza CW-C18  
 Unison UK-C18 Cadenza CD-C18

100 x 4.6 mm


HPLC COLUMN APPLICATION

Phenolphthalein (pH indicator)

フェノールフタレイン ( pH指示薬 )

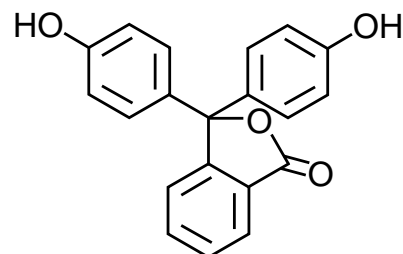

phenolphthalein (H<sub>2</sub>PP)  
 (solution color : colorless)

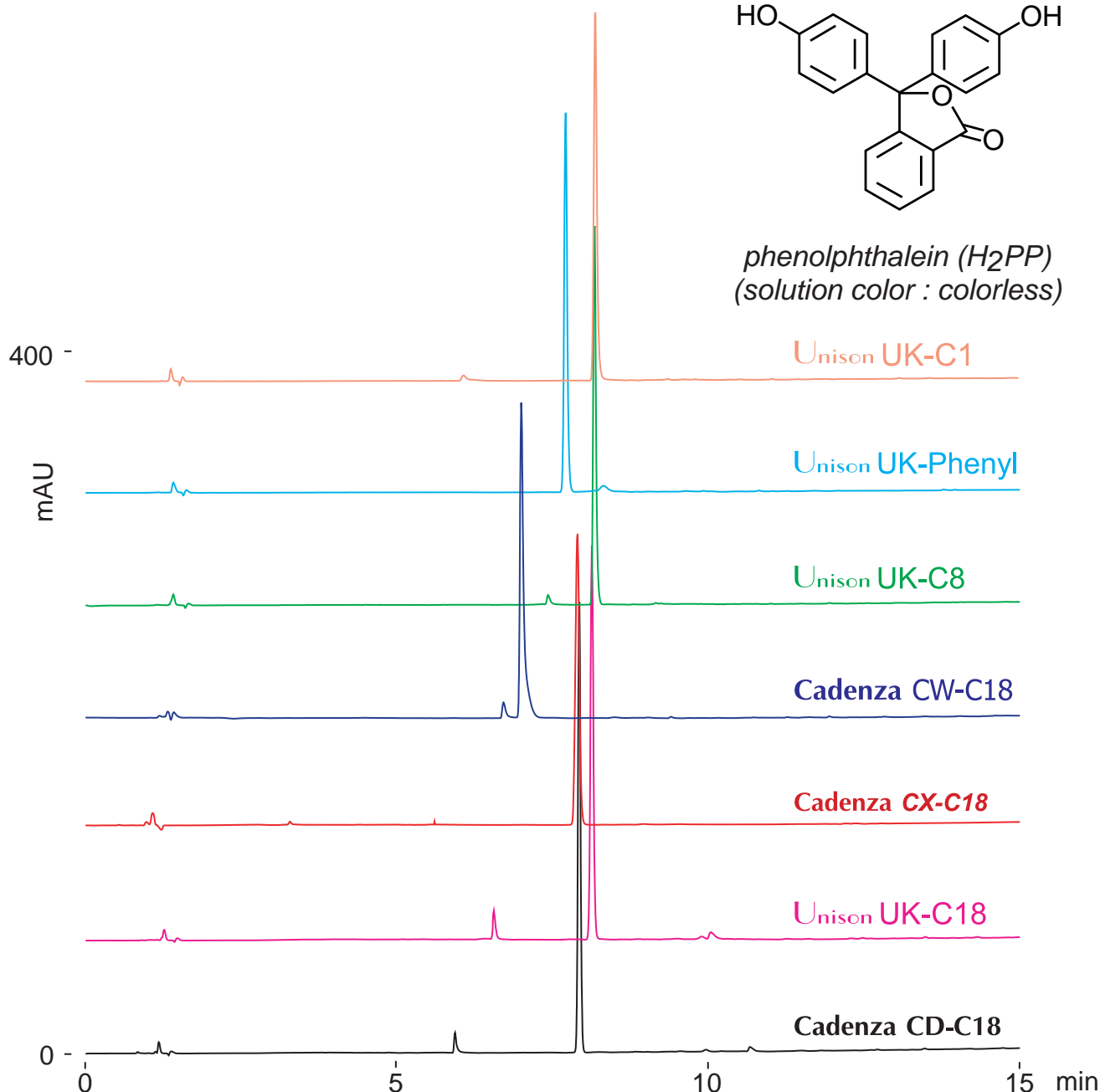

100 x 4.6 mm

A: water / acetonitrile / formic acid = 90 / 10 / 0.1

B: water / acetonitrile / formic acid = 10 / 90 / 0.1

0-100%B (0-15min)

1.0mL/min (11MPa), 37°C, 260nm

10µL (0.2mg/mL, water / acetonitrile / formic acid = 60 / 40 / 1)