



Unison UK-C1
 Unison UK-Phenyl Cadenza CX-C18
 Unison UK-C8 Cadenza CW-C18
 Unison UK-C18 Cadenza CD-C18

100 x 4.6 mm

HPLC COLUMN APPLICATION

Congo Red (pH indicator)

コンゴレッド (pH指示薬)

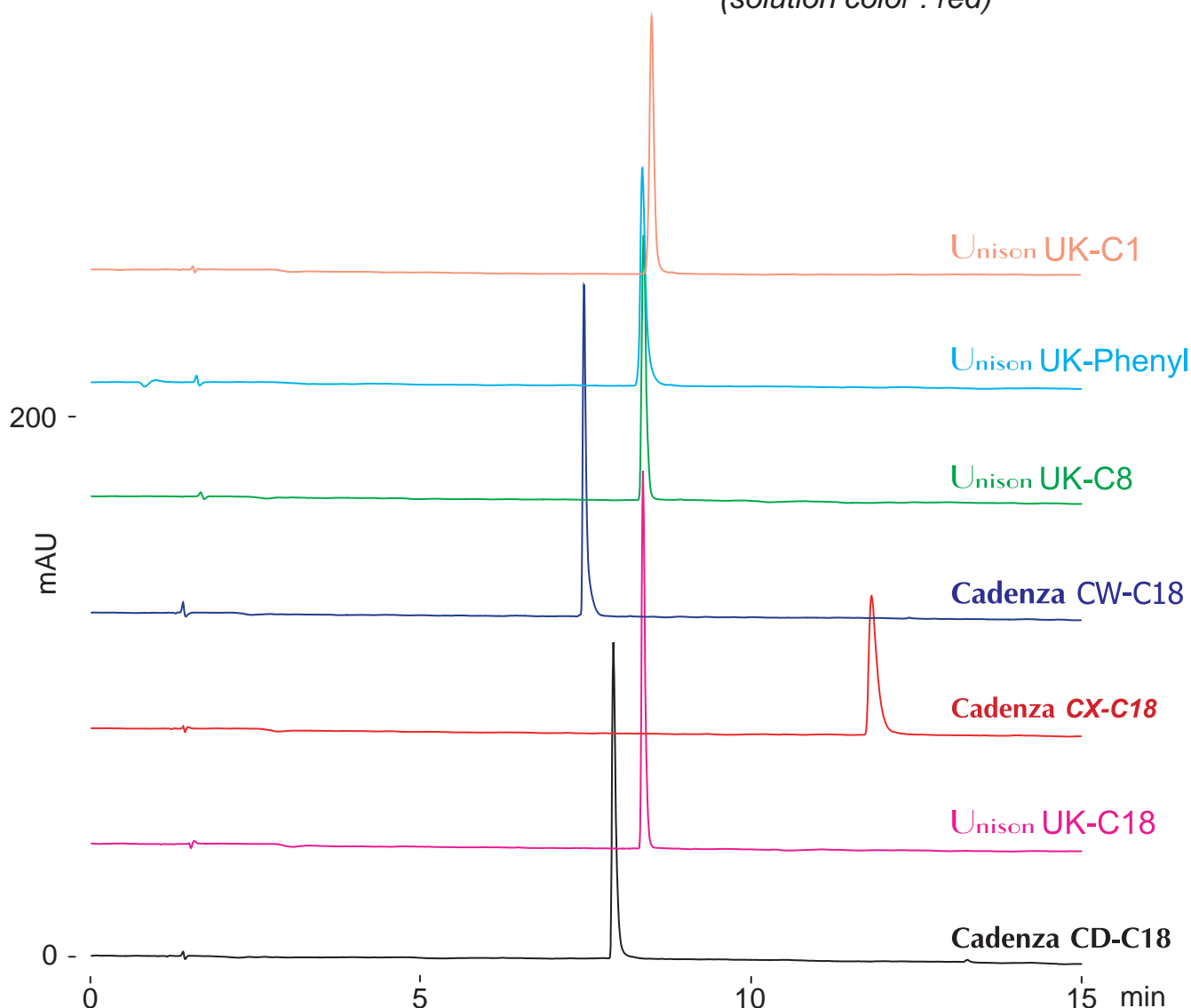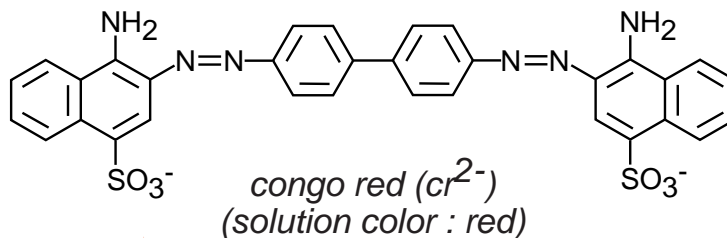

100 x 4.6 mm

A: 100mM ammonium formate / acetonitrile = 90 / 10

B: 100mM ammonium formate / acetonitrile = 10 / 90

0-70%B (0-15min)

1.0mL/min (14MPa), 37°C, 495nm

10µL (0.025mg/mL, 100mM ammonium formate)