


Scherzo SM-C18
 Scherzo SW-C18
 Scherzo SS-C18

150 x 3 mm

HPLC COLUMN APPLICATION

LC-MS: Guanidine hydrochloride

塩酸グアニジンの LC-MS 分析

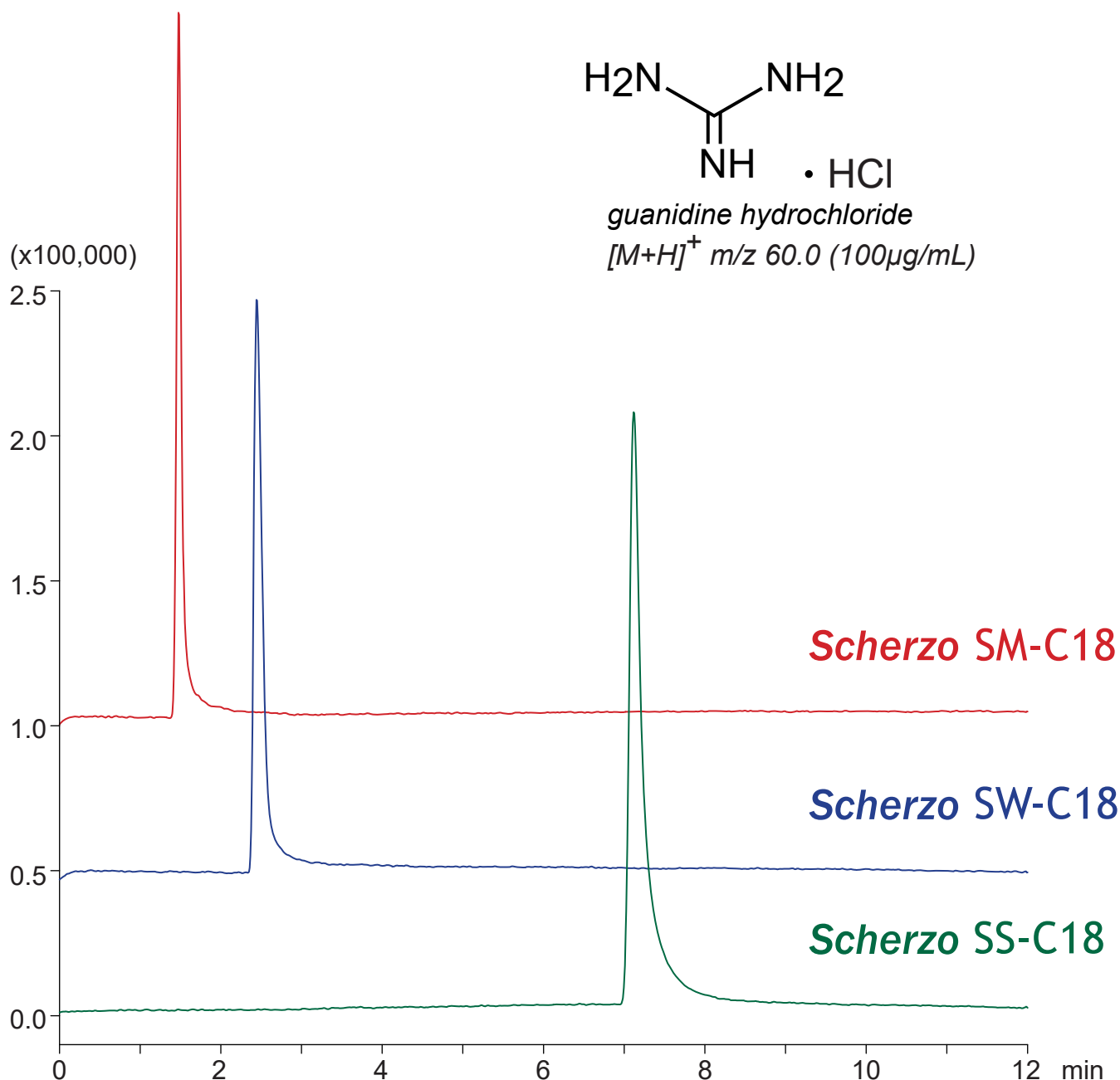

150 x 3 mm

A: water / formic acid = 100 / 0.1

B: 100mM ammonium formate / acetonitrile = 50 / 50

10-100%B (0-10min), 100%B (10-12min)

0.4mL/min (10-11MPa), 37°C, Single Quad. MS (ESI)

2µL (100µg/mL, water / formic acid = 100 / 1)