


**Nardis ND-RX**  
**Scherzo SM-C18**

150 x 3 mm

**HPLC COLUMN APPLICATION**

**LC-MS: Guanidine hydrochloride**

塩酸グアニジンの LC-MS 分析

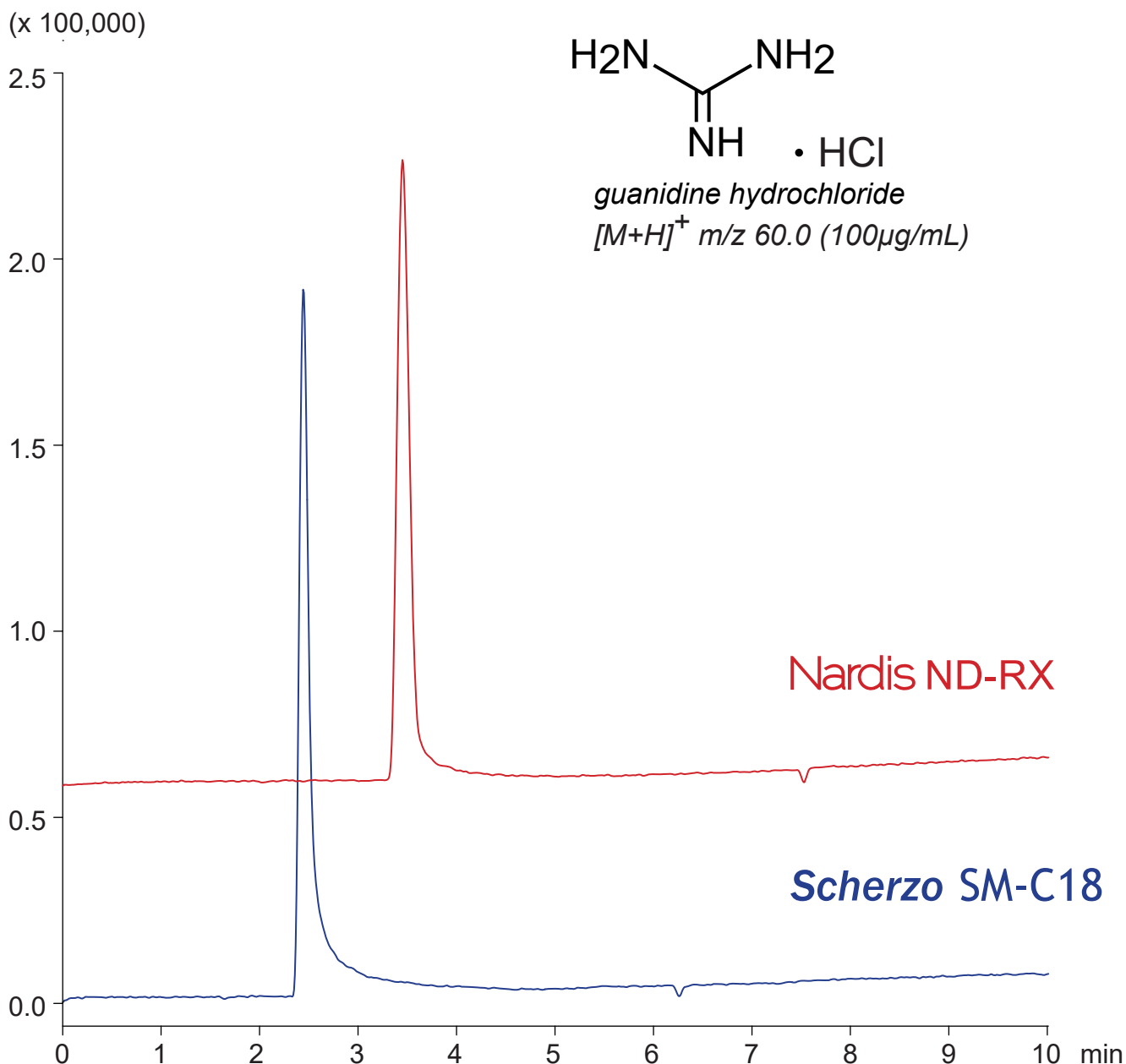

150 x 3 mm

A: acetonitrile / 5mM ammonium formate = 10 / 90

B: acetonitrile / water / formic acid = 20 / 80 / 0.5

0%B (0-1.5min), 0-50%B (1.5-10min)

0.4mL/min (10-11MPa), 37°C, Single Quad. MS (ESI)

2µL (100µg/mL, water / formic acid = 100 / 1)