


Scherzo SS-C18
 Scherzo SW-C18
 Scherzo SM-C18

150 x 4.6 mm

HPLC COLUMN APPLICATION

Guanidine hydrochloride (protein denaturant)

グアニジン塩酸塩 (タンパク質変性剤)

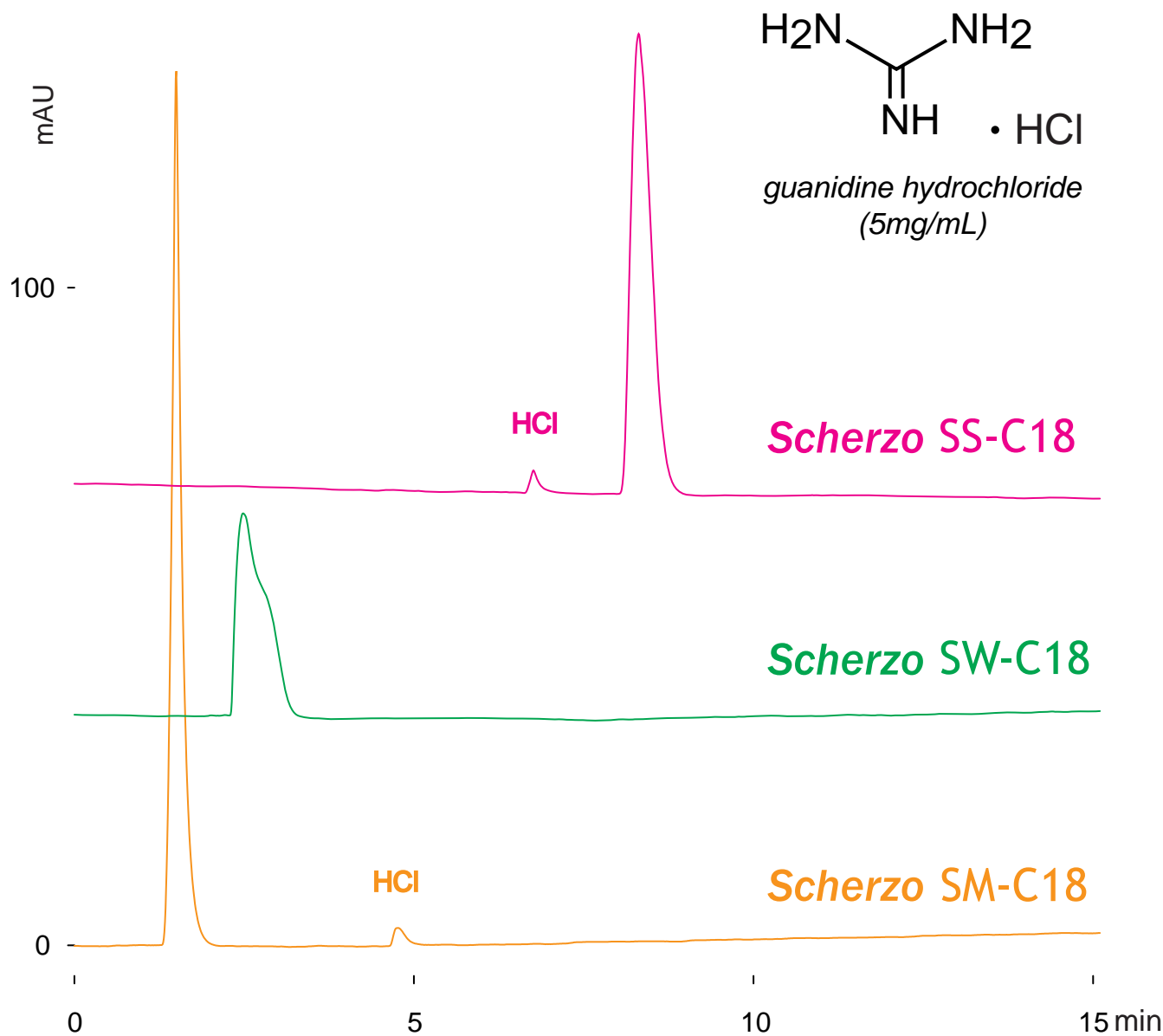

150 x 4.6 mm

A: water / formic acid = 100 / 0.1

B: 100mM ammonium formate / acetonitrile = 50 / 50

0-100%B (0-15min)

1.0mL/min (12-13MPa), 37°C

ELSD (spray chamber 70°C, drift chamber 100°C)

10µL (5mg/mL, water)