

Unison UK-Amino


Nardis ND-NX

150 x 4.6 mm

HPLC COLUMN APPLICATION

Idoxuridine (Antih herpesvirus drugs)

イドクスウリジン (抗ヘルペスウイルス薬)

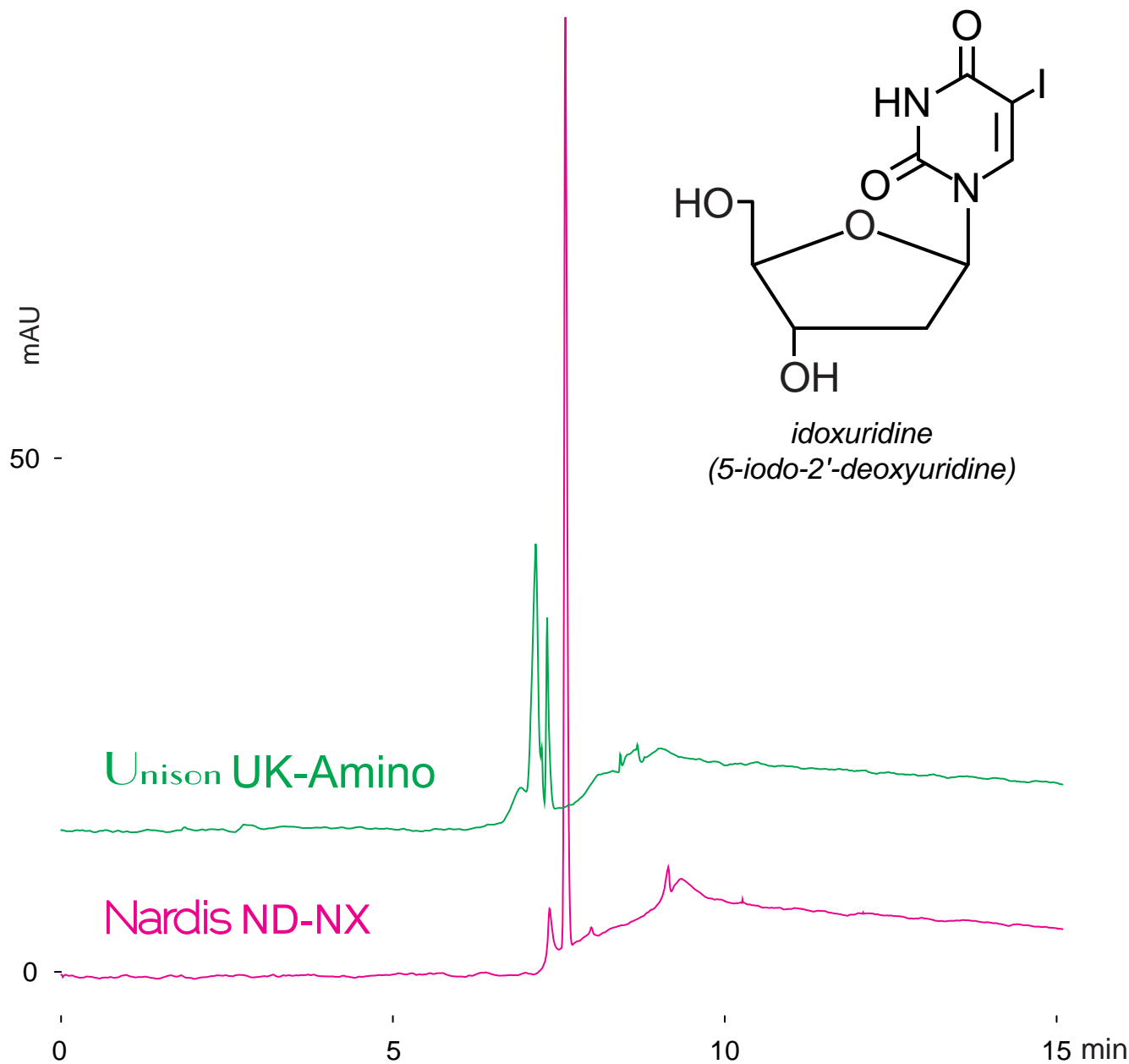

150 x 4.6 mm
 A: acetonitrile
 B: 100mM ammonium formate
 0%B (0-3min), 0-100%B (3-15min)
 1.0mL/min (5-7MPa), 37°C, 260nm
 10µL (0.1mg/mL, 1M sodium hydroxide)