






Cadenza CL-C18
Cadenza CD-C18
Unison UK-C18

250 x 4.6 mm

HPLC COLUMN APPLICATION

Fennel oil

ウイキョウ油

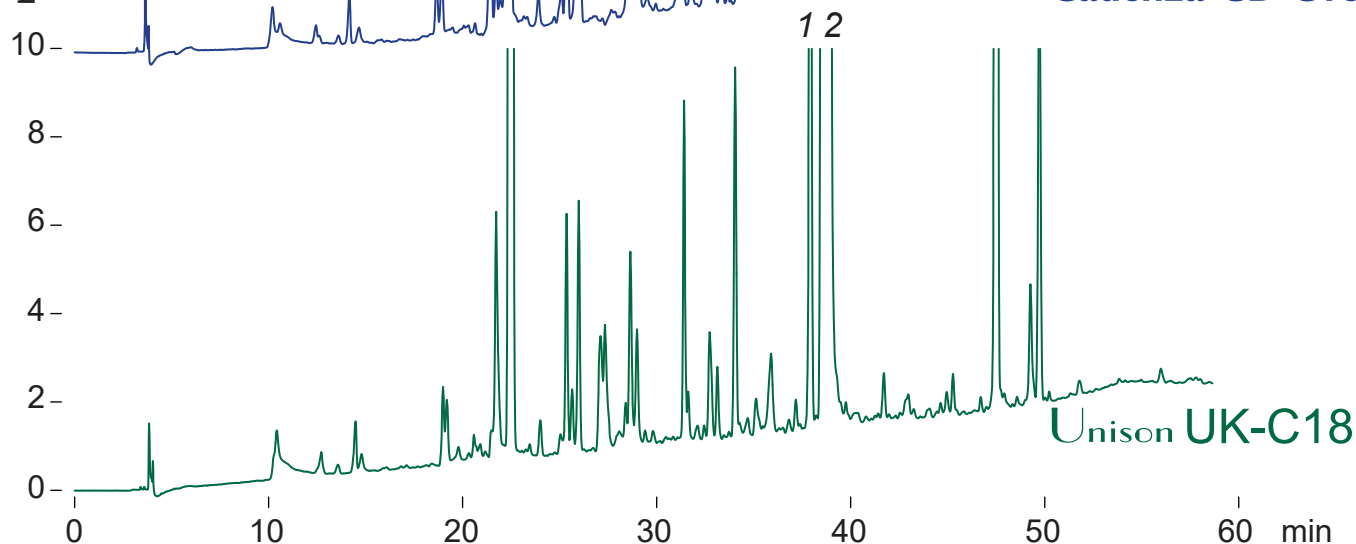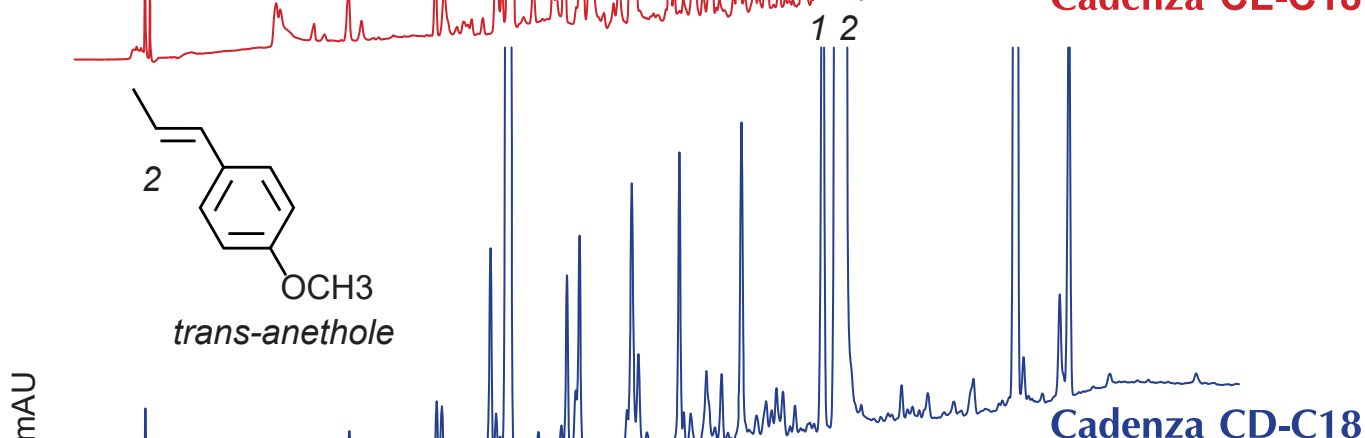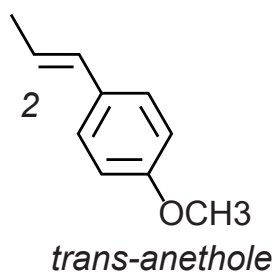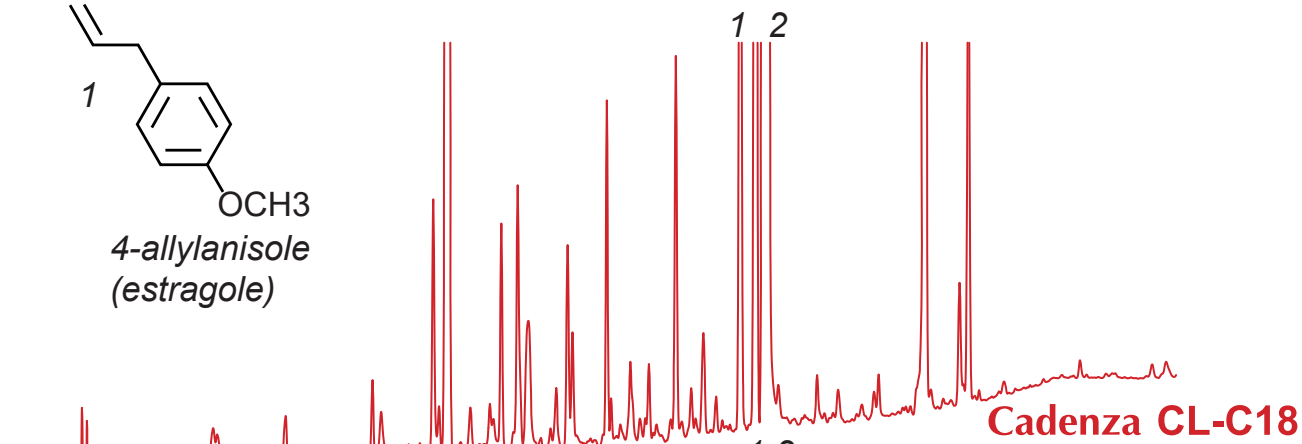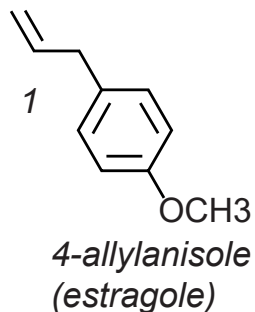

250 x 4.6 mm
 A: acetonitrile / water / formic acid = 10 / 90 / 0.1
 B: acetonitrile / water / formic acid = 95 / 5 / 0.09
 0-100%B (0-50min), 100%B (50-60min)
 0.8mL/min (14-15MPa), 37°C, 260nm
 10µL (10mg/mL, ethanol / formic acid = 100 / 1)