



Scherzo SM-C18 MF

100 x 2 mm

Application

LC-MS: Organophosphorus, quaternary ammonium and bis-guanidine herbicides / fungicides

有機リン系，四級アンモニウム系及びビスグアニジン系農薬の LC-MS 分析

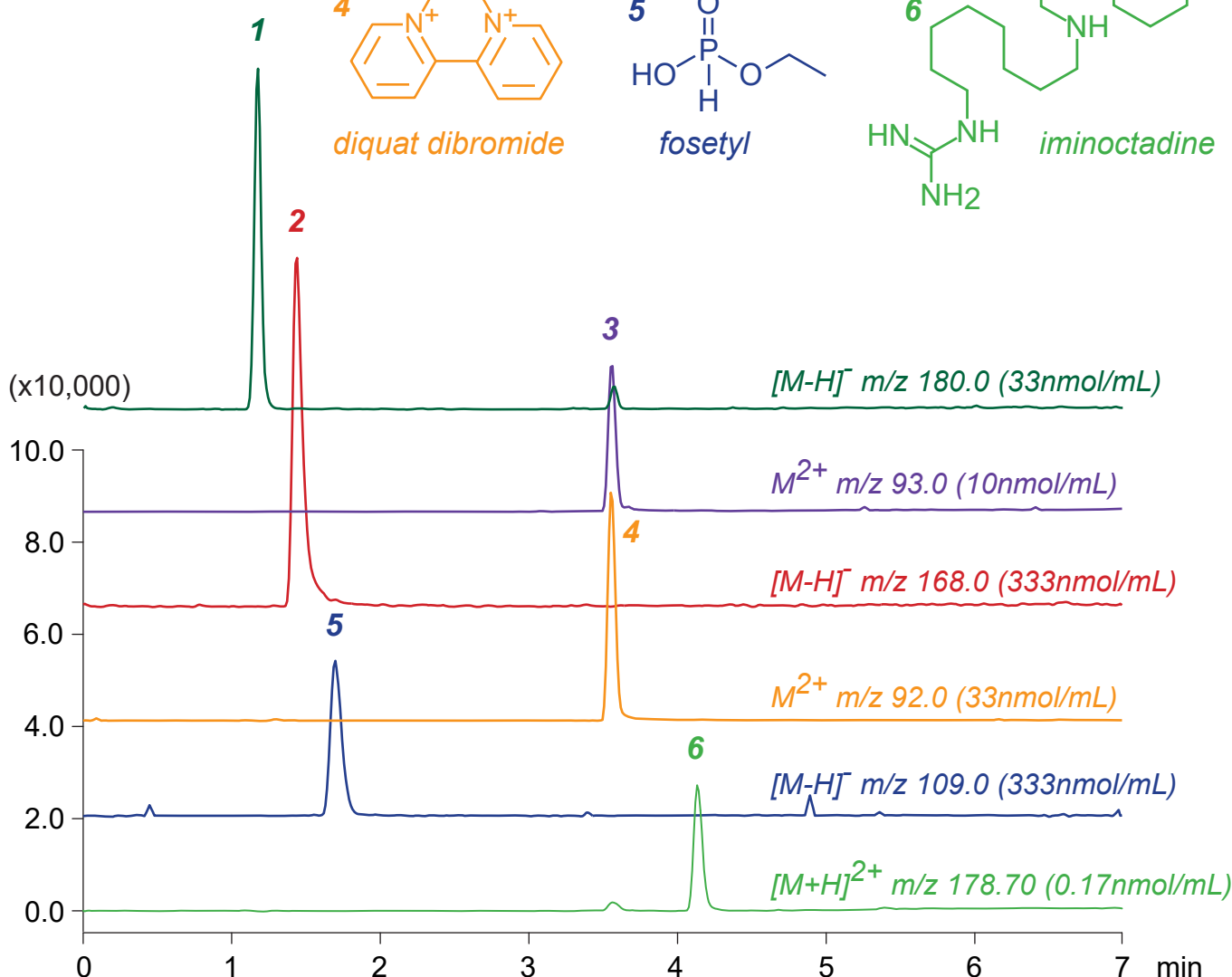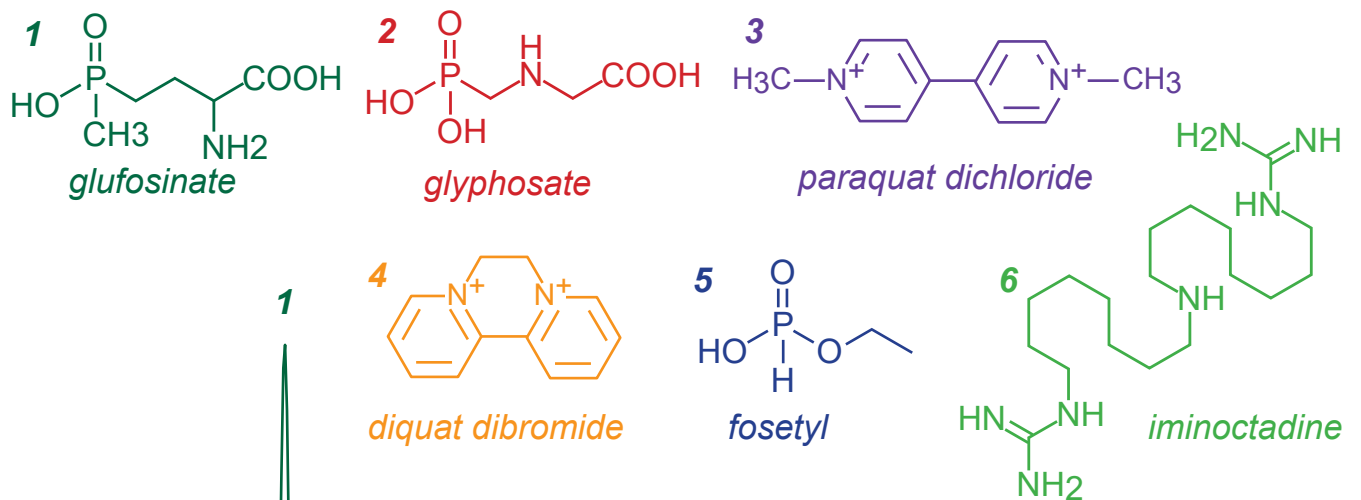

Scherzo SM-C18 MF, 100 x 2 mm
 A: 30mM ammonium formate
 B: acetonitrile / water / formic acid = 30 / 70 / 0.5
 0%B (0-1min), 0-100%B (1-4min), 100%B (4-7min)
 0.2 mL/min (7MPa), 37°C, Single Quad. MS (ESI)
 5 μL (0.17-333nmol/mL, 15mM HCOONH₄)