


Unison UK-C1
 Unison UK-Phenyl Cadenza CX-C18
 Unison UK-C8 Cadenza CW-C18
 Unison UK-C18 Cadenza CD-C18

100 x 4.6 mm


HPLC COLUMN APPLICATION

Variamine Blue B (Redox Indicator)

バリアミンブルー B (酸化還元指示薬)

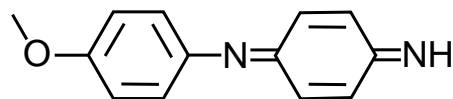

variamine blue b
(oxidized form, blue-purple)

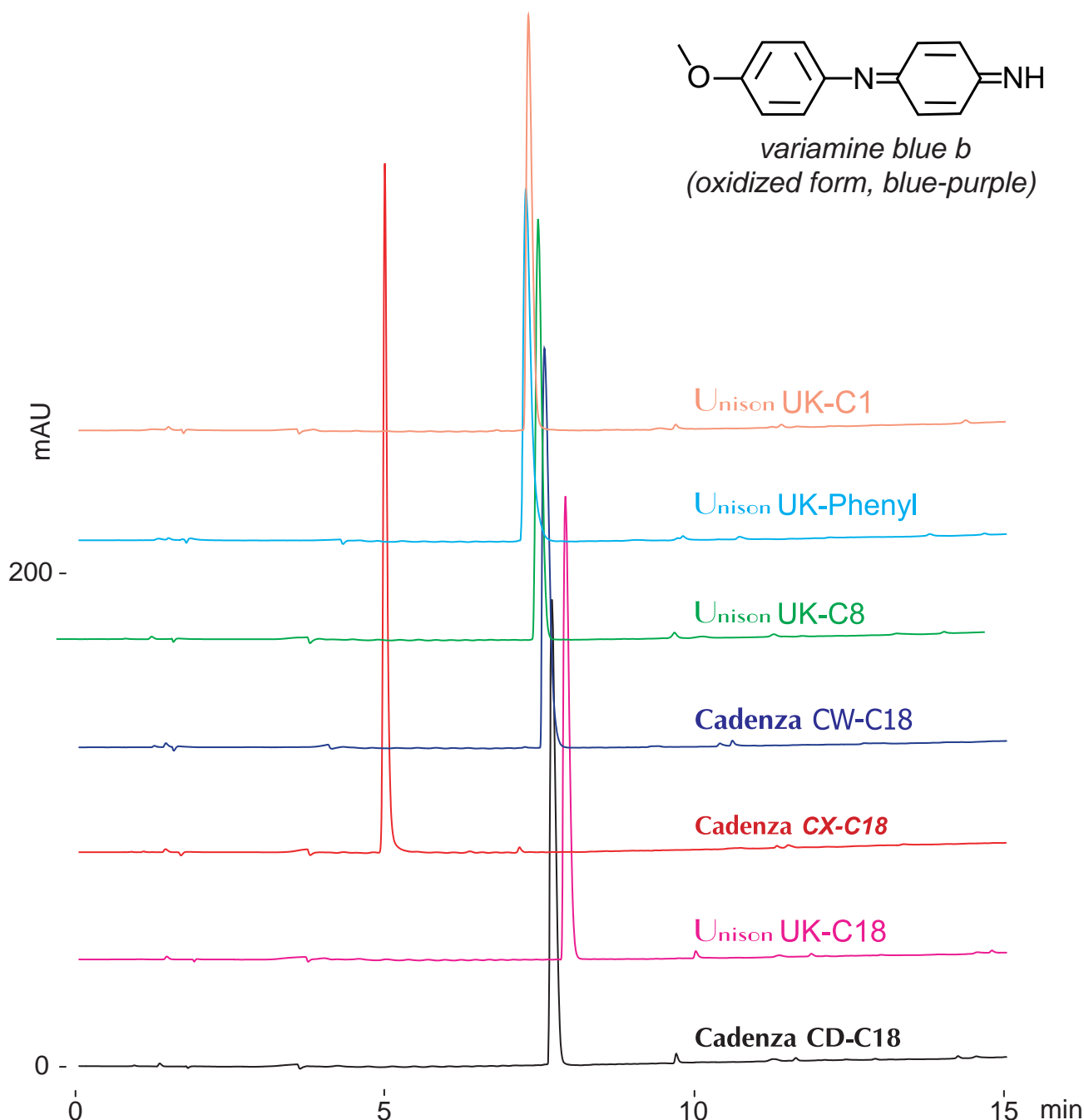

100 x 4.6 mm

A: water / formic acid = 100 / 0.5

B: water / acetonitrile / formic acid = 10 / 90 / 0.5

0%B (0-1min), 0-70%B (1-15min)

1.0mL/min (12-17MPa), 37°C, 260nm

10µL (0.1mg/mL, water / acetonitrile / formic acid = 80 / 20 / 1)