


Unison UK-C1  
 Unison UK-Phenyl Cadenza CX-C18  
 Unison UK-C8 Cadenza CW-C18  
 Unison UK-C18 Cadenza CD-C18

100 x 4.6 mm


HPLC COLUMN APPLICATION

## Variamine Blue B (Redox Indicator)

バリアミンブルー B (酸化還元指示薬)

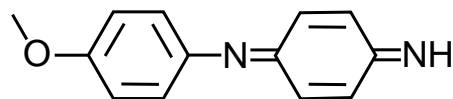

*variamine blue b*  
(oxidized form, blue-purple)

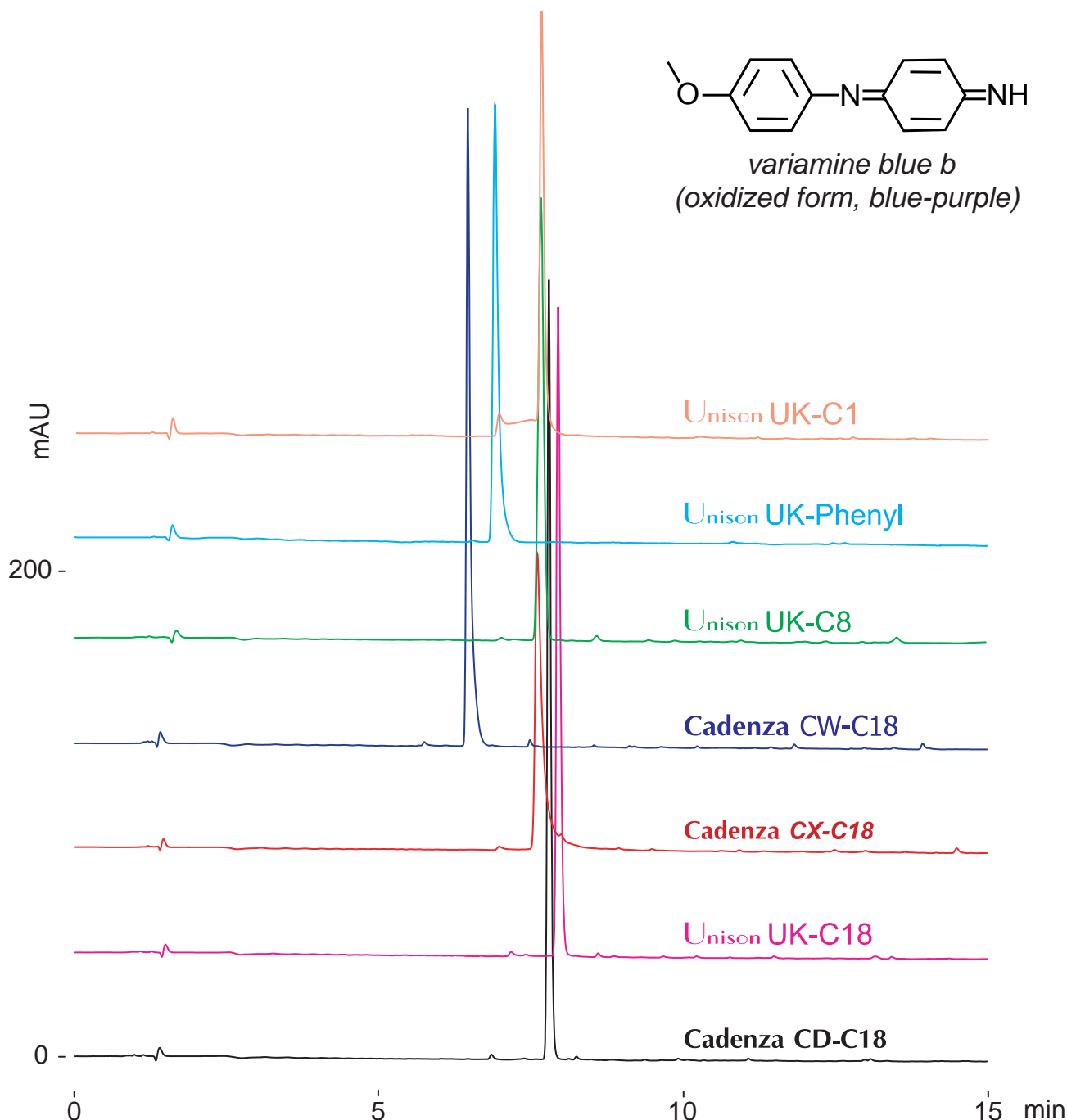

100 x 4.6 mm

A: 100mM ammonium formate / acetonitrile = 90 / 10

B: 100mM ammonium formate / acetonitrile = 10 / 90

0-100%B (0-15min)

1.0mL/min (11-19MPa), 37°C, 260nm

10µL (0.1mg/mL, 100mM ammonium formate / acetonitrile = 60 / 40)