


Unison UK-Phenyl Unison UK-C1  
 Unison UK-C8 Cadenza CW-C18  
 Unison UK-C18 Cadenza CD-C18

100 x 4.6 mm


HPLC COLUMN APPLICATION

**2-Hydroxy-1- (2-hydroxy-4-sulfo-1-naphthylazo)  
 -3-naphthoic acid (NN, Metal Indicator)**

2-ヒドロキシ -1-( 2-ヒドロキシ -4-スルホ -1ナフチルアゾ ) -3ナフトエ酸  
 ( NN, 金属指示薬 )

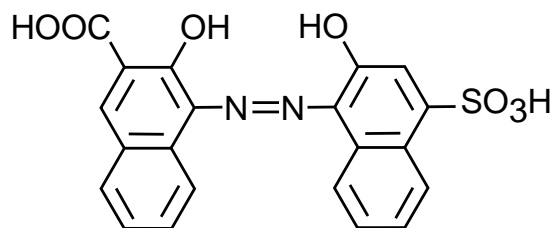

2-hydroxy-1-  
 (2-hydroxy-4-sulfo-1-naphthylazo)  
 -3-naphthoic acid (NN)  
 (solution color : purple)

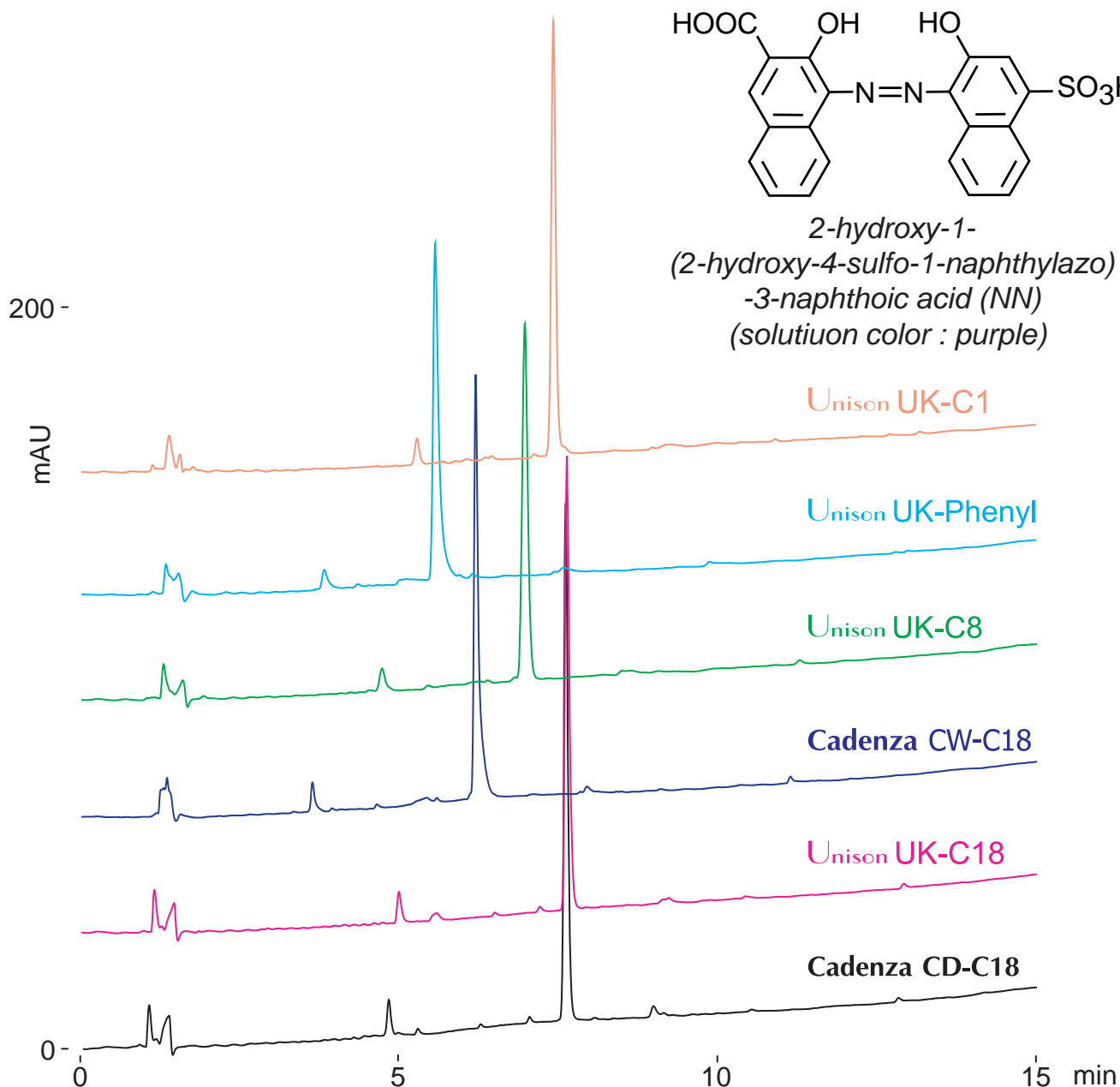

100 x 4.6 mm

A: water / acetonitrile / formic acid = 90 / 10 / 0.1

B: water / acetonitrile / formic acid = 10 / 90 / 0.1

10-100%B (0-15min)

1.0mL/min (13-19MPa), 37°C, 250nm

10µL (2mg/mL, water / methanol / formic acid = 50 / 50 / 1)