


Unison UK-C1
 Unison UK-Phenyl Cadenza CX-C18
 Unison UK-C8 Cadenza CW-C18
 Unison UK-C18 Cadenza CD-C18

100 x 4.6 mm

HPLC COLUMN APPLICATION


Xylenol Orange (Metal Indicator)

キシレノールオレンジ (金属指示薬)

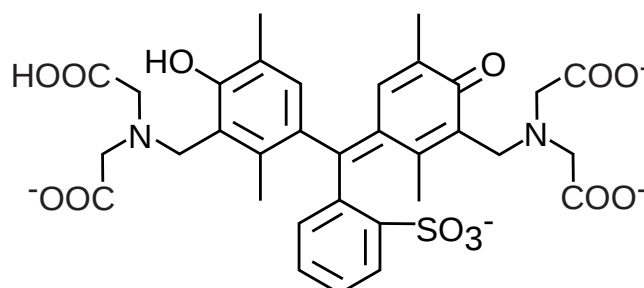

xylenol orange

(solution color : reddish purple)

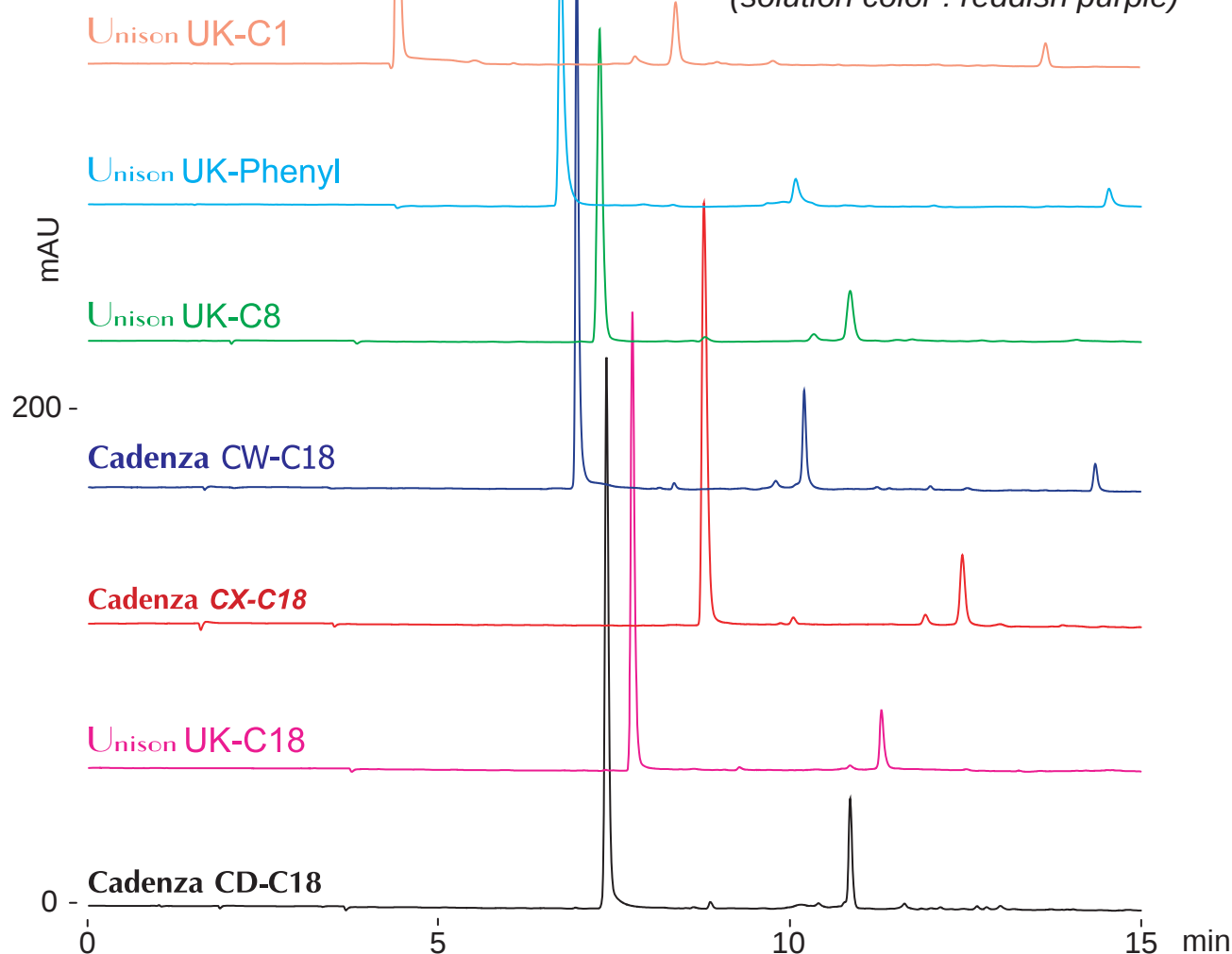

100 x 4.6 mm

A: 100mM ammonium formate

B: 100mM ammonium formate / acetonitrile = 50 / 50

0%B (0-1min), 0-60%B (1-15min)

1.0mL/min (12-18MPa), 37°C, 435nm

10μL (0.25mg/mL, 100mM ammonium formate)