


Unison UK-Phenyl Unison UK-C1
Unison UK-C8 Cadenza CW-C18
Unison UK-C18 Cadenza CD-C18

100 x 4.6 mm


HPLC COLUMN APPLICATION

Eriochrome Black T (Metal Indicator)

エリオクロムブラック T (金属指示薬)

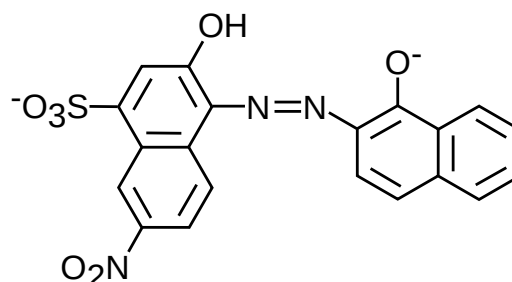

eriochrome black t
(solution color : blue)

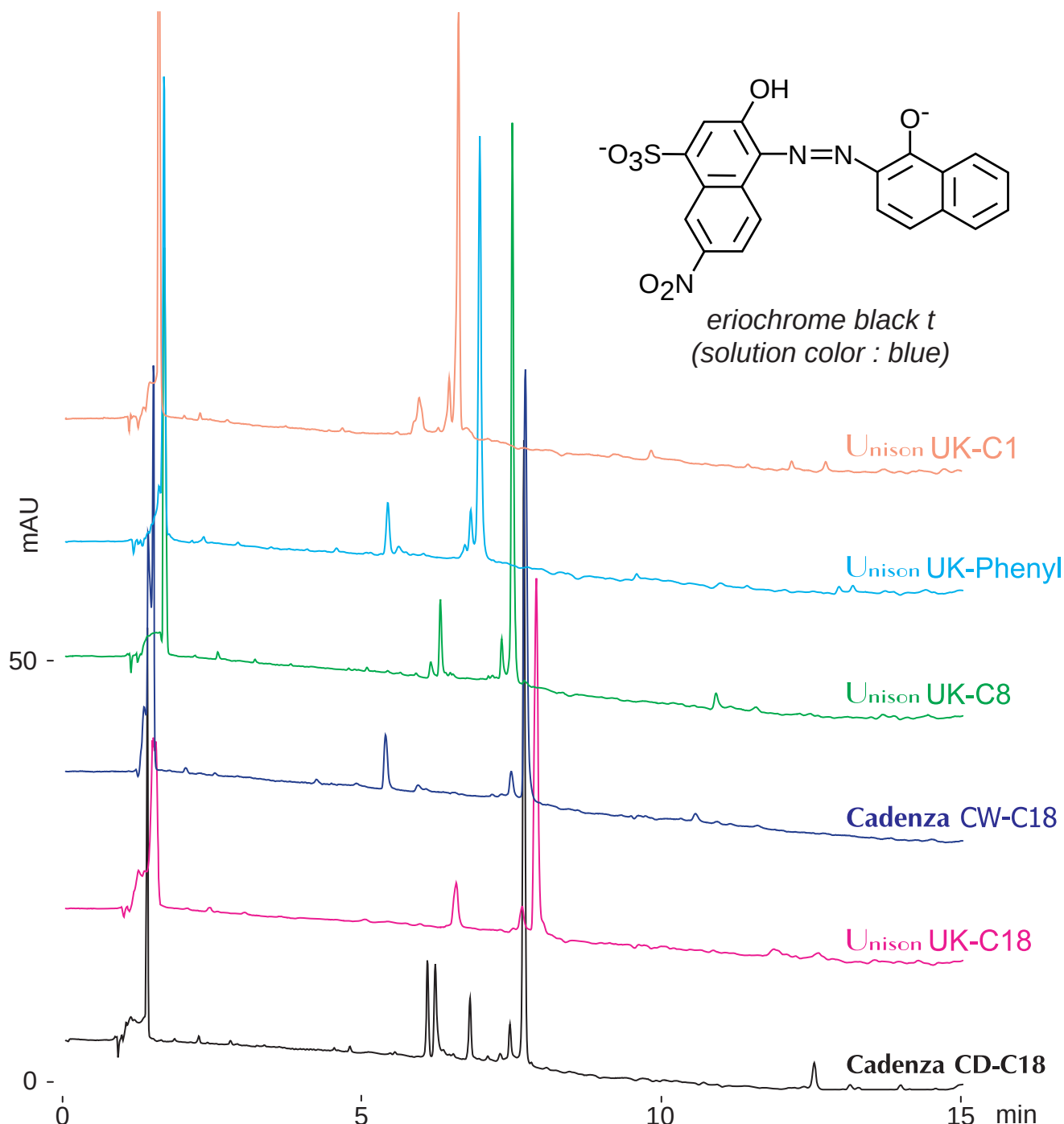

100 x 4.6 mm

A: 100mM ammonium formate / acetonitrile = 90 / 10

B: 100mM ammonium formate / acetonitrile = 10 / 90

15-100%B (0-15min)

1.0mL/min (10-13MPa), 37°C, 240nm

10μL (2mg/mL, 100mM ammonium formate / methanol = 50 / 50)