

Cadenza CD-C18 Unison UK-C8



Unison UK-C18 Unison UK-Phenyl

Cadenza CW-C18 Unison UK-C1 100 x 4.6 mm

HPLC COLUMN APPLICATION

Valsartan (Antihypertensive)

バルサルタン (高血圧治療薬)

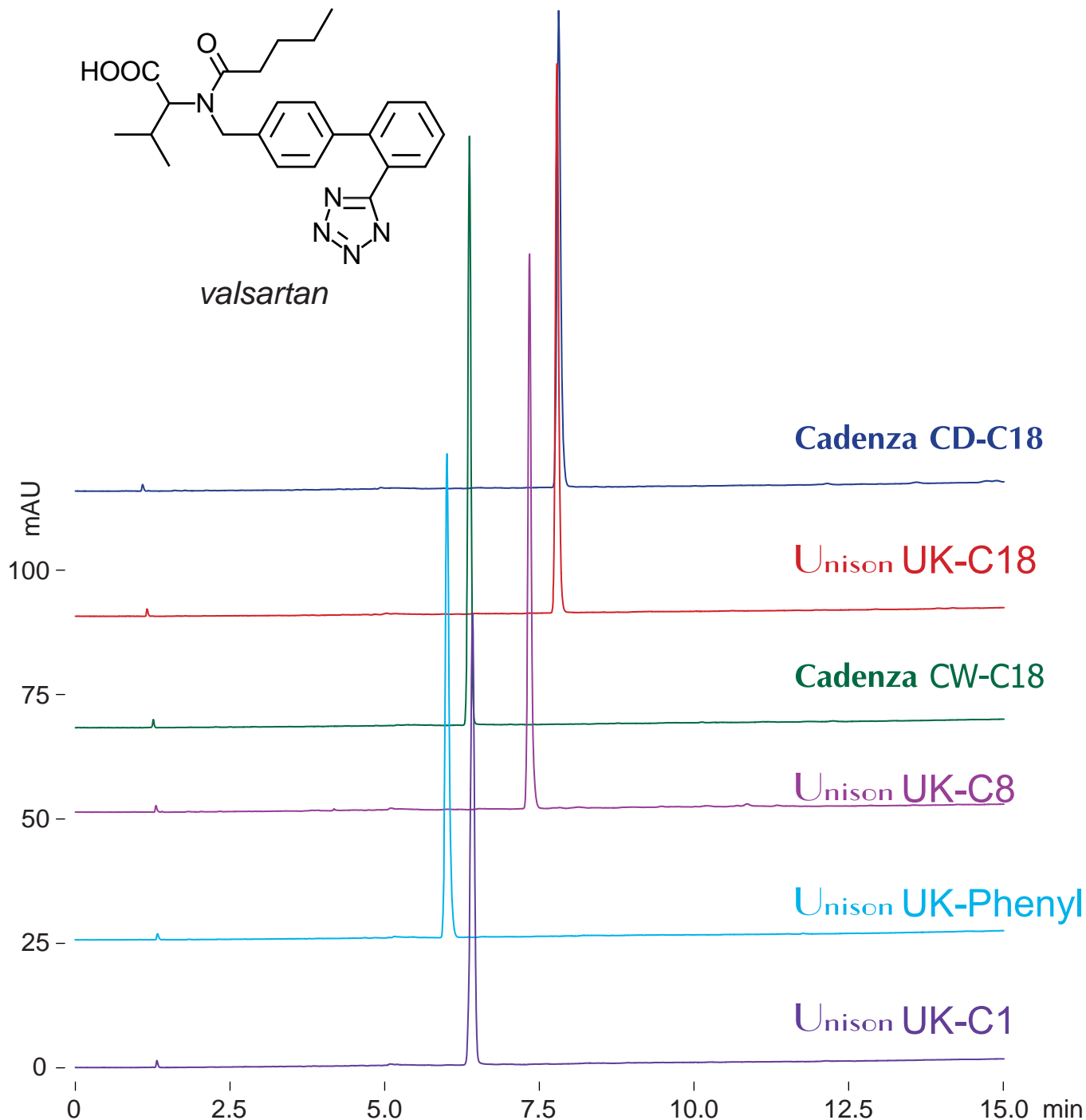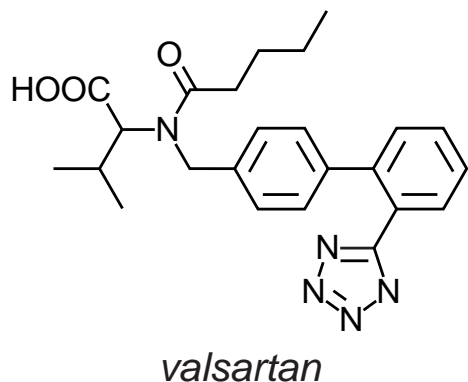

100 x 4.6 mm

A: acetonitrile / water / formic acid = 10 / 90 / 0.1

B: acetonitrile / water / formic acid = 90 / 10 / 0.1

20-100%B (0-15min)

1mL/min (12-13MPa), 37°C, 260nm

5µL (100µg/mL, acetonitrile / water / formic acid = 40 / 60 / 1)