

Cadenza CD-C18 Unison UK-C8


Unison UK-C18 Unison UK-Phenyl

Cadenza CW-C18 Unison UK-C1 100 x 4.6 mm


HPLC COLUMN APPLICATION

Valsartan (Antihypertensive)

バルサルタン (高血圧治療薬)

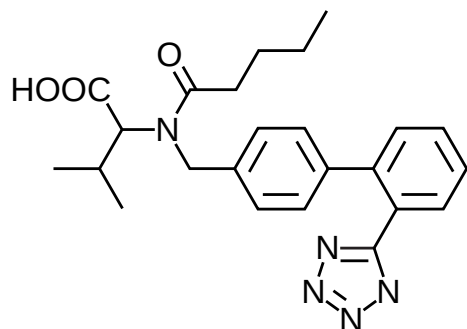

valsartan

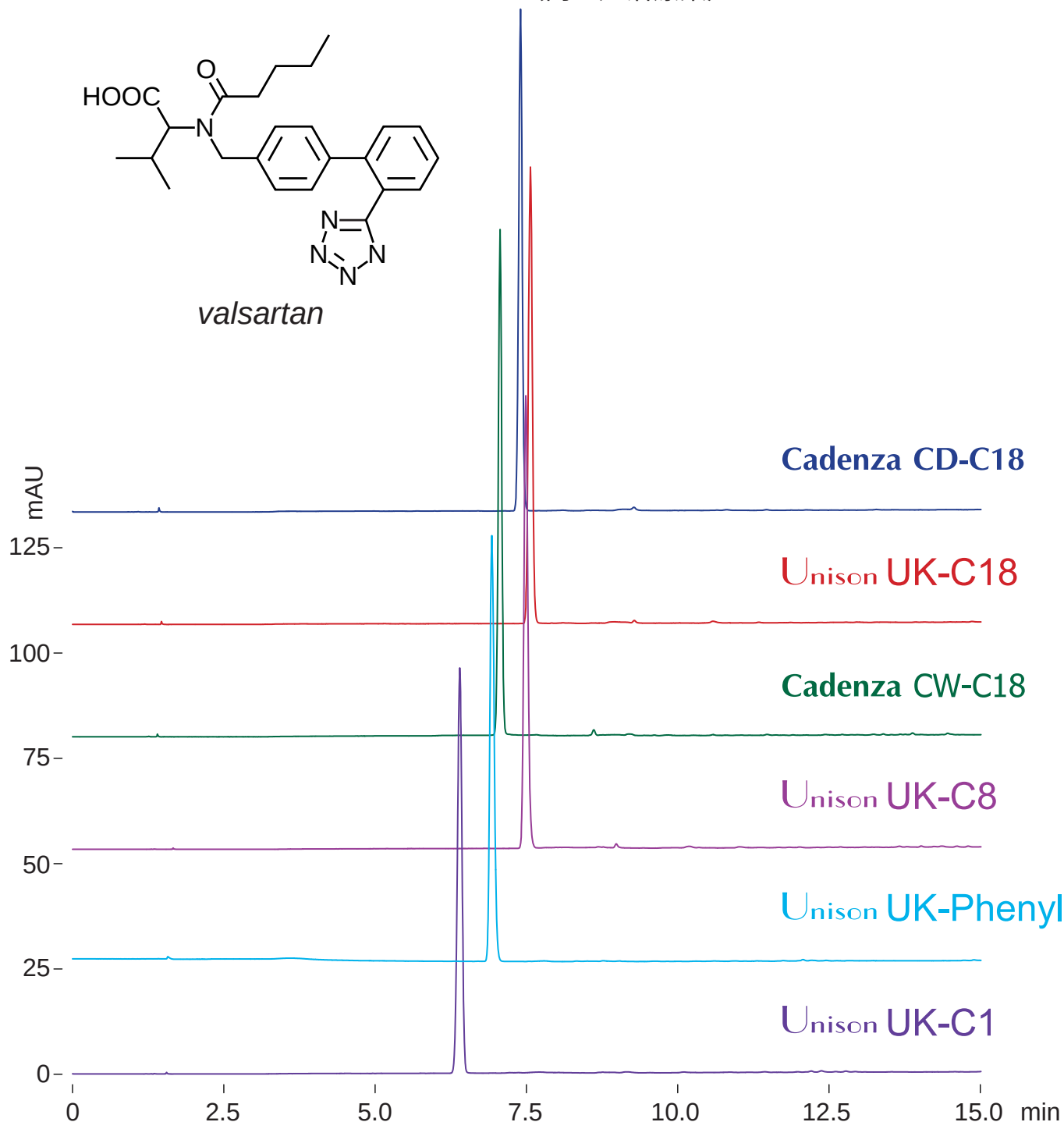

100 x 4.6 mm

A: acetonitrile / 100mM ammonium formate = 5 / 95

B: acetonitrile / 100mM ammonium formate = 90 / 10

0%B (0-1min), 0-90%B (1-15min)

1mL/min (10-13MPa), 37°C, 260nm

5μL (100μg/mL, acetonitrile / 100mM ammonium formate = 40 / 60)