

Cadenza CD-C18 Unison UK-C8


Unison UK-C18 Unison UK-Phenyl

Cadenza CW-C18 Unison UK-C1 100 x 4.6 mm


HPLC COLUMN APPLICATION

Cilostazol (Antithrombotic)

シロスタゾール (抗血栓薬)

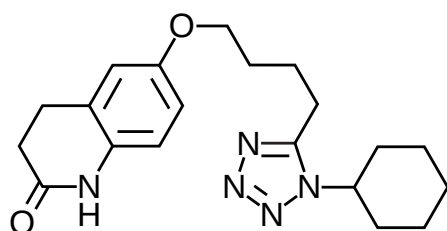

cilostazol

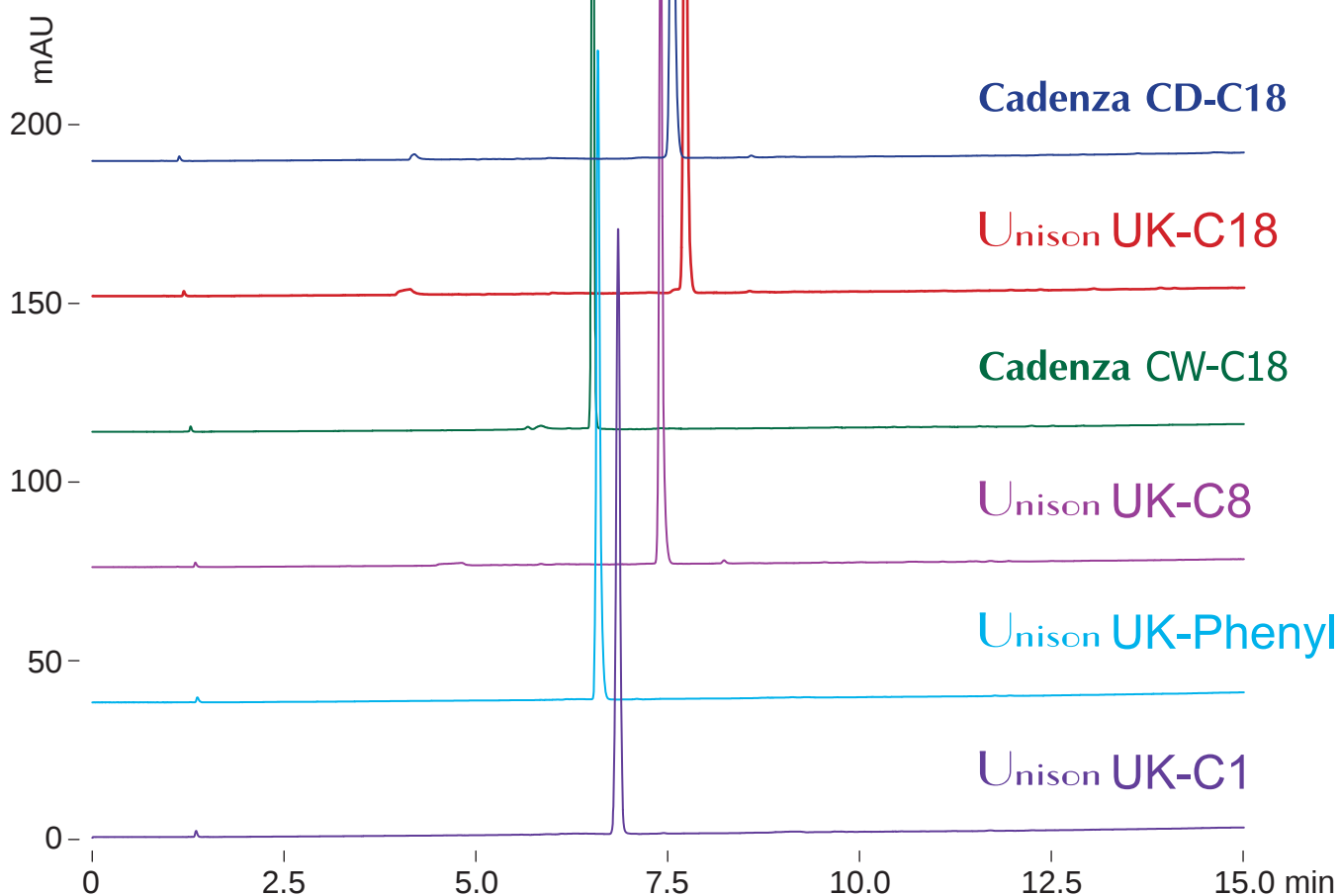

100 x 4.6 mm

A: acetonitrile / water / formic acid = 10 / 90 / 0.1

B: acetonitrile / water / formic acid = 95 / 5 / 0.1

10-100%B (0-15min)

1mL/min (9-12MPa), 37°C, 260nm

5μL (100μg/mL, acetonitrile / water / formic acid = 50 / 50 / 1)