

Unison UK-C8

Unison UK-C18


Cadenza CD-C18

100 x 4.6 mm


HPLC COLUMN APPLICATION

Sulfonazo III (Precipitation titration indicator)

スルホナゾIII (沈殿滴定指示薬)

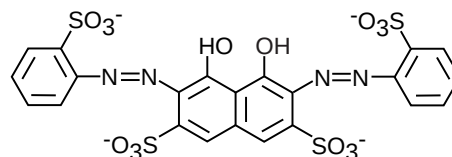

sulfonazo III
(solution color : red purple)

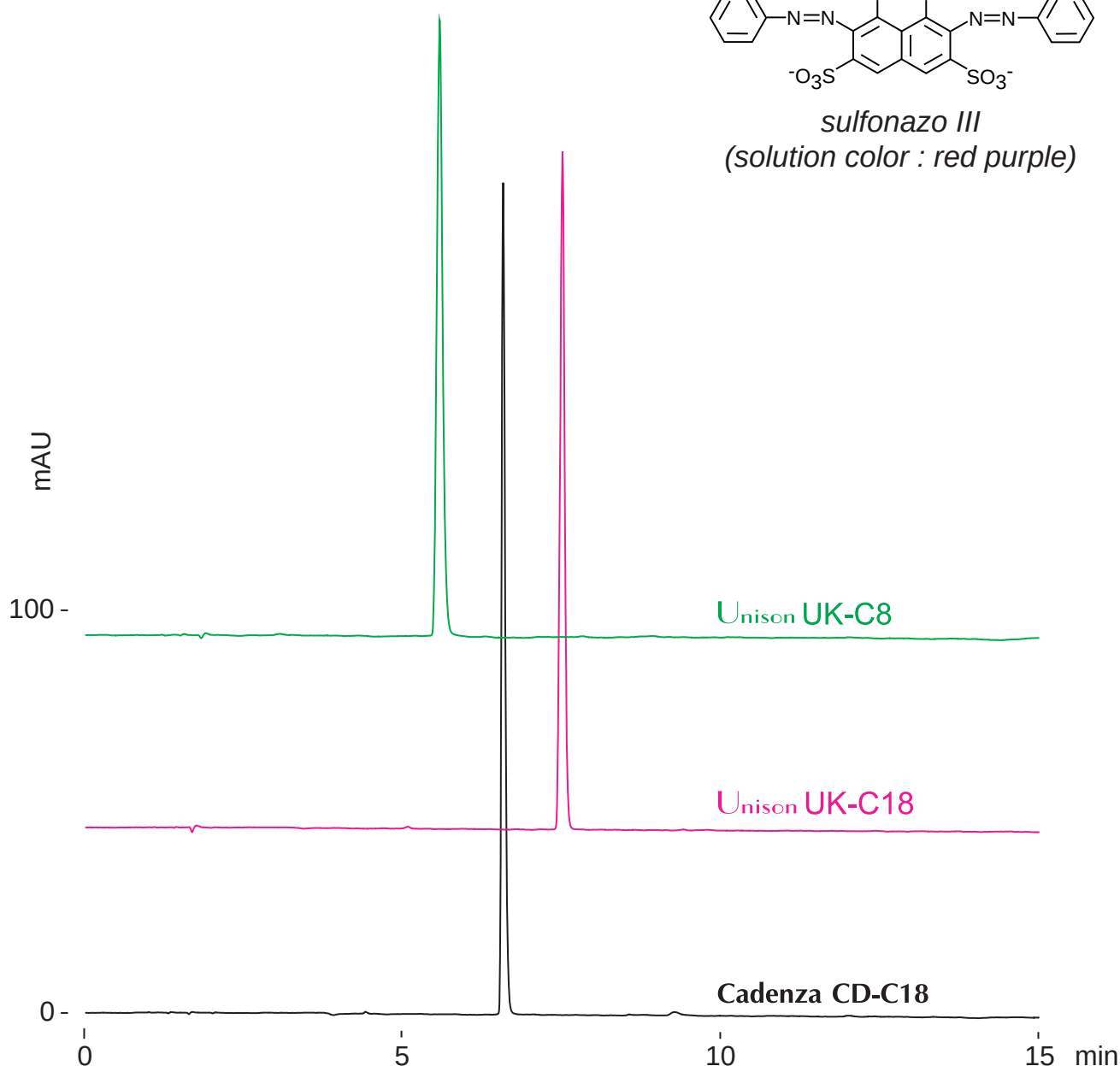

100 x 4.6 mm

A: 100mM ammonium formate / acetonitrile = 98 / 2

B: 100mM ammonium formate / acetonitrile = 50 / 50

0%B (0-1min), 0-40%B (1-15min)

1.0mL/min (11-14MPa), 37°C, 560nm

10μL (0.1mg/mL, 100mM ammonium formate)