

Unison UK-Amino



Nardis ND-NX

150 x 4.6 mm

HPLC COLUMN APPLICATION

Chlormequat (Plant growth retardant)

クロルメコート (植物成長調整剤)

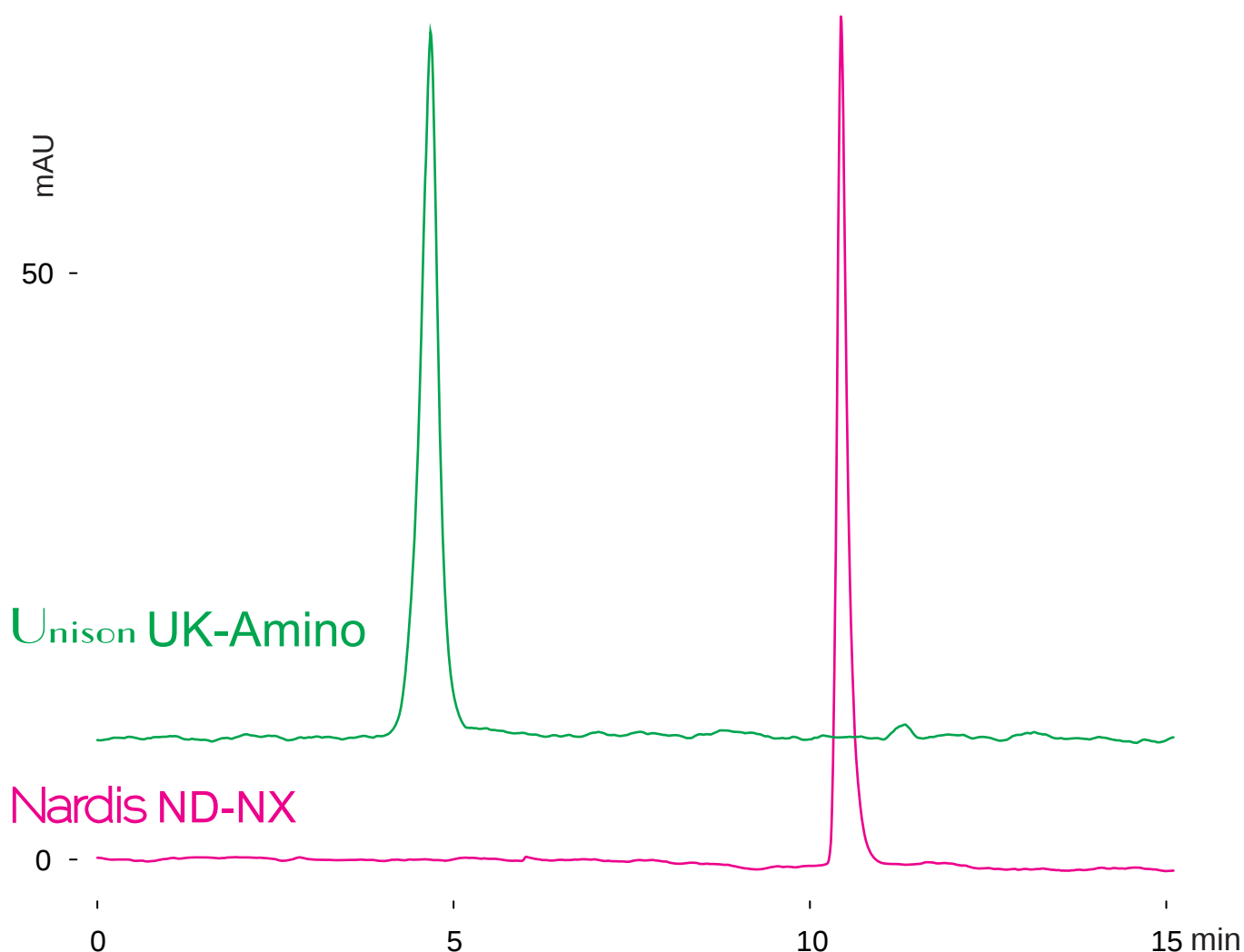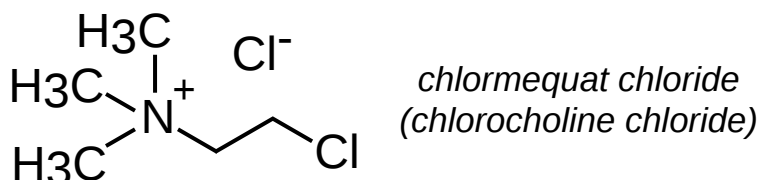

150 x 4.6 mm

A: acetonitrile

B: acetonitrile / 100mM ammonium acetate = 50 / 50

0-100%B (0-15min)

1.0mL/min (6-8MPa), 37°C

ELSD (spray chamber 70°C, drift chamber 100°C)

10μL (1mg/mL, water/acetonitrile/formic acid=50/50/0.1)