

Unison UK-Silica

Nardis ND-NX


Unison UK-Amino

150 x 4.6 mm

HPLC COLUMN APPLICATION

Oseltamivir phosphate (Anti-flu medicine)

オセルタミビル リン酸塩 (インフルエンザ治療薬)

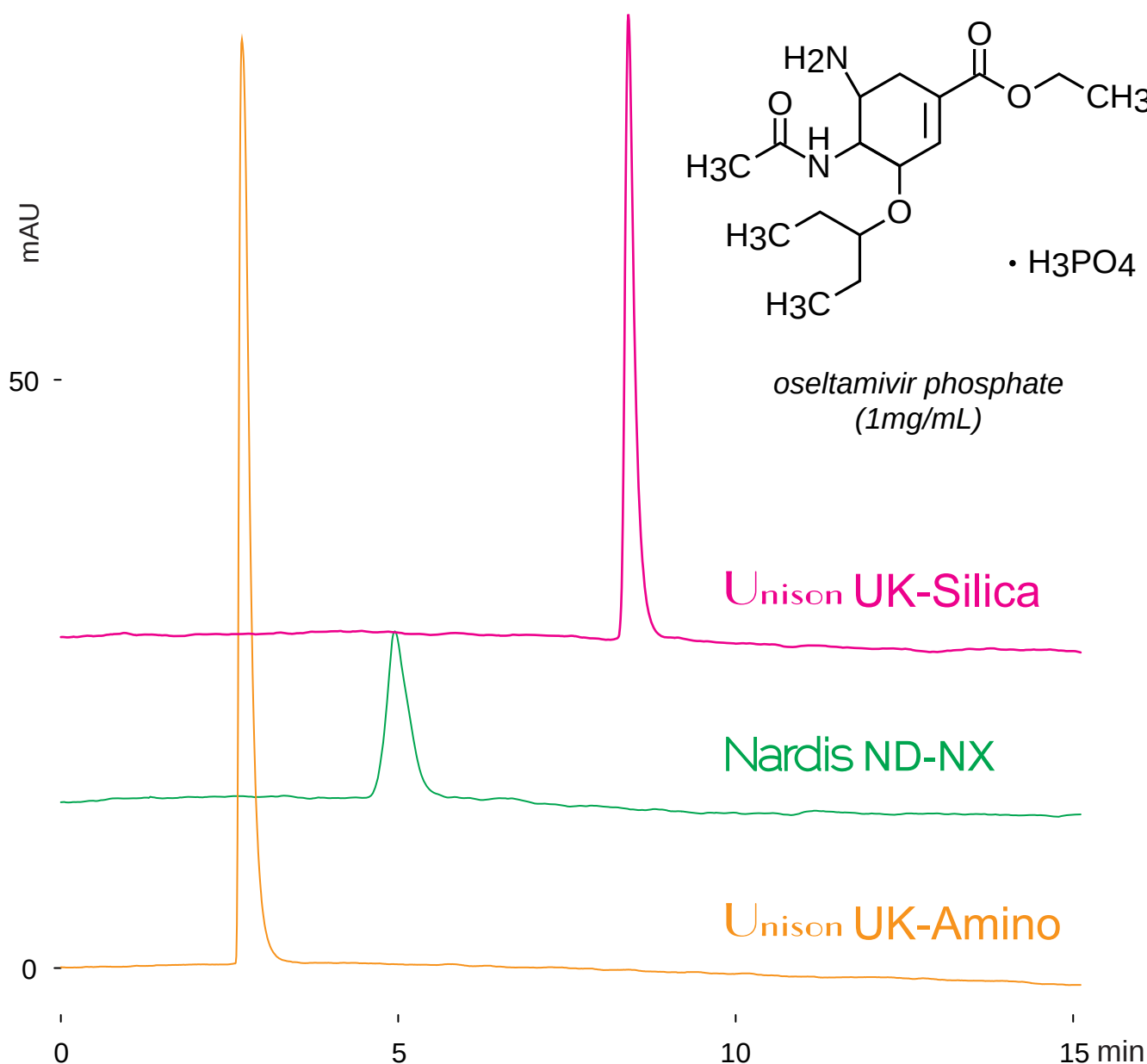

150 x 4.6 mm

A: acetonitrile / formic acid = 100 / 0.1

B: acetonitrile / water / formic acid = 50 / 50 / 0.5

0-100%B (0-15min)

1.0mL/min (5-7MPa), 37°C

ELSD (spray chamber 70°C, drift chamber 100°C)

10μL (1mg/mL, water)