

Unison UK-C1 Cadenza CW-C18

Unison UK-Phenyl Unison UK-C18


Unison UK-C8 Cadenza CD-C18

100 x 4.6 mm


HPLC COLUMN APPLICATION

2',7'-Dichlorofluorescein (Precipitation titration indicator)

2',7'-ジクロロフルオレセイン (沈殿滴定指示薬)

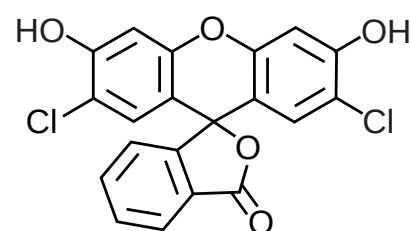

2',7'-dichlorofluorescein
(solution color : yellow)

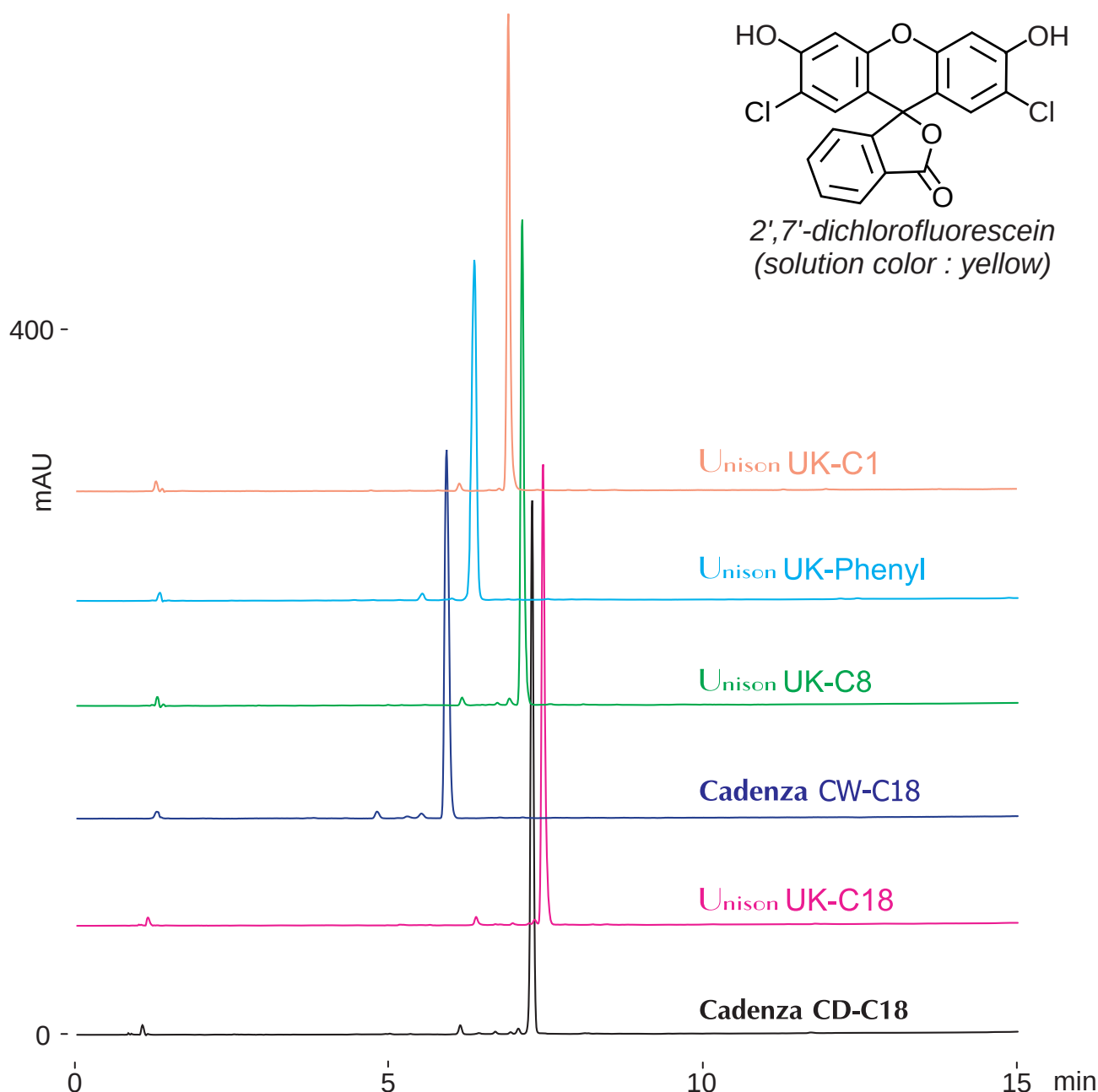

100 x 4.6 mm

A: water / acetonitrile / formic acid = 90 / 10 / 0.1

B: water / acetonitrile / formic acid = 10 / 90 / 0.1

25-100%B (0-15min)

1.0mL/min (9-16MPa), 37°C, 260nm

10μL (0.3mg/mL, water / acetonitrile / formic acid = 40 / 60 / 1)