

Unison UK-C1 Cadenza CW-C18

Unison UK-Phenyl Unison UK-C18


Unison UK-C8 Cadenza CD-C18

100 x 4.6 mm


HPLC COLUMN APPLICATION

2',7'-Dichlorofluorescein (Precipitation titration indicator)

2',7'-ジクロロフルオレセイン (沈殿滴定指示薬)

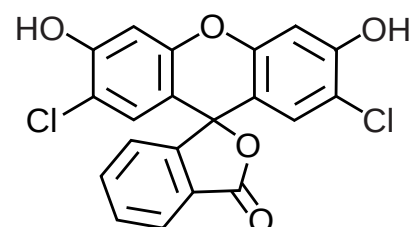

2',7'-dichlorofluorescein
(solution color : yellow)

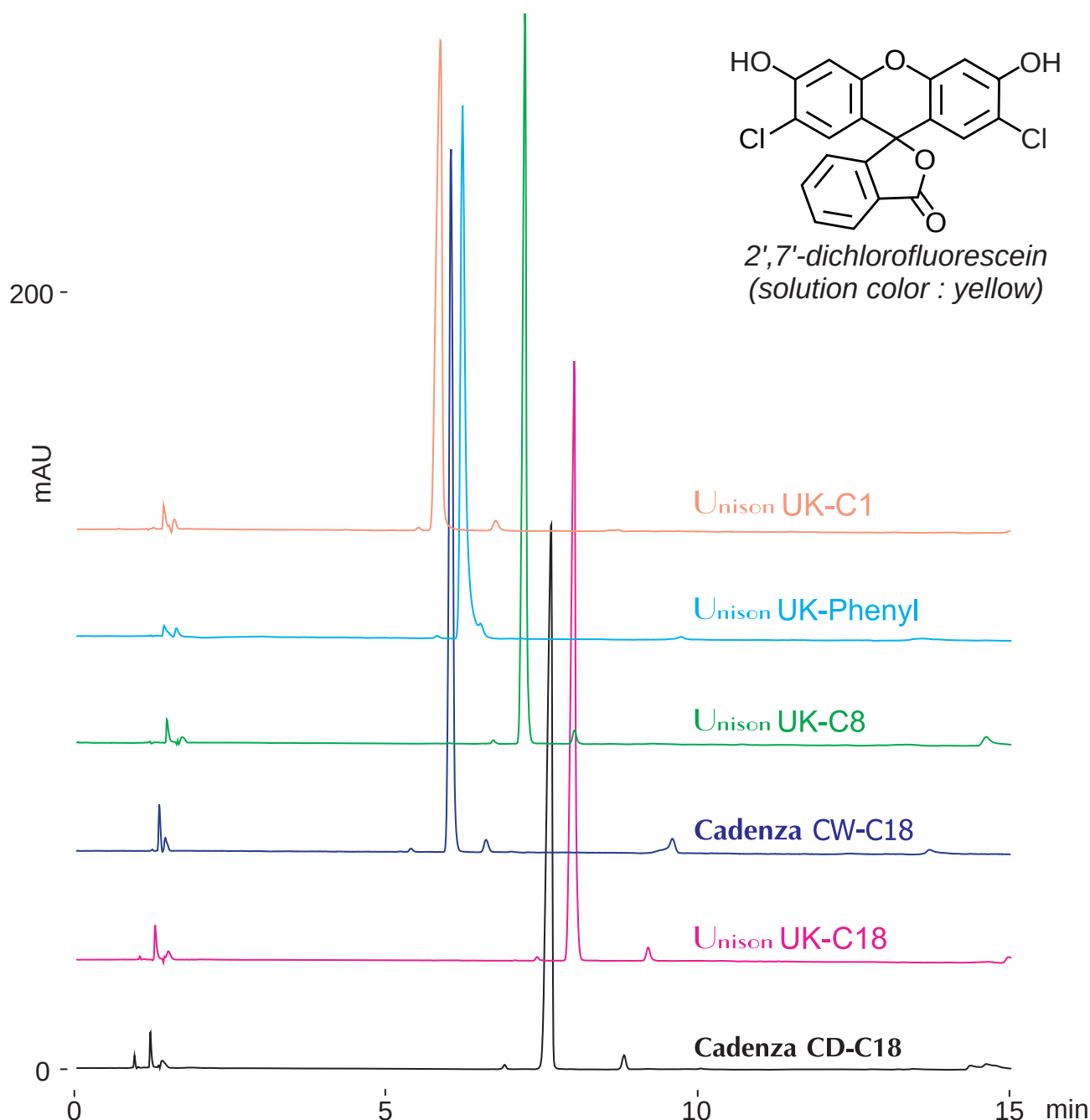

100 x 4.6 mm

A 100mM ammonium formate / acetonitrile = 90 / 10

B: 100mM ammonium formate / acetonitrile = 50 / 50

0-70%B (0-15min)

1.0mL/min, 37°C, 260nm

10μL (0.1mg/mL, water / acetonitrile / formic acid = 40 / 60 / 1)