

Cadenza CD-C18 Unison UK-C8


Unison UK-C18 Unison UK-Phenyl

Cadenza CW-C18 Unison UK-C1 100 x 4.6 mm


HPLC COLUMN APPLICATION

Losartan (Antihypertensive)

ロサルタン (高血圧治療薬)

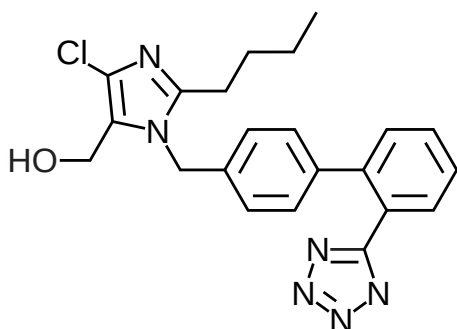

losartan

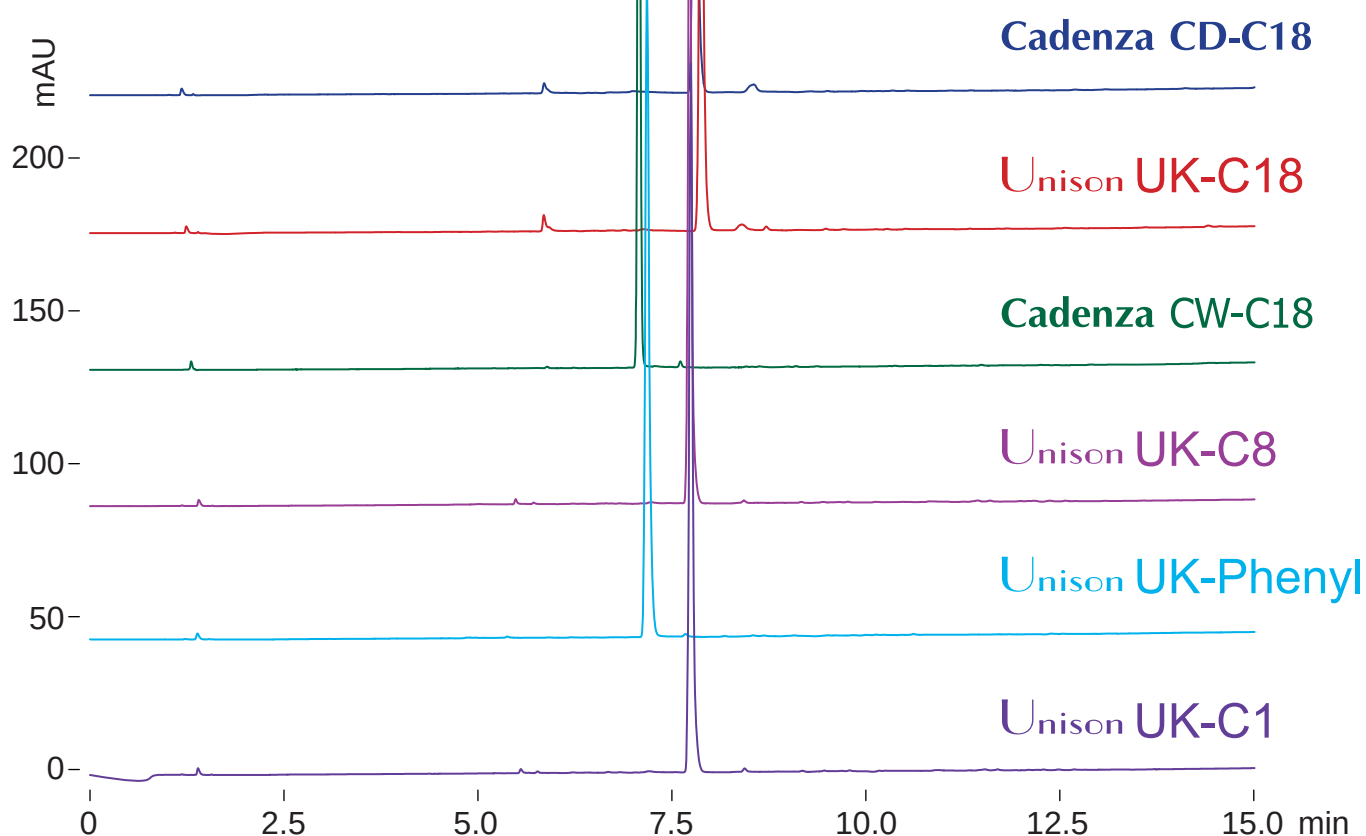

100 x 4.6 mm

A: acetonitrile / water / formic acid = 10 / 90 / 0.1

B: acetonitrile / water / formic acid = 95 / 5 / 0.1

0-100%B (0-15min)

1mL/min (10-13MPa), 37°C, 260nm

10μL (100μg/mL, acetonitrile / water / formic acid = 40 / 60 / 1)