

Cadenza CD-C18 Unison UK-C8


Unison UK-C18 Unison UK-Phenyl

Cadenza CW-C18 Unison UK-C1 100 x 4.6 mm


HPLC COLUMN APPLICATION

Losartan (Antihypertensive)

ロサルタン (高血圧治療薬)

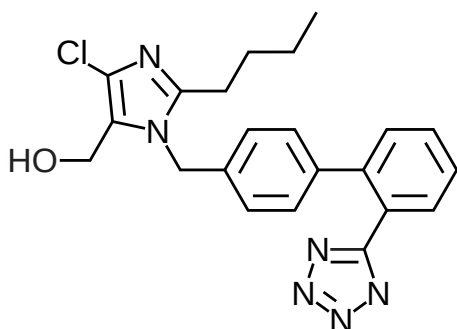

losartan

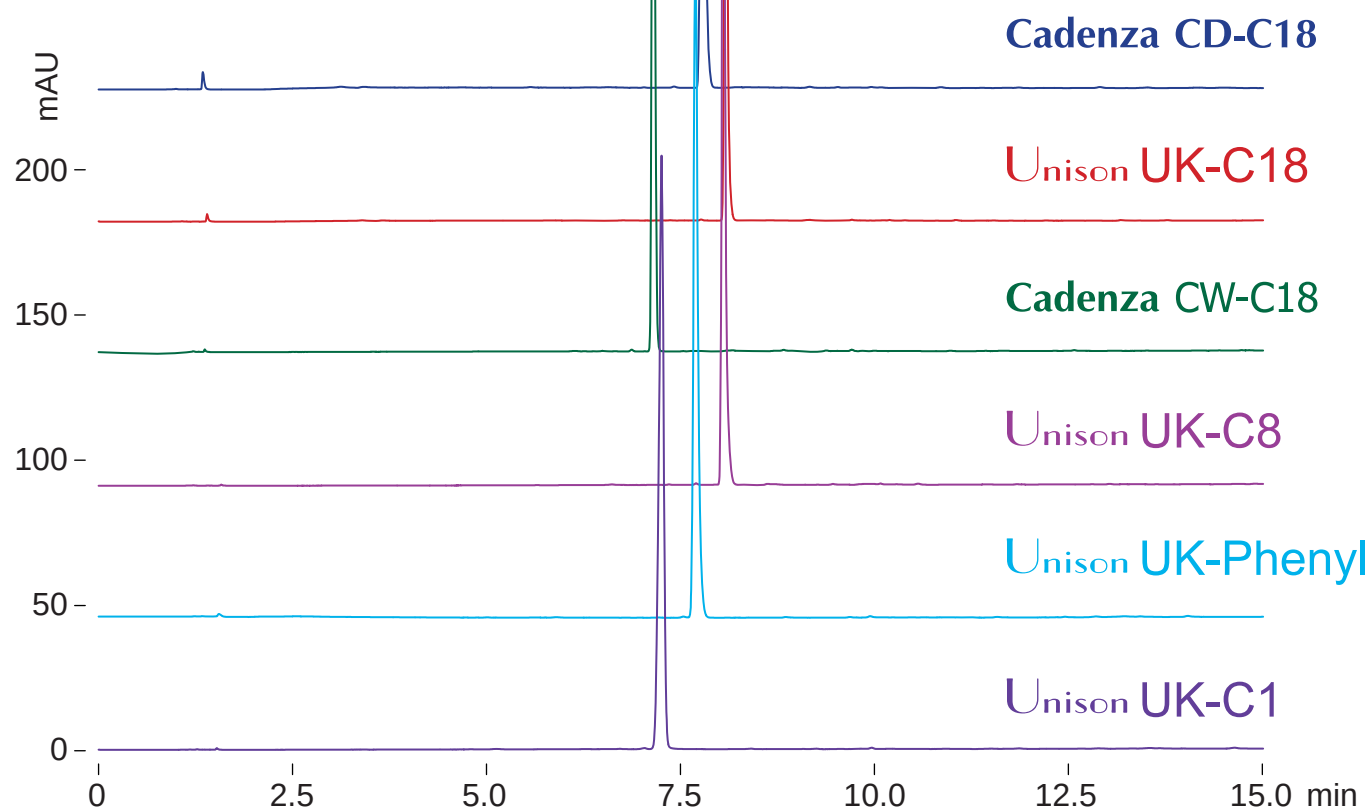

100 x 4.6 mm

A: acetonitrile / 100mM ammonium formate = 10 / 90

B: acetonitrile / 100mM ammonium formate = 95 / 5

0-70%B (0-15min)

1mL/min (10-14Pa), 37°C, 260nm

10µL (100µg/mL, ACN / 100mM HCOONH₄ = 40 / 60)